



## Laboratory Procedure Manual

*Analytes:* **Antimony, Barium, Beryllium, Cadmium, Cesium, Cobalt, Lead, Molybdenum, Platinum, Thallium, Tungsten, and Uranium**

*Matrix:* **Urine**

*Method:* **Urine Multi-Element ICP-DRC-MS**  
Renamed from "Inductively Coupled Plasma-Mass Spectrometry (ICP-DRC-MS)"

*Method No:* **ITU001B**

*Revised:* **February 16 2006**

*As performed by:* Inorganic Toxicology and Nutrition Branch  
Division of Laboratory Sciences  
National Center for Environmental Health

*Contact:* Dr. Kathleen L. Caldwell  
Phone: 770-488-7990  
Fax: 770-488-4097  
Email: [KCaldwell@cdc.gov](mailto:KCaldwell@cdc.gov)

Dr. Eric J. Sampson, Director  
Division of Laboratory Sciences

### **Important Information for Users**

The Centers for Disease Control and Prevention (CDC) periodically refines these laboratory methods. It is the responsibility of the user to contact the person listed on the title page of each write-up before using the analytical method to find out whether any changes have been made and what revisions, if any, have been incorporated.

## Procedure Change Log

Procedure: Urine Multi-Element ICP-DRC-MS DLS Method Code: ITU001B

Date	Changes Made	By	Rev'd By (Initials)	Date Rev'd
10/14/2005	RLT approach started; i.e. not diluting samples with results over the highest calibrator	JJ		
28/12/2005	"urblk" changed to "aqblk" on last line of Table 6.	JJ		
03/29/2006	Addition of CLIA cross reference index	DJK		
09May2006 ↓	<p><b><u>OLD (p.32, part 5.a.iii)</u></b>  <u>Regulator at ICP-DRC-MS:</u> Single stage "argon regulator filter kit" supplied with the ICP-DRC-MS. Set the delivery pressure to 52±1 psig. See Section 6.f. for part numbers</p> <p><b><u>New</u></b>  <u>Regulator at ICP-DRC-MS:</u> Single stage "argon regulator filter kit" supplied with the ICP-DRC-MS. Set the delivery pressure depending on the specifications for that model of ELAN ICP-DRC-MS instrument (see the PE Hardware Manual). This will be 52 ± 1 for instruments having a 0-60psi gauge and 60±1 for instruments having a 0-100psi gauge. See Section 6.f. for regulator part numbers.</p>	JJ ↓	KLC ↓	09May2006 ↓
8/14/2006	Corrected typos and small wording issues. Corrected page numbers in table of contents.	JJ		
2/8/2007	Updated section 8.b.viii (analyst evaluation of results) regarding boundary levels.	JJ	↓	03/01/2007
2/28/2007	Began using 1.5% (v/v) ethanol to both diluent and rinse solutions.	GX	↓	03/01/2007
3/5/2007	Clarified waste disposal procedures	JJ		

## Public Release Data Set Information

This document details the Lab Protocol for testing the items listed in the following table:

<b>Data File Name</b>	<b>Variable Name</b>	<b>SAS Label</b>
UHM_E	URXUBA	Barium, urine (ng/mL)
	URXUBE	Beryllium, urine (ng/mL)
	URDUCD	Cadmium, urine (ng/mL)
	URXUCO	Cobalt, urine (ng/mL)
	URXUCS	Cesium, urine (ng/mL)
	URXUMO	Molybdenum, urine (ng/mL)
	URXUPB	Lead, urine (ng/mL)
	URXUPT	Platinum, urine (ng/mL)
	URXUSB	Antimony, urine (ng/mL)
	URXUTL	Thallium, urine (ng/mL)
	URXUTU	Tungsten, urine (ng/mL)
	URXUUR	Uranium, urine (ng/mL)

## 1. Clinical Relevance & Summary of Test Principle

### A. Clinical Relevance:

This method is used to achieve rapid and accurate quantification of thirteen elements of toxicological and nutritional interest including Antimony (Sb), Arsenic (As), Barium (Ba), Beryllium (Be), Cadmium (Cd), Cesium (Cs), Cobalt (Co), Lead (Pb), Molybdenum (Mo), Platinum (Pt), Thallium (Tl), Tungsten (W), and Uranium (U). The method may be used to screen urine when people are suspected to be acutely exposed to these elements or to evaluate chronic environmental or other non-occupational exposure. [1-4].

### B. Test Principle:

Inductively coupled plasma dynamic reaction cell mass spectrometry (ICP-DRC-MS) is a multi-element analytical technique capable of trace level elemental analysis [1-4]. This ICP-DRC-MS method is used to measure the entire panel of 13 elements, or any subgroup of these.

Liquid samples are introduced into the ICP through a nebulizer and spray chamber carried by a flowing argon stream. By coupling radio-frequency power into flowing argon, plasma is created in which the predominant species are positive argon ions and electrons and has a temperature of 6000-8000 K. The sample passes through a region of the plasma and the thermal energy atomizes the sample and then ionizes the atoms. The ions, along with the argon, enter the mass spectrometer through an interface that separates the ICP (at atmospheric pressure, ~760 torr) from the mass spectrometer (operating at a pressure of  $10^{-5}$  torr). The ions pass through a focusing region, the dynamic reaction cell, the quadrupole mass filter, and finally are counted in rapid sequence at the detector allowing individual isotopes of an element to be determined. The dynamic reaction cell operates in one of two modes. In 'standard' mode the cell is not pressurized and ions pass through the cell to the quadrupole mass filter unaffected. In 'drc' mode the cell is pressurized with a gas which will collide or react with the incoming ions to either eliminate an interfering ion or change the ion of interest to a new mass which is free from interference. In this method the instrument is operated in drc mode when analyzing for cadmium and arsenic, but in standard mode when analyzing for all of the other analytes. For arsenic, the reaction cell is pressurized with a mixture of hydrogen (10%) and argon (90%) which causes the breakup of the  $^{40}\text{Ar}^{35}\text{Cl}^+$  ion which would otherwise interfere with detection of  $^{75}\text{As}$  at  $m/z$  75. When analyzing for cadmium, the reaction cell is pressurized with oxygen and the quadrupole mass filter in the reaction cell prevents  $^{98}\text{Mo}$  from entering the oxygen-rich environment of the cell. The  $^{98}\text{Mo}^{16}\text{O}^+$  ions which would normally interfere with detection of  $^{114}\text{Cd}$  at  $m/z$  114 react with the oxygen in the cell creating  $^{98}\text{Mo}^{16}\text{O}_2^+$  and  $^{98}\text{Mo}^{16}\text{O}_3^+$  at masses which no longer represent interference to  $^{114}\text{Cd}$  analysis. Since the quadrupole in the reaction cell prevents  $^{98}\text{Mo}^+$  from entering the cell, no additional interfering  $^{98}\text{Mo}^{16}\text{O}^+$  is formed in the oxygen-rich environment of the cell. Electrical signals

resulting from the detection of ions are processed into digital information that is used to indicate first the intensity of the ions and then the concentration of the element. This method was originally based on the method by Mulligan et al. [5]. The DRC portions of the method are based on work published by Tanner et al. [2, 3]. Urine samples are diluted 1+9 with 2% (v/v), double-distilled, concentrated nitric acid containing iridium (Ir), rhodium (Rh), and gallium (Ga) for multi-internal standardization. Nitric acid is used for the purpose of solubilizing and stabilizing metals in solution. Internal standards are a constant concentration in all blanks, calibrators and samples. Monitoring the instrument signal ratio of a metal to its internal standard allows correction for instrument noise and drift, and sample-to-sample matrix differences. Ethanol is used for the purpose of providing a constant amount of signal enhancement (carbon effect) for high ionization potential elements across all blanks, calibrators, and samples.

## 2. Safety Precautions

### A. General Safety

1. Observe all safety regulations as detailed in the Division (DLS) Safety Manual. Additional information can be found in your lab's chemical hygiene plan.
2. Observe Universal Precautions when working with urine.
3. Wear appropriate gloves, lab coat, and safety glasses while handling all solutions.
4. Exercise special care when handling and dispensing concentrated nitric acid. Add acid to water. Nitric acid is a caustic chemical that is capable of causing severe eye and skin damage. ***If nitric acid comes in contact with any part of the body, quickly wash the affected area with copious quantities of water for at least 15 minutes.***
5. Use secondary containment for containers holding biological or corrosive liquids.
6. The use of the foot pedal on the Micromedic Digiflex™ is recommended because it reduces analyst contact with work surfaces that have been in contact urine and also keeps the analyst's hands free to hold the specimen cups and autosampler tubes and to wipe off the tip of Micromedic Digiflex™.
7. Training will be given before operating the ICP-DRC-MS, as there are many possible hazards including ultraviolet radiation, high voltages, radio-frequency radiation, and high temperatures. This information is also detailed in the PerkinElmer ELAN® ICP-DRC-MS System Safety Manual.
8. Use flash arrestors on oxygen and argon / hydrogen gas cylinders and properly secure gas cylinders with safety harnesses.

9. Wipe down all work surfaces at the end of the day with bleach-rite spray or freshly prepared 10% (v/v) sodium-hypochlorite solution.

#### B. Radiation Safety

1. All personnel performing this method must successfully meet requirements of a CDC-OHS radiation worker (RW) due to the use of natural uranium in this method and observe all necessary radiation safety considerations indicated in the CDC Radiation Safety Manual [9].

C. Waste Disposal: Operators of this method should take the CDC-OHS Hazardous Chemical Waste Management Course (initial and yearly refreshers).

D. Waste to be Placed Into Biohazard Autoclave Bags & Pans:

1. All biological samples and diluted specimens (after analysis run).
2. All disposable plastic and paper which contact urine (autosampler tubes, gloves, etc.).
3. Used non-glass/quartz ICP-MS consumables (i.e. probes, tubing, cones, ion, and lenses).

E. Waste to be placed into Sharps Containers: Pipette Tips, broken glass or quartz instrument consumables (broken spray chambers, torches, nebulizers, etc. . .). Large broken glass which will not fit in the sharps container should be placed in a separate autoclave pan from other waste and labeled as "broken glass" (see the "Autoclaving" section of the CDC safety policies and practices manual located in the laboratory).

F. Waste to be Picked up by the Radiation Safety Office:

1. All liquid waste generated by the ICP-DRC-MS instrument is radiation waste due to the uranium content. When the waste jug is full, complete the radiation waste form, attach it to the jug, and place the jug in the designated radiation waste area for the lab. Record the waste production on the waste disposal log (kept at the waste storage location in the lab). Someone from the Radiation Safety Office (RSO) will test the waste and dispose of it appropriately.
2. Contact the laboratory radiation inventory person and the CDC Radiation Safety Office for disposal of any single element uranium standard, intermediate stock standard, or intermediate working standard solutions.

### **3. Procedures for Collecting, Storing, and Handling Specimens; Criteria for Specimen Rejection; Specimen Accountability and Tracking**

A. Procedures for Collecting, Storing, and Handling Specimens: Specimen handling conditions, special requirements, and procedures for collection and transport are

discussed in the division (DLS) Policies and Procedures Manual [8]. Copies are available in branch, laboratory, and special activities specimen-handling offices. An electronic copy is available at:

[http://intranet.nceh.cdc.gov/dls/pdf/policiesprocedures/Policy\\_and\\_Procedures\\_Manual.DLS.2002mod.pdf](http://intranet.nceh.cdc.gov/dls/pdf/policiesprocedures/Policy_and_Procedures_Manual.DLS.2002mod.pdf). In general,

1. No fasting or special diets are required before collection of urine.
2. Use sterile, lot screened collectors for specimen acquisition.
3. Urine specimens should be transported frozen (packed in dry ice during shipment is preferred when possible).
4. Once received, store long term at  $\leq -20^{\circ}\text{C}$  until time for analysis. Short-term storage at  $2-4^{\circ}\text{C}$  is acceptable. Refreeze at  $\leq -20^{\circ}\text{C}$  portions of the sample that remain after analytical aliquots are withdrawn. Thawing and refreezing samples has not been found to compromise sample results.
5. Acceptable containers for analytical aliquots include lot screened polypropylene (PP) cryovials or tubes (i.e. 5 mL cryogenic vial or 15 mL centrifuge tube).

B. Criteria for Specimen Rejection: Specimen characteristics that may compromise test results are indicated above. Reasons for rejection of a sample for analysis include

1. Low volume: Optimal amount of urine is 2 mL. The volume of urine used for one analysis is 0.5 mL.
2. Contamination: Improper collection procedures or collection devices can contaminate the urine by contact with dust, dirt, etc.

In all cases, request a second urine specimen.

C. Transfer or Referral of Specimens; Procedures for Specimen Accountability and Tracking: Location, status, and final disposition of the specimens will be tracked at least by paper document in the "Study Folder" (created before analysts receive the samples). Apart from this specimen tracking form, this folder will also contain the paper print outs of results from analysis of the specimens. Maintain records for a minimum of 3 years. Use only numerical identifiers for samples within the laboratory (e.g., case ID numbers) in order to safeguard confidentiality. Only the medical supervisor (MS) or project coordinator (PC) i.e. non CDC personnel should have access to the personal identifiers.

## **5. Procedures for Microscopic Examinations; Criteria for Rejection of Inadequately Prepared Slides**

Not applicable for this procedure.

## 6. Preparation of Reagents, Calibration (Standards), Controls, and All Other Materials; Equipment and Instrumentation

### A. Preparation of Reagent and Materials.

#### 1. Diluent

- a. Purpose: All samples (blanks, calibrators, QC, or patient samples) are combined with the diluent during the sample preparation step before analysis. This is where the internal standards are added which during the analysis will compensate for instrumental variations on the analyte signal.
- b. Contents: An aqueous solution of 10 microgram/L Rh, Ir, and Ga in 2% (v/v) double-distilled nitric acid and 1.5% (v/v) ethanol.
- c. Preparation (4L) & storage:
  - 1) Internal Standard Intermediate Mixture: Preparation of single intermediate solution containing all internal standards will simplify the addition of the internal standards into the final diluent solution. This solution can be purchased rather than prepared.
    - a) To prepare 200 mL of the Intermediate internal standard solution
    - b) Partially fill a 200 mL acid-washed volumetric flask (PP, PMP, or Teflon™) with  $\geq 18$  Mega-ohm·cm water (approximately 100-150 mL).
    - c) Carefully add 4 mL of double-distilled, concentrated nitric acid. Mix into solution.
    - d) Add 0.8 mL of 10,000 ug/mL Rh standard. If initial Rh standard concentration is different, adjust volume proportionally.
    - e) Add 0.8 mL of 10,000 ug/mL Ir standard. If initial Ir standard concentration is different, adjust volume proportionally.
    - f) Add 0.8 mL of 10,000 ug/mL Ga standard. If initial Ga standard concentration is different, adjust volume proportionally.
    - g) Fill to mark (200 mL) and mix thoroughly.
    - h) Label should include "Internal Standard Intermediate Mixture. 40 ug/mL Rh, Ir, and Ga. 2% (v/v) HNO<sub>3</sub>", "Store at room temperature", preparation date, expiration date 1 year from preparation date, and preparer's initials.
- d. Final Diluent Solution: This solution does not have to be made up in a volumetric flask. The important thing about the concentration of the internal standards is that they be consistent within all samples in one run. To prepare different volumes of diluent, add proportionally larger or smaller volumes of the solution constituents.
  - 1) Acid-rinse a 4 L container (material may be polypropylene (PP), polymethylpentene (PMP), or Teflon™).



- 2) Partially fill the 4 L container with  $\geq 18$  megaohm·cm water.
- 3) Carefully add 80 mL double-distilled, concentrated nitric acid and mix.
- 4) Carefully add 60 mL dehydrated 200 proof ethanol and mix.
- 5) Add 1 mL of the 40 ug/mL Rh, Ir, Ga internal standard solution. If other concentrations are used, the volume added should be adjusted proportionally.
- 6) Make up to volume (4 L) with  $\geq 18$  megaohm·cm water.
- 7) Store at room temperature and prepare as needed.
- 8) Label should include "10  $\mu\text{g/L}$  Rh, Ir, and Ga", "2% (v/v)  $\text{HNO}_3$ ", "1.5% (v/v) Ethanol", "Store at room temperature", preparation date, expiration date (1 year from prep), and preparer's initials.

## 2. Base Urine

- a. Purpose: This urine pool material will be mixed with the intermediate working calibrators just prior to analysis to matrix-match the calibration curve to the urine matrix of the unknown samples.
- b. Contents: A mixture of multiple urine sources collected from anonymous donors are used to approximate an average urine matrix.
- c. Preparation & Storage:
  - 1) Collect urine anonymously by placing screened containers and collection cups in the restrooms with a sign stating the reason the specimens are being collected, the name of the investigator to contact for additional information, and requesting that people provide a urine specimen (complete details can be found in CDC protocol #3994, ProTrack # DLSITN0313).
  - 2) Once the urine is collected from donors, it should be analyzed to ensure that concentrations of the analytes in this method are relatively low, so as to not interfere with the proper measurement of calibrators (see Table 2 for suggested maximum base urine concentrations).
  - 3) Once screened, mix the urine collections together in a larger container (i.e. acid washed polypropylene (PP), polymethylpentene (PMP), or Teflon™) and stir for 30+ minutes on a large stir plate (acid wash large Teflon™ stir bar before use).
  - 4) For short term storage, store at 2-4°C. For long-term storage, dispense into smaller-volume tubes (i.e., 50-mL acid-washed or lot screened polypropylene tubes) and store at  $\leq -20^\circ\text{C}$ .
  - 5) Labels on 50 mL tubes should include "Base Urine for Multi-element Method", "Store Long Term at  $\leq 20^\circ\text{C}$ ", "Store Short Term at 2-4°C", preparation date, expiration date 3 years from prep date, and preparer's initials.

## B ICP-DRC-MS Rinse Solution

1. Purpose: Pump this solution into the sample introduction system between samples to prevent carry-over of the analytes of interest from one sample measurement to the next.
2. Contents: A 0.002% Triton X-100™, 5% (v/v) double-distilled nitric acid solution and 1.5% (v/v) ethanol.
3. Preparation & Storage:
  - a. Intermediate Triton X-100 Solution: To avoid the process of dissolving pure Triton X-100 on a daily basis, prepare an intermediate 2% Triton X-100™ / 5% (v/v) double-distilled, nitric-acid solution for daily use.
    - 1) To prepare 2L of Intermediate Triton X-100 Solution:
      - a) Partially fill a 2 L acid-washed bottle (PP, PMP, or Teflon™) with  $\geq 18$  Mega-ohm·cm water (approximately 1-1.5 L). Use of volumetric flask is not required.
      - b) Add 20 mL of Triton X-100™ and stir until completely dissolved. Use a Teflon™ stir bar and stir plate if necessary (acid wash stir bar before use).
      - c) Carefully add 100 mL of double-distilled, concentrated nitric acid.
      - d) Fill to 2 L and stir thoroughly.
      - e) Label should include “2% Triton X-100™ / 5% (v/v) HNO<sub>3</sub>”, “Store at room temperature”, preparation date, expiration date 1 year from preparation date, and preparer’s initials.
    - b. Final Rinse Solution:
      - 1) To Prepare 4 L of the Final Rinse Solution:
        - a) Partially fill a 4 L acid-washed bottle (PP, PMP, or Teflon™) with  $\geq 18$  Mega-ohm·cm water (approximately 2-3 L). Use of volumetric flask is not required.
        - b) Add 4 mL of the 2% Triton X-100™ / 5% (v/v) double-distilled, nitric-acid intermediate stock solution and mix well.
        - c) Carefully add 200 mL of double distilled concentrated nitric acid and mix well.
        - d) Carefully add 60 mL dehydrated 200 proof ethanol and mix well.
        - e) Fill to 4 L using  $\geq 18$  Megaohm·cm water.
        - f) Store at room temperature and prepare as needed. To prepare volumes other than specified here, add proportionally larger or smaller volumes of the solution constituents.
        - g) Label should include “0.002% Triton X-100™ / 5% (v/v) HNO<sub>3</sub>, 1.5% (v/v) ethanol”, “Store at room temperature”, preparation date, expiration date one year from preparation date, and preparer’s initials.

## C. Standards and Calibrators

### 1. Multi-element Intermediate Stock Standard

- a. Purpose: This is the master solution from which all working calibrators will be prepared. It will be diluted to prepare intermediate working calibrators which are in turn diluted and included in each analytical run on the ICP-DRC-MS.
- b. Contents: An aqueous solution containing all 13 elements of interest for this method (does not include the internal standards). The concentrations of the 13 elements in the intermediate stock standard are listed in Table 3. The matrix is 2% (v/v) HNO<sub>3</sub> and 0.1% (v/v) HCL with traces of HF in  $\geq 18$  Mega-ohm-cm water.
- c. Preparation (Purchase) & Storage:
  - 1) Purchasing from vendors: The intermediate stock standard solution may be purchased as a custom mixture from any vendor which prepares multi-element solutions that are traceable to the National Institute for Standards and Technology (NIST) for their accuracy. Due to the uranium content, special procedures must be followed when ordering this solution.
  - 2) Current vendor & preparation process: Currently it is purchased from High Purity Standards (Charleston, SC, part number SM-2107-003). Details of the HPS preparation of the multi-element stock standard are as follows (per statement on their literature):

“Sub-boiled high purity acids were used to put the high purity metal, salts, or oxides into solution and to stabilize the standard. The solution matrix is 2% (v/v) nitric acid and 0.1% (v/v) HCl with traces of HF in  $\geq 18$  Mega-ohm-cm water. The standard was made gravimetrically by weighing the reference material to 5 significant figures. Volumetric glassware was calibrated gravimetrically to 5 significant figures.”
  - 3) In-house Preparation: If outside laboratories were not available to prepare the intermediate stock standard solution, it is also possible to make it in the laboratory from single element standards which are NIST traceable.
  - 4) Storage: Store the solution at room temperature. Due to the uranium content, and in keeping with the guidance of the CDC radiation safety manual [9], the intermediate stock standards must be kept in a lockbox. Label these bottles from HPS with additional information such as “store at room temperature”, date received, date opened, and initials of person to first open.

## 2. Multi-element Intermediate Working Standards

- a. Purpose: Use these solutions each day of analysis to prepare the final working calibrators that will be placed on the autosampler of the ELAN® ICP-DRC-MS.
- b. Content: The intermediate working standard solutions used in this method are a series of five aqueous dilutions of the multi-element intermediate

stock standard solution in 2% (v/v) double-distilled nitric acid and 1% (v/v) hydrochloric acid.

c. Preparation & Storage: To prepare different volumes of diluent, add proportionally larger or smaller volumes of the solution constituents.

1) Cleaning flasks: Acid-rinse three 100-mL, one 200-mL, one 500-mL PP, and one 2 L PP (or PMP) volumetric flasks. Check their cleanliness by comparing the counts observed on the ICP-DRC-MS for 1% (v/v) HNO<sub>3</sub> before and after contact with the flasks. Mark each of the flasks according to how they will be used. These flasks should be dedicated to this use in this method, and not used for other purposes.

2) HNO<sub>3</sub> & HCl Diluent Preparation: In the cleaned 2L flask, add 1-1.5L  $\geq 18$  Megaohm·cm water, 40 mL high purity concentrated HNO<sub>3</sub>, and 20 mL high purity concentrated HCl. Fill to the mark and mix thoroughly. Use this diluent to fill the remaining flasks during preparation of the intermediate working standards.

3) Dilutions & Storage:

a) Partially fill the 100 mL, 200 mL, and 500 mL flasks with the HNO<sub>3</sub> & HCl diluent (50-75% full).

b) Using the volumes listed (Table 4) pipette the appropriate volume of the multi-element intermediate stock standard solution into each of the five volumetric flasks. Dilute each solution to the mark with the HNO<sub>3</sub> & HCl diluent using a pipette for the final drops. Mix each solution thoroughly. The final concentrations of the 13 elements are listed in Table 4.

c) Once mixed, transfer to acid-cleaned, labeled, 50-mL containers (PP, PMP, or Teflon™) for storage. Labels should include information such as “Multi-element Urine Working Calibrators”, “2% (v/v) HNO<sub>3</sub>, 0.1% (v/v) HCl”, date of preparation, expiration date (1 year from date of preparation), “store at room temperature”, initials of preparer, and concentrations for each element.

d) Store at room temperature in a lock box (due to radiation safety protocol [9]).

### 3. Working Multi-element Calibrators

a. Purpose: The working multi-element calibrators are dilutions of the intermediate working standards. Analysis of these calibrators provides each run with a signal to concentration response curve for each analyte in the method. The concentration of an analyte in a patient urine sample dilution is determined by comparing the observed signal from the dilution of the patient urine sample to the response curve from the working multi-element calibrators.

b. Content: The working multi-element calibrators are 1:100 dilutions of the corresponding five intermediate working standards. The dilutions are

c. Preparation & Use: The working multi-element calibrators are made immediately prior to analysis when the intermediate working standards

are mixed with base urine (Section 7.b) and diluent (Section 7.a) using a Digiflex automatic pipetter. See Table 7 in section 8.b.2 for details of sample preparation.

#### 4. Internal Quality Control Materials (“Bench” QC)

- a. Purpose: Internal (or “bench”) quality control (QC) materials are used to evaluate the accuracy and precision of the analysis process, and to determine if the analytical system is “in control” (is producing results that are acceptably accurate and precise). They are included in the beginning and at the end of each analytical run.
- b. Content: The internal (or “bench”) quality control (QC) materials used in this method are pooled human urine, acidified to 1% (v/v) HNO<sub>3</sub>, and may have been spiked to reach a desired concentration. The analyte concentrations in the “low QC” are in the low-normal concentration range. The analyte concentrations in the “high QC” are in the high-normal concentration range.
- c. Preparation & Storage: Quality control materials can be either prepared by and purchased from an external laboratory or prepared within the CDC laboratories. Quality control must always be traceable to the National Institute for Standards and Technology (NIST). The CDC laboratory currently prepares its own bench QC materials using the following procedures:
  - 1) Collection of urine: Collect urine anonymously by placing screened containers and / or collection cups in the restrooms with a sign stating the reason the specimens are being collected, the name of the investigator to contact for additional information, and requesting that people provide a urine specimen (complete details can be found in CDC protocol #3994, ProTrack # DLSITN0313). Volume of urine to collect is dependent on the desired pool size. This write-up will assume a 10-L pool size for both the low and high bench QC.
  - 2) Screening Urine: Screen collected samples for metal content before mixing together to make 2 separate base urine pools (for preparing the low and high bench QC materials). Samples can be screened individually or after combining several together (reduces number of analyses).
    - a) Keep urine refrigerated whenever possible to minimize microbial growth.
    - b) Because this is only a quick screen of the metal content, the number of replicates in the urine method can be reduced to one in order to reduce analysis time.
    - c) Analyte concentrations in the final urine pool to be spiked for the low bench QC pool should be in the low-normal population range. Analyte concentrations in the final urine pool to be spiked for the high bench QC pool should be less than some pre-selected target concentration values in the high normal population range. See the Second National Report on Human Exposure to Environmental Chemicals for estimations of the

normal population ranges for metals  
(<http://www.cdc.gov/exposurereport/>).

- 3) Combining Collected Urine: Be attentive not to combine only diluted matrix urine samples into the low pool and only concentrated matrix urine samples into the high pool. The goal is for combining samples is to approach an 'average' matrix for each pool.
  - a) Graduate four acid-washed 10-L carboys (PP or PMP) in 0.5 L increments (two will be used for decanting into).
  - b) Combine collected urine samples into two separate acid-washed 10-L carboys (PP or PMP), according to their concentrations, for the low bench and high bench QC pools.
  - c) Mix each urine pool using large acid washed, Teflon™ coated stir bars and large stir plates. Keep urine refrigerated whenever possible.
  - d) Acidify each urine pool to 1% (v/v) HNO<sub>3</sub> by adding the appropriate volume of double distilled HNO<sub>3</sub>. Stir for 30+ min on large stir plates.
- 4) Settling out of solids:
  - a) Refrigerate the urine (no stirring) for 1-3 days to allow for settling out of solids.
  - b) For each urine pool, decant the urine into another of the acid-washed 10-L carboys to remove the urine from the solids settled out on the bottom of the carboy.
  - c) Repeat steps (i) and (ii) until minimal solids are left at the bottom of the carboy after sitting overnight.
- 5) Spiking of urine
  - a) Analyze a sample of each urine pool. Record these results for future recovery calculations.
  - b) Use these results to determine target analyte concentrations possible for the pools
  - c) Calculate the volume of single element standards needed to spike each pool to the desired concentrations.
  - d) While stirring the pools on large stir plates, spike each pool with calculated volumes of single element standards (all spiking standards used must be traceable to NIST).
  - e) Continue to stir pools for 30+ minutes after spiking, then reanalyze.
  - f) Repeat steps 4 and 5 until all analytes reach target concentrations keeping track of the total volume of spiking solution added to each urine pool.
- 6) Dispensing and Storage of urine

- a) Container Types: Dispense urine into lot screened containers (i.e. – 5 or 15 mL polypropylene tubes). If possible, prepare tubes of QC which have only enough volume for one typical run + 1 repeat analysis. This allows for one vial of QC to be used per day of analysis, reducing chances of contamination of QC materials due to multi-day use.
  - b) Labels: Place labels on vials after dispensing and capping if the vials are originally bagged separately from the caps. This minimizes the chance for contamination during the process. Include at least the name of QC pool (text and bar code), date of preparation, and a vial number on the labels.
  - c) Dispensing: Dispensing can be accomplished most easily using a Digiflex automatic pipetter in continuous cycling dispense mode. This process should be done in a clean environment (i.e. a class 100 clean room area or hood).
- 7) Dispensing
- a) Allow urine pool to reach room temperature before dispensing (to prevent temperature gradients possibly causing concentration gradients across the large number of vials being dispensed and to prevent condensation problems during labeling of vials). This may require leaving the carboy of urine at room temperature overnight before dispensing.
  - b) Replace the tubing attached to the dispensing syringe (left when looking at front of Digiflex) with a length of clean Teflon™ tubing long enough to reach into the bottom of the 10L carboy while it is sitting on the stir plate.
  - c) Check cleanliness of Digiflex before use by analyzing 1-2% (v/v) HNO<sub>3</sub> which has been flushed through the Digiflex with a portion of the same solution which has not been through the Digiflex.
  - d) Approximately one hour before dispensing begins,
    - i. With the large stir plate close to the left side of the Digiflex, begin stirring the urine pool to be dispensed.
    - ii. Also during this time, flush the Digiflex with urine from the pool to be dispensed. Place the ends of the tubing attached to both the sample and dispensing syringes into the carboy of urine so that urine won't be used up during this process. Be sure to secure both ends of tubing in the carboy with Parafilm so they will not come out during the flushing process.
  - e) After dispensing the urine into the vials, cap the vials and label them. Placing labels on vials after capping minimizes the chance for contamination during the process.
- 8) Homogeneity Testing: After dispensing, check homogeneity of analyte concentrations in pool aliquots by analysis of every Nth

sample dispensed (where  $N \sim 20 - 50$  depending on the pool size). Sample more heavily from the beginning and the ending portions of the tubes dispensed (these are the regions where most homogeneity problems occur). Keep samples pulled for homogeneity analysis in the sequence that they were dispensed for the purpose of looking for trends in concentrations. Once dispensed and homogeneity has been shown to be good throughout the tubes of a pool, store tubes at  $\leq -20^{\circ}\text{C}$  and pull tubes out as needed for analysis.

- 9) Storage: Urine pools should be stored long term at  $\leq -20^{\circ}\text{C}$ . Short term storage (several days) at refrigerator temperature ( $\sim 2-4^{\circ}\text{C}$ ).

## D. Instrument & Material Sources

### 1. Sources for ICP-MS Instrumentation

- a. ICP-MS: Inductively Coupled Plasma Dynamic Reaction Cell Mass Spectrometer (ELAN<sup>®</sup> 6100 DRC<sup>Plus</sup> or ELAN<sup>®</sup> DRC II) (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)).
- b. Recirculating chiller / heat exchanger for ICP-MS: Refrigerated chiller (PolyScience 6105PE for ELAN<sup>®</sup> 6100 DRC<sup>Plus</sup> instruments) or heat exchanger (PolyScience 3370 for ELAN<sup>®</sup> DRC II instruments) (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)).
- c. Autosampler: CETAC ASX 500 series autosampler. Models used have included ASX500, ASX510, and ASX510HS (Cetac Technologies, Omaha, Nebraska, [www.cetac.com](http://www.cetac.com)). See Section 6.b for source information on the PerkinElmer autosampler probe.

### 2. Sources for ICP-MS Parts & Consumables

**NOTE:** The minimum number of spares recommended before reordering (if owning one instrument) are listed as “# Spares =” in the descriptions below.

- a. Adapter, plastic: 1/4-28 female threads on one side, 1.8 mm barb adapter on the other. Connects 1/4-28 nut at flanged tubing connection to 0.045” i.d. peristaltic pump tubing. Use part # B019-3342 (“Type A” adapter, PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 4.
- b. Adapter, PEEK: Securely connects 1.6 mm O.D. PFA tubing to 0.03” I.D. peristaltic tubing. Composed of three PEEK parts.
  - 1) Female nut for 1.6 mm O.D. (1/16”) tubing. Like part P-420 (Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)).
  - 2) PEEK ferrule. Like part P-260x (10 pk SuperFlangeless ferrule, Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)).



- 3) Conical Adapter Body. Like part P-692 (Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)).
- c. Coolant, for Polyscience chiller or heat exchanger: Only PerkinElmer part # WE01-6558 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) is approved for use by PerkinElmer. # Spares = 6.
  - d. Cone, sampler (nickel): PerkinElmer part # WE021140 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Part # SC2011-Ni (Testing has also found Spectron, Ventura, CA, [www.spectronus.com](http://www.spectronus.com) cones to be comparable). # Spares = 4.
  - e. Cone, skimmer (nickel): PerkinElmer part # WE021137 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Part # SC2012-Ni (Testing has also found Spectron, Ventura, CA, [www.spectronus.com](http://www.spectronus.com) cones to be comparable) # Spares = 4.
  - f. Connector (for tubing): Use to connect 1/8" I.D. PVC tubing to 0.125" I.D. peristaltic pump tubing. Use part # 3140715 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 4.
  - g. Counterweight: Teflon counterweight with 1/16" hole through middle which holds flanged end of PFA tubing at the bottom of the rinse solution container. Use part #B0191059 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 2.
  - h. Detector, electron multiplier: Like part # N8125001 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Available direct from manufacturer (part # 14210, SGE Incorporated, Austin, Texas, <http://www.etpsci.com>) or various distributors. # Spares = 1.
  - i. Hose, for connection to chiller: Push on hose. I.D. = 1/2", O.D. = 3/4". Use part # PB-8 (per inch, Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. Do not normally need spare hose (unless moving instrument into a new location).
  - j. Hose, for exhaust of ELAN: Available as part of ELAN installation kit from Perkin Elmer (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Available direct from manufacturer as part # S-LP-10 air connector (Thermafex, Abbeville, SC, [www.thermafex.net](http://www.thermafex.net)). Equivalent part may be substituted. # Spares = 10 feet of 4" diameter and 10 feet of 6" diameter hose.
  - k. Injector, quartz with ball joint: I.D. = 2.0 mm. PerkinElmer part # WE023948 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Available direct from manufacturer as part # 400-30 (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)) or from various distributors. # Spares = 2.
  - l. Injector support (for pass-through injector): PerkinElmer part # WE023951 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Available direct from manufacturer as part # 400-37 (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)) or from various distributors. # Spares = 2.

- m. Ion Lens: PerkinElmer part # WE018034 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). # Spares = 3.
- n. Nebulizer, quartz concentric: Initial work using this method has used the standard Type A, 3 mL/min nebulizer. Alternatively, the Type C, 1mL/min nebulizer may be used to improve sensitivity and precision. The ELAN supplies 30 psi argon to the nebulizer. Variations of these nebulizers may be substituted with or without quick connects for the gas and liquid ports. Quartz nebulizers are used to avoid potential contamination from borosilicate glass (i.e. barium, uranium). # Spares = 2.
- o. Type A, Standard ELAN 3 mL/min nebulizer: PerkinElmer part # WE024371 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Available directly from manufacturer as part # TQ-30-A3 (Meinhard Glass Products, Golden, CO, [www.meinhard.com](http://www.meinhard.com)) or from various distributors. The flangeless nut and ferrule assembly has been used for liquid sample back-end connection to this nebulizer.
- p. Type C, 1 mL/min nebulizer with quick disconnects for liquid and gas ports: One example is part # 500-70QQDAC (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)). This nebulizer is designed to use quick disconnects part # 500-QD (liquid) and # 500-AC (argon).
- q. Nebulizer Connections (gas): (for nebulizer argon side-arm).
  - 1) If *not* using quick disconnection fitting, insert nebulizer argon side-arm into the 1/8" i.d. vinyl tubing and secure the connection with a hose clamp for 1/4" o.d tubing (like part # EW-06832-01, Cole Palmer Instrument Company, Vernon Hills, Illinois, [www.colepalmer.com](http://www.colepalmer.com)). # Spares = 2.
  - 2) Quick disconnection fitting: Like part # 500-AC (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)). # Spares = 2.
- r. Nebulizer Connections (liquid): (for nebulizer 4mm o.d. liquid sample backend). Can use quick disconnect or flangeless nut and ferrule assembly.
  - 1) Quick Disconnect: Like Part # 500-QD (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)). # Spares = 2.
  - 2) Flangeless nut and ferrule assembly: An assembly such as part # FIT KIT 3 (Meinhard Glass Products, Golden, CO, [www.meinhard.com](http://www.meinhard.com)) or equivalent. Individual pieces of FIT KIT #3 can be purchased as follows.
    - a) Nut, flangeless, 1/16", 1/4-28, Delrin<sup>®</sup> (Acetal), red. Part # P202x (10pk, Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)). # Spares = 10.
    - b) Ferrule, flangeless, 1/16", Tefzel<sup>®</sup> (ETFE), blue. Part # P-200x (10pk, Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)). # Spares = 10.

- c) Adapter, 1/4-28 internal to 5/16-24 internal, PEEK™. Part # P-135 (Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)). # Spares = 2.
  - d) Nut, 4mm ID PEEK. Part of fit kit 3 for concentric nebulizers. Part # S-1050 (Meinhard Glass Products, Golden, CO, [www.meinhard.com](http://www.meinhard.com)). # Spares = 2.
  - e) Ferrule, 4mm ID green Delrin. Part of fit kit 3 for concentric nebulizers. Part # S-1121 (Meinhard Glass Products, Golden, CO, [www.meinhard.com](http://www.meinhard.com)). # Spares = 2.
- s. Nut: (for flanged connections of 1.59 mm (1/16") o.d. PFA tubing) Flanged, for 1/16" o.d. tubing, 1/4-28 threads. Use part # P-406x (pkg. of 10, Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)) or equivalent. Use a Teflon-coated Viton o-ring with this nut instead of the stainless steel washer that comes with part # P-406x). # Spares = 10.
  - t. Nut: (for bottom port of autosampler rinse station) 10-32 UMC threads for 1/16" tubing. Such as part # M653x (Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)) or equivalent. # Spares = 2.
  - u. Nut and Ferrule set, 1/8" Swagelok: Such as part # SS-200-NFSET (stainless steel) or part # B-200-NFSET (brass) (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. For part numbers listed here a quantity of 1 means 1 nut, 1 front ferrule, and 1 back ferrule. Spares = 20.
  - v. Nut and Ferrule set, 1/4" Swagelok: Such as part # SS-400-NFSET (stainless steel) or part # B-400-NFSET (brass) (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. For part numbers listed here a quantity of 1 means 1 nut, 1 front ferrule, and 1 back ferrule. Spares = 20.
  - w. Oil, Welch Director Gold: For roughing pumps. Available direct from manufacturer as part # 8995G-15 (1 gallon, Welch Rietschle Thomas, Skokie, IL, [www.welchvacuum.com](http://www.welchvacuum.com)) or from various distributors. Equivalent oil may be substituted. # Spares = 4.
  - x. O-ring: (for sampler cone) PerkinElmer part # N8120511 (pkg. of 5, PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 20 o-rings.
  - y. O-ring: (for skimmer cone) PerkinElmer part # N8120512 (pkg. of 5, PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 20 o-rings.
  - z. O-ring: (for flanged connections of 1.59mm (1/16") o.d. PFA tubing) Teflon-coated Viton o-ring, i.d. = 1/16", thickness = 1/16", o.d. = 3/16". Such as part # V75-003 (O-rings West, Seattle, WA, [www.oringswest.com](http://www.oringswest.com)) or equivalent. # Spares = 20.
  - aa. O-ring: (for injector support).
    - 1) Internal o-rings: ID = 1/4", OD = 3/8", thickness = 1/16". Need 2 o-rings per injector support setup. PerkinElmer part # N8122008 (PerkinElmer,

Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent (such as part # V75-010, O-rings West, Seattle, WA, [www.oringswest.com](http://www.oringswest.com)). # Spares = 20.

2) External o-rings: ID = 3/8", OD = 1/2", thickness = 1/16". Need 2 o-rings for each injector support setup. PerkinElmer part # N8122009 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent (such as part # V75-012, O-rings West, Seattle, WA, [www.oringswest.com](http://www.oringswest.com)). # Spares = 20.

- bb. O-ring: (for inside spray chamber at nebulizer port) Such as part # 120-56 (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)). Additional o-rings can sometimes be obtained free of charge or at reduced price when acquired while purchasing spray chambers. # Spares = 20.
- cc. O-ring: (for inside of torch mount): Part # WE017284 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Do not substitute. The PerkinElmer o-ring is especially metal impregnated to minimize RF leakage though the torch mount. # Spares = 2.
- dd. Photon Stop: PerkinElmer part # WE018278 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). # Spares = 1.
- ee. Plugs, Quick Change for Roughing Pump Oil: These plugs will only work on the Varian roughing pumps which come standard on ELAN DRC II ICPMS instruments. These plugs will not fit the Leybold pumps which come standard on the ELAN DRC Plus instruments. Part # W1011013 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). No spares typically needed.
- ff. Probe, for Autosampler: Probe with screw fitting. Such as part # B3000055. (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 2.
- gg. RF coil. PerkinElmer part # WE02-1816 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 2.
- hh. Rinse Station, modified (for CETAC ASX 500 series autosampler): Order part # SP5337 (CETAC Technologies, Omaha, Nebraska, [www.cetac.com](http://www.cetac.com)) with bottom port plug not inserted (especially not glued). A special nut will be inserted in this port to allow connection between the rinse station and the rinse solution container with PFA tubing. # Spares = 2.
- ii. Screw, for Torch Mount: PerkinElmer part # WE011870. (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 3.
- jj. Spray chamber, quartz concentric: PerkinElmer part # WE025221 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. Available direct from manufacturer as part # 400-20 (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)) or from various distributors. # Spares = 2.
- kk. Torch, quartz: PerkinElmer part # N812-2006 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. Available direct from manufacturer as part # 400-10 (Precision Glass Blowing, Centennial, CO,

[www.precisionglassblowing.com](http://www.precisionglassblowing.com)) or various distributors. Damaged torches can often be repaired for substantially lower cost than purchasing a new one by companies such as Wilmad LabGlass (Buena, NJ, [www.wilmad-labglass.com](http://www.wilmad-labglass.com)) or Precision Glass Blowing (Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)). # *New Spares* = 2.

- ll. Tubing, main argon delivery to instrument: I.D. = 1/8", O.D. = 1/4". Such as part # C-06500-02 (pkg. of 100 ft, polypropylene, Fisher Scientific International, Hampton, NH, [www.fishersci.com](http://www.fishersci.com)) or equivalent. # *Spares* = 50 ft.
- mm. Tubing, PFA: I.D. = 0.5 mm, O.D. = 1.59mm (1/16"). Used to transfer liquid between probe and peristaltic tubing, between rinse solution jug and peristaltic pump tubing, between autosampler rinse station bottom port and peristaltic pump tubing and possibly used between nebulizer and peristaltic pump tubing (if quick connection is not used for liquid sample delivery)  
  
The Perfluoroalkoxy (PFA) copolymer is a form of Teflon<sup>®</sup>. Such as part # 1548 (20 ft length, Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)) or equivalent. # *Spares* = 20 ft.
- nn. Tubing, peristaltic, 0.045" i.d. (rinse station feed): Standard PVC, 2-stop (red / red) peristaltic pump tubing, i.d. = 0.045". PerkinElmer part # N0680375, (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # *Spares* = 6 packs of 12 tubes.
- oo. Tubing, peristaltic, 0.03" i.d. (sampling): Standard PVC, 2-stop (black / black) peristaltic pump tubing, i.d. = 0.03". PerkinElmer part # 09908587 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # *Spares* = 6 packs of 12 tubes.
- pp. Tubing, peristaltic, 0.125" i.d. (spray chamber drain): Standard PVC, 2-stop (black / white) peristaltic pump tubing, i.d. = 0.125". PerkinElmer part # N812-2012 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # *Spares* = 6 packs of 12 tubes.
- qq. Tubing, PVC, i.d. = 1/8", o.d. = 3/16". Used to transfer liquid 1)between spray chamber waste port and peristaltic pump and 2)between peristaltic pump and liquid waste jug  
  
Like part # 14-169-7A (pkg. of 50 ft, Fisher Scientific International, Hampton, NH, [www.fishersci.com](http://www.fishersci.com)) or equivalent. # *Spares* = 20 ft.
- rr. Tubing, Stainless Steel, o.d. = 1/8", wall thickness = 0.028": Used to connect DRC gas cylinders to ELAN DRC gas ports. Also used to replace plastic tubing in the DRC gas path within the ELAN. Like part # SS-T2-S-028-20 (20 ft, Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. *Spares* = 20 ft.
- ss. Tubing, Teflon, corrugated, 1/4" o.d.: Connects to the auxiliary and plasma gas side-arms of the torch. Part # WE015903 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). # *Spares* = 2.
- tt. Tubing, Tygon, i.d. = 3/16", o.d. = 5/16": Used to transfer liquid between rinse station drain port and liquid waste jug. Like part # EW-06409-15

(50 ft, Cole Parmer, Vernon Hills, Illinois, [www.coleparmer.com](http://www.coleparmer.com)) or equivalent. # Spares = 20 ft.

- uu. Tubing, vinyl (argon delivery to nebulizer): Vinyl Tubing, 1/8" ID x 1/4" OD. Like part # EW-06405-02 (Cole Parmer, Vernon Hills, Illinois, [www.coleparmer.com](http://www.coleparmer.com)) or equivalent. Equivalent tubing material may be substituted. # Spares = 10 ft.
- vv. Union Elbow, PTFE 1/4" Swagelok: Connects argon tubing to torch auxiliary gas sidearm. Like part # T-400-9 (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. Spares = 2.
- ww. Union Tee, PTFE, 1/4" Swagelok: Connects argon tubing to torch plasma gas sidearm and holds igniter inside torch sidearm. Like part # T-400-3 (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. Spares = 2.

b. Sources for ICP-MS Maintenance Equipment & Supplies

- 1) Anemometer: Like digital wind-vane anemometer (Model 840032, SPER Scientific LTD., Scottsdale, AZ, [www.sperscientific.com](http://www.sperscientific.com)) or equivalent. Use to verify adequate exhaust ventilation for ICP-MS (check with hoses fully disconnected).
- 2) Pan, for changing roughing pump oil: Like part # 53216 (United States Plastics Corporation, Lima, OH, [www.usplastic.com](http://www.usplastic.com)) or equivalent. # On hand = 1.
- 3) Container, to hold acid baths for glassware: Polypropylene or polyethylene containers with lids (must be large enough for torch, injector, or spray chamber submersion). May be purchased from laboratory or home kitchen supply companies. # On hand = 4.
- 4) Cotton swabs: Any vendor. For cleaning of cones and glassware.
- 5) Cutter (for 1/8" o.d. metal tubing): Terry tool with 3 replacement wheels. Like part # TT-1008 (Chrom Tech, Inc., Saint Paul, MN, [www.chromtech.com](http://www.chromtech.com)) or equivalent.
- 6) Getter Regeneration Kit: Part # WE023257 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Use this as needed (at least annually) to clean the getter in the pathway of channel A DRC gas.
- 7) Magnifying glass: Any 10x + pocket loupe for inspection of cones and other ICP-MS parts. Plastic body is preferred for non-corrosion characteristics. Like part # 5BC-42813 (Lab Safety Supply, Janesville, WI, [www.labsafety.com](http://www.labsafety.com)).
- 8) Screw Driver, for Ion Lens Removal: Screw driver with long, flexible shaft, and 2 mm ball-Allen end for removal of ion lens screws, part # W1010620. Extra 2 mm bits, part # W1010598 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)).
- 9) Toothbrush: Any vendor. For cleaning ion lens and glassware.
- 10) Ultrasonic bath: Like ULTRASONIK™ Benchtop Cleaners (NEYTECH, Bloomfield, CT, [www.neytech.com](http://www.neytech.com)) or equivalent.

c. Sources for General Laboratory Consumable Supplies

- 1) Bar Code Scanner: Like Code Reader 2.0 (Code Corporation, Draper, UT, [www.codecorp.com](http://www.codecorp.com)) or equivalent. For scanning sample IDs during analysis setup. Any bar code scanner capable of reading Code 128 encoding at a 3 mil label density can be substituted.
- 2) Carboy (for preparation of urine quality control pool and waste jug for ICPMS sample introduction system): Polypropylene 10-L carboy (like catalog # 02-960-20C, Fisher Scientific, Pittsburgh, PA, [www.fischersci.com](http://www.fischersci.com)) or equivalent. Carboys with spouts are not advised due to potential for leaking.
- 3) Containers for diluent and Rinse Solution: Two liter Teflon™ containers (like catalog# 02-923-30E, Fisher Scientific, Pittsburgh, PA., [www.fishersci.com](http://www.fishersci.com)) and 4L polypropylene jugs (like catalog# 02-960-10A, Fisher Scientific, Pittsburgh, PA, [www.fishersci.com](http://www.fishersci.com)) have both been used. Acid rinse before use. Equivalent containers may be substituted.
- 4) Cups for urine collection: Like polypropylene 4.5 oz cup, catalog # 354013 (Becton Dickinson Labware, Franklin Lakes, NJ, [www.bd.com](http://www.bd.com)) or equivalent. Each lot of cups used must be lot screened (tested to be free of trace metal contamination). Clear plastics tend to have lowest trace metal contamination.
- 5) Gloves: Powder-free, low particulate nitrile (like Best CleaN-DEX™ 100% nitrile gloves, any vendor). Equivalent nitrile or latex gloves may be substituted.
- 6) Paper towels: For general lab use, any low-lint paper wipes such as KIMWIPES®EX-L Delicate Task Wipers or KAYDRY®EX-L Delicate Task Wipers (Kimberly-Clark Professional, Atlanta, GA, [www.kcprofessional.com](http://www.kcprofessional.com)). For sensitive applications in clean rooms, a wipe designed for clean room use may be desired such as the Econowipe or Wetwipe (Liberty, East Berlin, CT, [www.liberty-ind.com](http://www.liberty-ind.com)).
- 7) Pipette (for preparation of urine dilutions to be analyzed): Micromedic Digiflex-CX Automatic™ pipette equipped with 10.0-mL dispensing syringe, 2 mL sampling syringe, 0.75-mm tip, and foot pedal (Titertek, Huntsville, AL, <http://www.titertek.com/>).
- 8) Pipettes (for preparation of intermediate stock working standards & other reagents): Like Brinkmann Research Pro Electronic pipettes (Brinkmann Instruments, Inc., Westbury, NY, <http://www.brinkmann.com/home/>). 5-100 µL (catalog #4860 000.070), 20-300 µL (catalog #4860 000.089), 50-1000 µL (catalog #4860 000.097), 100-5000 µL (catalog #4860 000.100). Note: pipette catalog numbers are without individual chargers. Can purchase individual chargers (pipette catalog numbers will differ) or a charging stand that will hold four pipettes (catalog #4860 000.860). When purchasing pipette tips (epTips), purchase one or more boxes, then “reloads” for those boxes after that: 5-100 µL (box catalog # 22 49 133-4, reload catalog # 22 49 153-9), 20-300 µL (box catalog # 22 49 134-2, reload catalog # 22 49 154-7), 50-1000 µL (box catalog # 22 49 135-1, reload catalog # 22 49 155-5), 100-5000 µL (box catalog # 22 49 138-5, reload catalog # 22 49 198-9, bulk bag catalog # 22 49 208-0). Equivalent pipettes and tips can be substituted.

- 9) Tubes for sample analysis (for autosampler): Like polypropylene 15-mL conical tubes, BD Falcon model #352097 (Becton Dickinson Labware, Franklin Lakes, NJ, [www.bd.com](http://www.bd.com)). Equivalent tubes may be substituted which are shown by lot screening to be free of trace metal contamination. Clear plastics tend to have lowest trace metal contamination. Blue colored caps have also been used successfully for this method.
- 10) Tubes for storage of intermediate working stock standards: Like polypropylene 50-mL conical tubes, BD Falcon model #352098 (Becton Dickinson Labware, Franklin Lakes, NJ, [www.bd.com](http://www.bd.com)). For use in storage of intermediate working stock standards. Equivalent tubes may be substituted which are shown by lot screening to be free of trace metal contamination. Clear plastics tend to have lowest trace metal contamination. Blue colored caps have also been used successfully for this method.
- 11) Vortexer: Like MV-1 Mini Vortexer (VWR, West Chester, PA, [www.vwr.com](http://www.vwr.com)). Used for vortexing urine specimens before removing an aliquot for analysis. Equivalent item can be substituted.
- 12) Water purification system: Like NANOpure Diamond Ultrapure Water System (Barnstead International, Dubuque, Iowa, [www.barnstead.com](http://www.barnstead.com)). For ultra-pure water used in reagent and dilution preparations. An equivalent water purification unit capable of producing  $\geq 18$  Mega-ohm-cm water may be substituted.

#### **d. Sources of Chemicals, Gases, and Regulators**

- 1) Acid, Hydrochloric acid: Veritas™ double-distilled grade, 30-35% (GFS Chemicals Inc. Columbus, OH, [www.gfschemicals.com](http://www.gfschemicals.com)). This is referred to as “concentrated” hydrochloric acid in this method write-up. For use in preparation of intermediate working stock standards. An equivalent hydrochloric acid product may be substituted, but it must meet or exceed the purity specifications of this product for trace metals content.
- 2) Acid, Nitric acid: Veritas™ double-distilled grade, 68-70% (GFS Chemicals Inc. Columbus, OH, [www.gfschemicals.com](http://www.gfschemicals.com)). For use in diluent, rinse solution, intermediate working stock standards, and QC pool preparations. This is referred to as “concentrated” nitric acid in this method write-up. An equivalent nitric acid product may be substituted, but it must meet or exceed the purity specifications of this product for trace metals content.
- 3) Ethanol (EtOH): USP dehydrated 200 proof (Pharmco Products, Inc.) or equivalent.
- 4) Argon Gas (for plasma & nebulizer) and Regulator: High purity argon (>99.999% purity, Specialty Gases Southeast, Atlanta, GA, [www.sgsgas.com](http://www.sgsgas.com)) for torch and nebulizer. Minimum tank source is a dewar of liquid argon (180-250L). Bulk tank (1500<sup>+</sup>L is preferred).
  - a) Regulator for argon (at dewar): Stainless steel, single stage, specially cleaned regulator with 3000 psig max inlet, 0-100 outlet pressure range, CGA 580 cylinder connector, and needle valve shutoff on delivery side terminating in a ¼” Swagelok connector. Part number KPRAFPF415A2AG10



- (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)). An equivalent regulator from an alternate vendor may be substituted. # Spares = 1.
- b) Regulator for argon (between bulk tank and PerkinElmer filter regulator): Single Stage 316SS Regulator, with 0-300 psi Inlet Gauge, 0-200 psi Outlet Gauge, Outlet Spring Range, 0-250 psi, 1/4" Swagelok Inlet Connection, 1/4 turn Shut off Valve on Outlet with 1/4" Swagelok Connection and Teflon Seals. Part number KPR1GRF412A20000-AR1 (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)). An equivalent regulator from an alternate vendor may be substituted. # Spares = 1.
  - c) Regulator for argon (PerkinElmer filter regulator on back of ELAN): Argon regulator filter kit. Catalog number N812-0508 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)).
  - d) Argon / hydrogen: Argon (90%) / hydrogen (10%) for DRC channel A. Initial purity of argon = 99.9997% ("Research grade 5.7"). Initial purity of hydrogen = 99.9999% ("Research Grade 6.0"). Mixture is typically purchased in cylinder size 35 (6"x24") (Airgas South, Atlanta, GA, [www.airgas.com](http://www.airgas.com)).
  - e) Regulator for argon / hydrogen: Stainless steel, two stages, specially cleaned regulator with 3000 psig max inlet, 0-25 outlet pressure range, CGA 350 cylinder connector, and needle valve shutoff on delivery side terminating in a 1/4" Swagelok connector. Like part number KCYADPF412A2AD10 (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)). An equivalent regulator from an alternate vendor may be substituted. # Spares = 1.
  - f) Flash Arrestor (Stainless steel): Like part # 6104 (Matheson Tri Gas, Montgomeryville, PA, [www.mathesontrigas.com](http://www.mathesontrigas.com)) or equivalent.
- 5) Disinfectant, for work surfaces: Bleach-rite sprays (any distributor). On-site dilutions of bleach (1part bleach + 9 parts water) may be substituted, but must be re-made daily.
- 6) Oxygen: Oxygen ("Research Grade Research Grade 5.0", 99.9999% purity) for DRC channel B. Typically purchased in cylinder size 300 (9.5" x 54") (Airgas South, Atlanta, GA, [www.airgas.com](http://www.airgas.com)).
- a) Regulator for oxygen: High purity brass body with monel trim, two stage regulator. Stainless steel is not used for this application due to safety concerns of working with oxygen at high pressure [10]. For one regulator, order the following parts, and ask that they be tested and assembled (Engineered Specialty Products, Kennesaw, GA, [www.espgauges.com](http://www.espgauges.com)).
    - i. Tescom part # 44-3410S24-555  
*Regulator body:* Brass bar stock, two stage, Monel trim, TFE seats, Elgiloy diaphragms, CV=0.05, 3000 psig max inlet, 1-25 psig outlet range, 1/4 FNPT inlet / outlet / gauge Vports, O<sub>2</sub> cleaned to ASTM G93 and CGA4.1.
    - ii Tescom part # 60500-3000N  
*Inlet pressure gauge:* 2" diameter, 0-3000 psig range, O<sub>2</sub> cleaned, 1/4" MNPT bottom, brass.
    - iii Tescom part # 60500-0015N

*Delivery pressure gauge:* 2" diameter, 0-15 psig range , O<sub>2</sub> cleaned, ¼" MNPT bottom, brass.

- iv Tescom part # 63842-540-B  
*NPT to CGA Adaptor:* ¼" NPT to CGA 540 adapter, brass.
- v Swagelok part # B-200-1-4:  
*Adapter:* Brass male connector, ¼" MNPT to 1/8" Swagelok (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)).

An equivalent regulator from an alternate vendor may be substituted. #  
*Spares = 1.*

- b) Flash Arrestor (brass): Like part # 6103 (Matheson Tri Gas, Montgomeryville, PA, [www.mathesontrigas.com](http://www.mathesontrigas.com)) or equivalent.
- 7) Standard, Gallium: Like 1,000 mg/L, item # PLGA2-2Y. (SPEX Industries, Inc., Edison, NJ, [www.spexcsp.com](http://www.spexcsp.com)). Used as an internal standard in diluent. Any vendor whose standards are traceable to the National Institute for Standards and Technology may be substituted. The standard must have low trace metal contamination.
  - 8) Standard, Iridium: Like 1,000 mg/L iridium, item # PLIR3-2Y (SPEX Industries, Inc., Edison, NJ, [www.spexcsp.com](http://www.spexcsp.com)). Used as an internal standard in diluent. Any vendor whose standards are traceable to the National Institute for Standards and Technology may be substituted. The standard must have low trace metal contamination.
  - 9) Standard, Multi-element intermediate stock standard: Item number SM-2107-003 (High Purity Standards, Charleston, SC, <http://www.hps.net/>). This is a custom mix solution (see Table 3 for concentrations). This solution is diluted to prepare the intermediate stock working standards, which are in turn diluted to prepare the working calibrators. This solution can be prepared in-house from NIST traceable single element stock solutions if necessary.
  - 10) Standard, Rhodium: Like 1,000 mg/L, item # PLRH3-2Y. (SPEX Industries, Inc., Edison, NJ, [www.spexcsp.com](http://www.spexcsp.com)). Used as an internal standard in diluent. Any vendor whose standards are traceable to the National Institute for Standards and Technology may be substituted. The standard must have low trace metal contamination.
  - 11) Standard, single element stock standards for preparation of urine quality control pools: National Institute of Standards and Technology (NIST) Standard Reference Materials (SRMs) 3103a (As), 3105a (Be), 3113 (Co), 3134 (Mo), 3108 (Cd), 3102a (Sb), 3111a (Cs), 3104a (Ba), 3163 (W), 3128 (Pb), 3140 (Pt), 3158 (Ti), and 3164 (U) (National Institute of Standards and Technology (NIST), Office of Standard Reference Materials, Gaithersburg, MD, [www.nist.gov](http://www.nist.gov)). Other sources of standards can be used if they are NIST traceable.
  - 12) Triton X-100™ surfactant: Like "Baker Analyzed" TritonX-100™ (J.T. Baker Chemical Co., [www.jtbaker.com](http://www.jtbaker.com)). Another source may be substituted, but it must be free of trace-metal contamination.

## 7. Calibration and Calibration-Verification Procedures

### A. Calibration Curve

Generate a simple linear calibration curve for each of the 12 elements in this method using a series of external standards whose concentrations are defined in the calibration page of the quantitative analysis method software. The ratio of the analyte isotopes listed in Table 6 versus the internal standards is calculated as the net intensities of the analytes. Blank subtraction is performed after the analyte / internal standard ratio is calculated.

### B. Calibration Verification

In order to verify that the calibration of this test system is accurate throughout reportable range, use external reference materials such as NIST SRM 2670, NIST SRM 2670A, and / or proficiency testing samples (such as those from the Center of Toxicology of Quebec) at least once every 6 months. Verification needs to be done at three points across the calibration range (low, medium, and high concentrations). If external reference materials such as these are not available for each analyte in the method, follow CLIA instructions of inter-laboratory sample exchanges.

1. Bi-annual tests as defined in the DLS Policy and Procedures manual: CLIA requires the verification of accuracy of instrument response to analyte concentration be completed at least every 6 months. NIST traceable calibrators are analyzed in each run to define this response up to the concentration of the highest calibrator in the run. To verify accuracy of instrument response at concentrations higher than the highest calibrator in each run, analyze a NIST traceable standard with very high concentrations (see Table 8 in the Appendix for concentrations) at least every 6 months. Prepare the Calibration Verification Standard for analysis just as a working calibrator is prepared. Use the "Urine Blank" as the blank when it is analyzed. If the observed concentrations for the Calibration Verification Standard are not within 10% of the target value (see Table 8 in the Appendix) the lab supervisor should be notified and the issue should be investigated. Do not substitute external reference materials (i.e. biological samples from a PT program) for the Calibration Verification Standard when performing this. Solutions needed for the Calibration Verification checks can be purchased from standards vendors (i.e. SPEX, High Purity Standards, etc) or prepared in-house from NIST traceable single element standards. Always verify that normal background levels have been re-achieved through adequate rinse time following analysis of elevated standards for calibration verification.
  - a. As-needed confirmations (per supervisor discretion): When a sample result is greater than the highest calibrator in the run, the supervisor may request that the result be confirmed in an analysis run which includes a standard or external reference material with equivalent (within 10%) or greater concentration than the sample. In order to avoid needless contamination of the instrument with high concentrations of analytes, the analyst should use the lowest appropriate calibration verification solution concentrations to meet the need.

For *infrequent* verification needs, the calibration verification stock solutions can be used to prepare verification standards to appropriate concentrations. This will, however, introduce elevated concentrations of all elements in the method to the sample introduction system. Frequent measurement of these very high concentrations can result in high background levels in the instrument which are difficult to rinse out and which may limit the ability to measure low concentrations.

For frequent verification needs (i.e. when certain studies have many elevated results on particular elements) or when a concentration higher than those shown in Table 8 needs to be verified, use NIST-traceable single element stock standards to prepare single element verification standards. This will limit the exposure of the instrument to elevated concentrations of only the elements needing verification.

Always verify that normal background levels have been re-achieved through adequate rinse time following analysis of elevated standards for calibration verification.

An external reference material (i.e. historical proficiency testing sample) can be substituted in place of the Calibration Verification Standard sample in these situations if the target value has been assigned by an external source (i.e. NIST, or the proficiency testing program).

- b. The concentration of the external reference material is within 10% or is higher than the concentration of the material you need it to confirm.
- c. There is confidence that there is no contamination of previously used external reference material.
- d. A note to file is made that this was done.
- e. If the observed concentrations are not within 10% of the target value the lab supervisor should be notified and the issue should be investigated.

## **8. Operating Procedures; Calculations; Interpretation of Results (see Section 6 for details on hardware used, including sources)**

### **A. Instrumentation & Equipment Setup:**

ICP-DRC-MS: Inductively Coupled Plasma Dynamic Reaction Cell Mass Spectrometer ELAN<sup>®</sup> 6100 DRC<sup>Plus</sup> or ELAN<sup>®</sup> DRC II.

#### **1. Modifications made to ICP-DRC-MS**

- a. Plastic tubing for between mass flow controllers and dynamic reaction cell have been replaced with stainless steel. Stainless steel tubing is preferred between the reaction gas cylinder / regulator and the back of the ICP-DRC-MS instrument.
- b. A second mass flow controller has been added (channel B).

2. Sample introduction system setup:  
*(Quartz items are preferred over glass to prevent background issues related to certain metals such as barium and uranium.)*
  - a. Concentric quartz nebulizer (quick connect arrangement for liquid and gas connections available from some vendors).
  - b. Quartz injector, 2 mm ID, ball joint end.
  - c. Quartz cyclonic spray chamber.
3. Configuration of tubing for liquid handling:  
(See Section 6.f. for part numbers and ordering details.)
  - a. Tubing for liquid sample uptake:
    - 1) Probe-to-peristaltic pump tubing: Flanged 0.5 mm x 1.59mm PFA tubing attaches to top of probe by a nut and rubber o-ring (Teflon coated). Angle-cut the opposite end of PFA tubing before inserting into end of black / black peristaltic pump tubing.
    - 2) Nebulizer-to-peristaltic pump tubing: Quick connect fitting fits inside back side of nebulizer. Angle-cut opposite end of tubing (PFA or Tefzel) before inserting into end of black / black peristaltic pump tubing. (Alternative: Hold square-cut end of 0.5 mm x 1.59mm PFA tubing against the inside tapered nebulizer capillary using a flangeless nut and ferrule assembly.)
  - b. Tubing for autosampler rinse solution:

The connections to the CETAC rinse station are modified to accept PFA tubing. This is done to avoid the use of PVC tubing in the rinse station feed line since PVC tubing can be a source of metals contamination in the tubing system.

    - 1) Peristaltic pump tubing – to – Rinse solution container: In the following order, assemble
      - a) Screw connector adaptor (PerkinElmer # B0193342). The barbed end of the adapter fits into the end of the red / red peristaltic pump tubing.
      - b) Flanged end of 0.5 mm x 1.59mm PFA tubing is connected to the screw connector adapter using a nut.
      - c) Opposite end of PFA tubing goes through an inert, plastic counterweight and is held on by flanged tubing. Place counterweight into the rinse bottle and cap the bottle with Parafilm to avoid contamination and crimping the PFA tubing.
    - 2) Peristaltic pump tubing – to – Autosampler Rinse Station: In the following order, assemble

- a) Screw connector adaptor (PerkinElmer # B0193342). Insert barbed end of the adapter into end of the red / red peristaltic pump tubing.
- b) One flanged end of 0.5 mm x 1.59mm PFA tubing is connected to the screw connector adapter using a nut.
- c) The other end of the PFA tubing should go through a nut with 10-32 threads and be flanged on the end. Use epoxy glue to glue the back side of the flange to the front of the nut (this is to prevent leakage around the microbore tubing since the flange will not be pressed against a flat surface by the nut). The nut and tubing screws into the base of the autosampler rinse station.

### 3) Autosampler Rinse Station to Waste Jug:

Connect ends of 3/16" i.d. x 5/16" o.d. PVC tubing directly onto the barbed fitting at the top of the autosampler rinse station and the top of the waste jug. Because the rinse station drains by gravity to the waste jug, the tubing must go directly down into the waste jug with no upward bends or flooding of the autosampler tray will result.

#### a) Configuration of tubing for spray chamber waste removal:

- i. Chamber-to-peristaltic pump tubing: Connect 1/8" i.d. x 1/4 inch o.d. PVC tubing directly to the waste port on the spray chamber. Connect other end of PVC tubing to the white / black peristaltic pump tubing using a tubing connector (PerkinElmer item # B3140715).
- ii. Waste Jug-to-peristaltic pump tubing: Connect 1/8" i.d. x 1/4" o.d. PVC tubing to the white / black peristaltic pump tubing using a tubing connector (PerkinElmer item # B3140715). Place the free end of the PVC tubing through the lid of the waste jug (be sure it is secure). Waste jug should be sitting in a secondary containment tray in case of overflow.

### 4. Cones used

Nickel cones have been used and tested to be comparable in performance from either PerkinElmer or Spectron.

### 5. Gases & Regulators setup:

- a. Argon: Argon stored as liquid in a dewar (180-250L) or bulk tank. Gaseous argon used for plasma and nebulizer.

- 1) Regulator for argon source (if a dewar): Keep the inlet pressure (headspace pressure of liquid argon dewar) above 100 psi. Set delivery pressure to 60-100psi to allow for pressure drop across tubing

- that stretches to the instrument. See Section 6.f. for part numbers and details.
- 2) Step down regulator (if source of argon is a bulk tank): Place this single stage regulator in the lab so that incoming argon pressure can be monitored and adjusted. Set delivery pressure to 60-100 psig. See Section 6.f. for part numbers and details.
  - 3) Regulator at ICP-DRC-MS: Single stage “argon regulator filter kit” supplied with the ICP-DRC-MS. Set the delivery pressure depending on the specifications for that model of ELAN ICP-DRC-MS instrument (see the PE Hardware Manual). This will be  $52 \pm 1$  psi for instruments having a 0-60psi gauge and  $60 \pm 1$  for instruments having a 0-100psi gauge. See Section 6.f. for regulator part numbers.
- b. Argon (90%) / hydrogen (10%) mixture for DRC channel A.
- 1) Regulator for Ar / H<sub>2</sub> gas mixture: Set delivery pressure to 5-7 psig. See Section 6.f. for part numbers and details.
  - 2) Flash arrestor: Stainless steel flash arrestor is used on outlet side of regulator. See Section 6.f. for part numbers and details.
- c. Oxygen (99.999+ %) gas for DRC channel B.
- 1) Regulator for O<sub>2</sub> gas: Set the delivery pressure = 5-7 psig. See Section 6.f. for part numbers and details.
  - 2) Flash arrestor: Brass flash arrestor is used on outlet side of regulator. See Section 6.f. for part numbers and details.
- d. Chiller / Heat Exchanger: Refrigerated chiller (for ELAN<sup>®</sup> 6100 DRC<sup>Plus</sup> instruments) or heat exchanger (for ELAN<sup>®</sup> DRC II instruments). For refrigerated chiller, set temperature control to 18°C.
- e. Computer: Dell Optiplex GX150, GX270, or GX280 have all been used. Processors used have included Pentium III (1 GHz) through Pentium IV (2.8 GHz). Recommend 512Mb - 1Gb RAM. External hard disk drive for nightly backups of data connects via USB port. Software used includes Windows XP Professional, service pack 2 and ELAN v3.3.
- f. Autosampler: CETAC ASX 500 series autosampler (500, 510, 510 HS have all been used) is substituted for the standard PerkinElmer AS-93plus autosampler. However, the corrosion resistant PerkinElmer autosampler probe is used instead of the Teflon probe from CETAC to minimize splashing of solution off the end of the probe tip. The rinse port has been modified to accept 1/16” PFA tubing at the inlet port (see Section 5.a.i.3 for details).
6. Parameters for Instrument and Method: See Table 1 for a complete listing of the instrument and method parameters.

## B. Daily Analysis of Samples

### 1. Preparation of the Analytical Equipment

*For further details on any part of this description, see the ITN Daily Startup SOP for ELAN ICPMS instruments.*

2. Power on the computer, printer, peristaltic pump, and autosampler, and log into the operating system.
3. Peristaltic pump: Set up the peristaltic pump tubing with proper tension for the sample rinse station.
4. ELAN software: Start the ELAN® ICPMS software from Windows™.
5. Daily Pre-Ignition Maintenance Checks: Perform daily maintenance checks as described in the ITN Daily Startup SOP for ELAN instruments (i.e., Ar supply pressure, interface components cleanliness and positioning, interface pump oil condition, vacuum pressure, etc.). Make appropriate notes in the Daily Maintenance Checklist and Instrument Log Book.
6. Start the Plasma: In the INSTRUMENT window of the software (or on the front of the ELAN), press the “Start” button to ignite the plasma.
7. Send Probe to Rinse Station: Through the METHOD/SAMPLING window in the software, press the “Probe” button, and then the “Go to Rinse” button to lower the autosampler probe into the rinse solution.
8. Start the peristaltic pump:
  - a. *If using an external peristaltic pump, after lighting the plasma go to the DEVICES window of the software and press the “Connect” button to establish communication between the computer and the autosampler. Next, start the peristaltic pump by pressing the appropriate arrow in the DEVICES window (make sure that the rotational direction is correct for the way the tubing is set up in the peristaltic pump). Set the pump speed to 10 rpm in the DEVICES window.*
  - b. *If using the on-board ICP-MS peristaltic pump, start the peristaltic pump by pressing the appropriate arrow in the DEVICES window (make sure that the rotational direction is correct for the way the tubing is set up in the peristaltic pump). Set the pump speed to 10 rpm in the DEVICES window.*
9. Warm-up time: Allow at least 45 minutes warm-up time for the ICP-DRC-MS after igniting the plasma. This warm-up time is for the RF generator. There will be another “Stability time” for the DRC later in this procedure.
10. Optimizations and Daily Performance Check: After this warm-up time, perform a daily performance check and any optimizations necessary (as described in the ITN Daily Startup SOP for ELANs). Include Be (m/z 9) in the daily performance check. Fill in the Daily Maintenance Checklist according to the optimization procedures performed.
  - a. Magnesium (<sup>24</sup>Mg) may have high RSDs due to the use of Triton-X100 in the rinse solution. Avoid this problem by either temporarily using non-Triton-containing rinse solution during the daily check, or repeating the daily check multiple times in succession with no rinse time between.



- 1) Saving the Files: Save new tuning (mass calibration) parameters to the file "default.tun." Save new optimization parameters (i.e., detector voltages, autolens values, nebulizer gas flow rate) to the file "default.dac." Monthly or any time large changes are made in optimization parameters, save a separate copy of these optimization files under a different name (i.e. – default\_070703.dac).

#### 11. Software setup for Analysis:

- a. Workspace (files & folders): Click on "Open Workspace" from the "File" menu. Select the workspace file "CDC\_urine multi-element.wrk" (or one customized for user preferences). Select "Review Files" from the "File" menu. Verify & set up the correct files and data directories for your analysis (See Table 1 "File Names & Directories").
- b. Samples / Batch Window: Update the window to reflect the current sample set. The only fields which need to be filled in include the autosampler location, sample identification (id), measurement action, method, sample flush time, sample flush speed, read delay time, read delay & analysis speed, wash time, wash speed. Use a bar code scanner to input data whenever possible. See Table 1 for times and speeds. Save the Sample window file and re-use it on other days by simply replacing the sample IDs for the patient samples.

- 1) DRC Stability Time: Best analyte-to-internal standard ratio stability is obtained after 1-1.5 hrs of analysis of urine samples using the DRC method. Analyze enough "dummy" urine sample dilutions prior to any DRC analysis run to fill 1-1.5 hours of analysis time (not necessary if analyzing only a subgroup of the method containing no DRC analytes). If analyzing the full set of method analytes, 10 samples will be sufficient. See Table 5 for example of setup in the Samples / Batch window.

#### 2) Urine vs. Aqueous Method Files:

- a) The difference: There are two method files for this one method (see Table 1). It is necessary to use both to accomplish each run because the current PerkinElmer software will not allow for more than one blank per method file. The ONLY DIFFERENCE between these two files is on the Sampling tab where one lists the autosampler positions of the urine blank and urine calibrators (the "urblk" method file) and the other lists the autosampler position of the aqueous blank (the "aqblk" method file).
- b) Use: The ONLY TIME when it matters which of these files is used is when the measurement action *includes* "Run blank" or "Run standards". When the measurement action is only 'run sample', it does not matter whether the "urblk" or "aqblk" method file is used. Analysts typically follow the pattern below, however, for the sake of consistency and as a reminder of which blank must be used for which type of sample. See Table 6.

- i. The “*urblk*” method file: Use to analyze the initial urine blank (blank for the calibration curve), the urine calibrators, and the urine blank checks (*urblkchk1* & *urblkchk2*) at the very beginning of the run. The urine blank method (set up for a CETAC ASX500 series autosampler, tray B) defines the urine blank in autosampler location 11 and the urine calibration standards 1-5 in autosampler locations 12-16, respectively.
  - ii. The “*aqblk*” method file must be used to analyze all QC materials and patient samples. The aqueous blank method (set up for a CETAC ASX500 series autosampler, tray B) defines the aqueous blank in autosampler location 19.
- 3) Notation of Dilutions: To designate an extra dilution of a sample, edit the sample ID to reflect the level of dilution being performed (i.e., A 1:2 dilution of sample 1 would be reflected in the sample ID “sample 1 (2x dilution)”. This sample ID will be edited during the data-import process to the database so that it is recognized as the appropriate sample. Do not use the ELAN® software to automatically correct for sample dilutions. Extra dilution is performed on urine samples whose concentration is greater than the concentrations listed in Table 8 in the Appendix (linearity of the method has been documented up to these concentrations).

#### C. Preparation of Samples for Analysis (See Table 7)

1. Thaw the frozen urine specimens; allow them to reach ambient temperature.
2. DRC stability “dummy urine matrix”. Prepare 50+mL of standard 2 or standard 3 to be analyzed for 1-1.5 hr before the beginning of the run. This can be prepared using 50 mL polypropylene tubes or a wide-mouth bottle (which can be put on the autosampler in place of one of the tube trays).
3. Set up a series of 15-mL polypropylene tubes corresponding to the number of blanks, standards, QCs, and patient samples to be analyzed.
4. Prepare the following solutions in the 15-mL falcon tubes using the Micromedic Digiflex™ (see Table 3 for a summary).
  - a. *Aqueous Blank*: Prepare two aqueous blanks consisting of 1,000  $\mu\text{L}$  of  $\geq 18$  Mega-ohm·cm water and 9,000  $\mu\text{L}$  of diluent. One will be the actual aqueous blank and the other will be a backup (“Aqueous Blank Check”) in case the original aqueous blank gets contaminated.
  - b. *Urine Blank*: Prepare two urine blank dilutions consisting of 900  $\mu\text{L}$  of base urine (same material used to prepare the urine calibration standards), 100  $\mu\text{L}$  of  $\geq 18$  Mega-ohm·cm water, and 9,000  $\mu\text{L}$  of diluent. One of these urine blanks will be the blank for the calibration standards; the other will be analyzed twice after standard 5 as *UrBlkChk1* and *UrBlkChk2*, respectively. Results from the *UrBlkChks* will be used to determine the method limit of detection.

- c. *Calibrators*: Prepare the working calibration standards as 100  $\mu\text{L}$  of the appropriate aqueous intermediate working calibration standard, 900  $\mu\text{L}$  of base urine, and 9,000  $\mu\text{L}$  of diluent.
  - d. *Patient & QC Samples*: Before taking an aliquot for analysis, mix the sample so that no particulates remain on the bottom of the tube. Prepare urine sample dilutions as 4,500  $\mu\text{L}$  of diluent and 500  $\mu\text{L}$  of the urine sample.
  - e. Cap all of the blanks, standards, and samples and mix them well. Uncap them and place them in the autosampler of the ELAN<sup>®</sup> ICPMS in the order that was entered in the Samples / Batch window of the ELAN software.
- D. Specimen Storage and Handling During Testing: Specimens may be left at room temperature during analysis in case confirmation analyses must be made. Take stringent precautions to avoid external contamination by the metals to be determined. Specimens may be stored short term at refrigerated temperatures, but should be stored long term (>4 weeks) at  $\leq -20\text{ }^{\circ}\text{C}$ .
- E. Starting the Analysis: To begin analysis, highlight (click and drag with the mouse) the table rows of the samples that should be included in the run, and then click on "Analyze Batch."
- F. Monitoring the Analysis: Initiate work in a timely manner so that the run may be monitored. Make every effort to complete analysis within the work day so that the entire run can be monitored. If it is not possible to complete the analysis by the end of the work day, the run may be left to complete itself unattended as long as appropriate planning is made for either overnight operation or Auto Stop (see below).

Monitor the analysis for the following:

1. *DRC stability (analyte / internal standard ratio stability)*

After the analysis of the DRC stability "dummy" samples, review these results to determine if sufficient stability of the analyte-to-internal standard ratio has been reached before beginning analysis. Importing data into an MS Excel template file is useful to simplify this procedure.

2. *Proper operation of the instrument.*

3. *Contaminated blanks.*

4. *Linear calibration curves.*

a. Typical correlation coefficients will be 0.999 to 1.000.

b. The ELAN software generates a "simple linear" calibration curve (using a least squares calculation) for each of the 13 elements in this method. The curves are generated using the results from analysis of the urine blank and the 5 external urine calibrators whose concentrations are defined in the Calibration tab of the Method file. Specifically, the software plots the "net intensity" (y-axis) versus the analyte concentration (x-axis). The "net intensity" is the blank subtracted *ratio* of the measured intensity for the analyte to the measured intensity of the associated internal standard and is calculated as follows:

$$\text{net intensity} = \frac{\text{Analyte Meas Intensity}_{\text{sample}}}{\text{Internal Std Meas. Intensity}_{\text{sample}}} - \frac{\text{Analyte Meas Intensity}_{\text{Blank}}}{\text{Internal Std Meas Intensity}_{\text{Blank}}}$$

5. *Bench QC results within the acceptable limits.*

If an analyte result for the beginning QC material(s) falls outside of the 99% limits, then the following steps are recommended:

- a. If a particular calibration standard is obviously in error, remake a new dilution at the Digiflex of that working calibrator, reanalyze it, and reprocess the sample analyses using this new result as part of the calibration curve.
- b. Prepare a fresh dilution of the failing QC material and reanalyze it.
- c. Prepare fresh dilutions at the Digiflex of all of the calibration standards (working urine multi-element standards) and reanalyze the entire calibration curve using the freshly prepared standards.

If these three steps do not result in correction of the out-of-control values for QC materials, consult the supervisor for other appropriate corrective actions. Do not report analytical results for runs that are not in statistical control.

6. *Good precision among replicates.*

7. *Consistent measured intensities of the internal standards.*

Some sample-to-sample variations are to be expected. However the intensities should be within a few percent of one another, and should fluctuate around an average value (not drift continuously in one direction).

8. *Elevated patient results.*

G. Records of Results: Run results will be documented daily in both electronic and paper form.

1. Electronic Records:

- a. Transfer of Results to the Laboratory Information System / Database: Transfer data electronically between computers or software to reduce errors. When keyboard entry must be used, proofread transcribed data after entry.
- b. Long-Term Storage of ELAN software files: Files used and produced by the ELAN software in analyzing samples will be backed up long term on compact disk and kept a minimum of three years.

2. Paper Records: The paper copy of the results from the run should be put into the study folder(s) and should include

- a. A summary of the calibration curve statistics.
- b. A printout of analysis of each measurement made during the run.
- c. Optional, but helpful, is a printout of the DRC stability check measurements in graphical form.
- d. On the front sheet of the printed records, write the following

- 1) Analyst initials
  - 2) Instrument ID
  - 3) Date of Analysis
  - 4) Run # for the day on this instrument
  - 5) Study ID and Group Number
  - 6) Database batch ID (Not known until the run is imported into the database)
- H. Transfer of Results to the Laboratory Database: Every analytical run performed for the analysis of patient samples should be entered into the laboratory results database unless the run is not useable for obvious reasons (i.e. the run is stopped for some reason before ending QC is analyzed, no internal standard spiked into the diluent, etc. . . ).
- I. Data Export Process (from ELAN® software to .TXT file): If the data file was not created during the initial analysis, reprocess the data of interest either with “original conditions” option, or by loading the files and folders used during the analysis. In the ELAN® ICP-DRC-MS software, select “Review Files” from the “File” menu. From this window, you must open the files and directories that were used when collecting the data of the run that you wish to export. (If the analysis has just ended, all of these files and directories will still be open.) NOTE: A second copy of the ELAN® software can be run as an Edit/Reprocess copy without affecting an ongoing analysis by the first copy of the software running in Windows. After you open the relevant files, go to the “Report” page in the METHOD window. Deselect the box that prints a paper copy of data and select the box that sends data to a file. Select the “Report Options Template” named “CDC\_Database Output.rop” and type in a report filename using a format such as “2005-0714a\_group55.txt” to designate data from analysis of group 55 from July 14, 2005, run #1. Under “Report Format”, choose the “Use Separator” option, and under the “File Write” section choose “Append.” Finally, reprocess the data of interest. (See PerkinElmer ELAN® ICPMS Software Manual.) Make sure you apply the aqueous blank to all sample and quality control material analyses.
- J. Data Import Process (from .TXT file to Microsoft Access™ database):
1. Move the .TXT file to the appropriate subdirectory on the network drive where exported data are stored. Directories for data storage are named according to instrument \ year \ month\, such as I:\Instruments\ELANDRC2A\2005\07\.
  2. Using the ITN Database Frontends, import the instrument file into the database. On the GoTo window, click on “Add Sample Results to Database”, then “Import Instrument Data File”.
  3. Enter the appropriate information to identify the instrument, assay, analysis date & time, run number, analyst, calibrator lot number and prep date used (use the “IS Lot Number” field) and study. If other than default values for Method LOD, High Calibrator, Rep Delta Limit, and units were used in the run, document what was used by clicking on the “View/Set Batch Parameters” button, changing the appropriate values, and then clicking “Back”.

4. Press the “Import” button, and then browse to the correct network folder to select the file which contains the results from the run. Select the file and click “OK”.
5. In the “Import Instrument Results” table, pressing the “Find X’s” button will show only those samples whose sample ID is not recognized as a valid QC pool ID or sample ID for this study. (Sample IDs are set up when the study is logged into the database.) Corrections to sample IDs and dilution factors can be made in this table (e.g., correction of transcription errors and adjustment for level of dilution). If samples were diluted for analysis, both the sample ID and the dilution factor need to be edited in this table before the values are transferred to the database (the Replace command under the Edit window is helpful in this case). When corrections to sample IDs are made, press the “Check IDs” button to re-evaluate the sample IDs. Any sample or analyte row marked “Not Recognized” will not be transferred to the database when the “Transfer” button is pressed. Once transferred into the database, the data should be evaluated for QC pass / fail, then set with the appropriate settings for QC accept / reject, final value status, and comment(s). See the database programmers for more detail on working in the database.

#### K. Analyst Evaluation of Run Results:

1. Bench Quality Control: After completing a run, and importing the results into the database, export the QC results to the SAS program where the run will be judged to be in or out of control. The QC limits are based on the average and standard deviation of the beginning and ending analyses of each of the bench QC pools, so it will not be possible to know if the run is *officially* accepted or rejected until it is completed.
  - a. Quality Control Rules: The SAS program applies the division QC rules to the data as follows:
    - 1) If both QC run means (low & high bench QC) are within 2Sm limits and individual results are within 2Si limits, then accept the run.
    - 2) If 1 of the 2 QC run means is outside a 2Sm limit - reject run if:
      - a) Extreme Outlier – Run mean is beyond the characterization mean +/- 4Sm
      - b) 1 3S Rule - Run mean is outside a 3Sm limit
      - c) 2 2S Rule - Both run means are outside the same 2Sm limit
      - d) 10 X-bar Rule – Current and previous 9 run means are on same side of the characterization mean
    - 3) If one of the 4 QC individual results is outside a 2Si limit - reject run if:
      - a) R 4S Rule – Within-run ranges for all pools in the same run exceed 4Sw (i.e., 95% range limit)

Note: Since runs have multiple results per pool for 2 pools, the R 4S rule is applied within runs only.

Abbreviations:

$S_i$  = Standard deviation of individual results (the limits are not shown on the chart unless run results are actually single measurements).

$S_m$  = Standard deviation of the run means (the limits are shown on the chart).

$S_w$  = Within-run standard deviation (the limits are not shown on the chart).

L. Patient Results:

1. Elevated Results:

a. Boundaries Requiring Confirmatory Measurement:

- 1) Results Greater than the First Upper Boundary (1UB): Concentrations observed greater than the “first upper boundary” (defined in the laboratory database as the “1UB”) should be confirmed by repeat analysis of a new sample preparation. The concentration assigned to the 1UB for an element is determined by study protocol but default concentrations are in Table 9 in the Appendix. Report the original result, as long as the confirmation is within 10% of the original. Continue repeat analysis until a concentration can be confirmed.
  - 2) Results Greater Than Highest Calibrator: When a sample result is greater than the highest calibrator in the run, the supervisor may request that the result be confirmed in an analysis run which includes a standard or external reference material with equivalent (within 10%) or greater concentration than the sample.
  - 3) Results Greater Than Range of Linearity Tested: Perform an extra dilution on any urine sample whose concentration is greater than those listed in Table 8 in the Appendix (the linearity of the method has been documented up to these concentrations). See Table 7 for description of sample preparation with extra dilution.
  - 4) Uranium Isotope Ratio Measurement for Elevated Uranium Concentrations: A uranium 235/238 isotope ratio analysis is performed for all urine uranium samples where the urine total uranium concentration is greater than the 2UB boundary (see Table 9).
2. Inadequate Precision in Confirmation of a Measurement: If a sample is reanalyzed to obtain a confirmation of an initially elevated result, the confirmation should be within 10% of the original result.
3. Analyst Reporting of Elevated Results: Concentrations observed greater than the “second upper boundary” (defined in the laboratory database as the “2UB”) should be reported to the QC reviewer as an “elevated result”. The concentration assigned to the 2UB for an element is determined by study protocol but default concentrations are in Table 9 in the Appendix. The analyst should report any patient results confirmed to be greater than the second upper boundary to the QC reviewer as an “elevated result”. There is no routine notification for elevated levels for the metals determined in this

method. The protocol for supervisors reporting elevated results to medical personnel is defined according to the study protocol.

4. Inadequate Precision Within One Measurement: If the range of the three replicate readings (maximum replicate concentration value - minimum replicate concentration value) for a single sample analysis is greater than the criteria listed in Table 9 in the Appendix (“>Lim Rep Delta” in the database) **and** the range of the three replicate readings is greater than 10% of the observed concentration, do not use the measurement for reporting. Repeat the analysis of the sample.
- M. Submitting final work for Review: Once results have been imported, reviewed, and set as final in the database by the analyst,
1. Submit an email to the QC reviewer informing them of the readiness of the data for final review. The email should include
    - a. Instrument ID, run Date, run number, study ID, group ID.
    - b. Any bench QC failures (include reasons if known).
    - c. Any patient sample result greater than the 2UB boundaries (see Table 9 in the Appendix).
    - d. Anything out of the ordinary about this analytical work which could have a bearing on the availability (i.e. insufficient sample to analyze), accuracy, or precision of the results.
  2. Include all items called for by the study folder cover sheet in the study folder (i.e. printouts from the ICP-MS, bench QC evaluation) together in the study folder before submitting the folder for review when analysis is complete.
- N. Overnight operation or Using Auto Stop: Make every effort to complete analysis within the work day so that the entire run can be monitored. If it is not possible to complete the analysis by the end of the work day, the run may be left to complete itself unattended as long as appropriate planning is made for either overnight operation or Auto Stop.
1. *24 hrs / day operation in DRC mode:*
    - a. To reduce startup time in the mornings, the analyst is encouraged to operate the ELAN in DRC mode 24 hrs/day during the work week. This eliminates the need for daily 45 minute RF generator warm-up, and possibly the need for DRC stability time (if the DRC gas is not off for extended periods of time before analysis). To maintain the instrument in DRC mode when not analyzing patient samples, setup multiple sample rows in the Samples / Batch window with autosampler position n zero (rinse station of autosampler) and wash time of 1800s (30 minutes). Repeat this sample row enough times to keep the instrument in analysis mode overnight (1 sample with 15 minute wash will take ~ 25 minutes).
    - b. *AutoStop:* If 24 hrs / day ELAN operation is not desired, the instrument can shut the plasma off unattended after analysis. Setup this as follows:
    - c. On the “Auto Start / Stop” tab of the Instrument window, enable the Auto Stop feature.



- d. Press the “Change” button within the Auto Stop box and set the Delayed shutdown time to 5 minutes. This will rinse the sample introduction system of urine matrix before turning off the plasma.
  - e. It will be necessary to replace the sample peristaltic pump tubing the next day since it will have been clamped shut overnight.
- O. Equipment Maintenance: Analysts are expected to follow a 4-day analysis / 1-day maintenance schedule in the laboratory.
1. ICPMS Maintenance: On the maintenance day, perform all maintenance per the Inorganic Toxicology and Nutrition Branch ELAN ICP-MS Weekly Maintenance SOP. All equipment maintenance should be documented in the instrument logbook.
  2. Data Backup: Data on the ELAN computer will be backed up via two backup routines.
    - a. Daily Backups to External Hard Drive: Automatic backups of the “elandata” directory and all subdirectories should be programmed to occur each night onto an external hard disk.
    - b. Weekly Backup to CD: Backup all files in the active “elandata” directory and all subdirectories onto one recordable compact disc during the weekly maintenance SOP. When the active “elandata” directory on the ICP-DRC-MS computer hard drive becomes too large to fit onto a single recordable compact disc, the oldest data can be removed from the computer to make it easier to backup the entire directory weekly. This can usually be done annually.
      - 1) Backup the oldest data on the hard drive to two duplicate compact disks and verify that the files on the CD are readable
      - 2) Label them with the name of the instrument, the date range of the data, the current date, your name, and “Copy 1 of 2” or “Copy 2 of 2”
      - 3) After verifying that the CDs are readable, the oldest, backed up data can be deleted from the ICP-MS computer hard drive.
      - 4) It is best to not store duplicate copies in the same location.

## 9. Reportable Range of Results

Urine multi-element values are reportable in the range between the method LOD (see Appendix, Table 8 in the Appendix) and the highest concentration verified accurate by bi-annual calibration verification tests (see Appendix, Table 8 in the Appendix). For example, if a urine cadmium value is less than the method LOD of 0.042, report it as < 0.042 µg/L). Above the highest concentration verified, extra dilutions are made of the urine sample to bring the concentration within the verified range.

## 10. Quality Control (QC) Procedures

The Inorganic Toxicology and Nutrition Branch uses the method described in this protocol for environmental and occupational health screening studies.

This analytical method uses two types of QC systems: With one type of the QC system, the analyst inserts bench QC specimens two times in each analytical run (a set of consecutive assays performed without interruption) so that judgments may be made on the day of analysis. With the other type of QC system, “blind” QC samples are placed in vials, labeled, and processed so that they are indistinguishable from the subject samples. The supervisor decodes and reviews the results of the blind specimens. With both systems, taking these samples through the complete analytical process assesses all levels of the analyte concentrations. The data from these materials are then used to estimate methodological imprecision and to assess the magnitude of any time-associated trends. The bench QC pools used in this method comprise two levels of concentration spanning the “low-normal” and “high-normal” ranges of the analyte of interest. Both of these pools are analyzed after the calibration standards are analyzed but before any patient samples are analyzed so that judgments on the calibration curves may be made before analysis of patient samples. These bench QCs should be analyzed again at the end of the run (approximately 20 patient samples total). If more patient samples are analyzed on the same calibration curve after the second run of the bench QC (after approximately 20 patient samples), both the low-normal and high-normal bench QC should be reanalyzed both before and after the additional samples. For example, the following schemes shown below are both acceptable ways to analyze more than one set of patient samples in one day. CDC typically uses scheme 2:

<p><b>Scheme 1.</b> Calibration Standards Low Bench QC High Bench QC patient samples Low Bench QC High Bench QC</p> <p>Low Bench QC High Bench QC patient samples Low Bench QC High Bench QC</p>	<p><b>Scheme 2.</b> Calibration Standards Low Bench QC High Bench QC patient samples Low Bench QC High Bench QC</p> <p>Calibration Standards Low Bench QC High Bench QC patient samples Low Bench QC High Bench QC</p>
--	--

Before bench QC materials can be used, QC limits must be established for each pool. This is done by performing statistical calculations on each pool from the concentration results observed in at least 20 characterization runs. During the 20 characterization runs, previously characterized QCs or pools with target values assigned by outside laboratories should be included to evaluate the analysis. In addition to providing QC limits, the characterization runs can also serve to establish homogeneity of the pools.

A. Precision and Accuracy (see Appendix, Table 8)

B. Method Calculations

1. Method Limit of Detection (LODs): The detection limits for elements in urine specimens are based on 3 times the concentration standard deviation of urine blanks (named UrBlkChk1 or UrBlkChk2) analyzed in at least 20 separate runs. Method LODs are re-evaluated periodically.
2. Method Limit of Quantitation (LOQ): The Division of Laboratory Sciences does not currently utilize limits of quantitation in regards to reporting limits [8].
3. QC Limits: Quality control limits are calculated based on concentration results obtained in at least 20 separate runs. It is preferable to perform separate analyses on separate days and using multiple calibrator lot numbers, instruments, and analysts to best mimic real-life variability. The statistical calculations are performed using the SAS program developed for the Division of Laboratory Sciences (DLS\_QC\_compute\_char\_stats.sas).

## 11. Remedial Action If Calibration or QC Systems Fail to Meet Acceptable Criteria

- A. Implications of QC Failures: If the division SAS program declares the run out of control” for any analyte, use the following to determine the implications on usability of the data from the run.
1. For 1 or 2 analytes: ONLY the analytes which were “out of control” are invalid for reporting from the run. Set all run results for those 1 or 2 analytes as “QC Rejected” in the database.
  2. For 3 or more analytes: All results, regardless of analyte, are invalid for reporting from the run. Set all run results for all analytes as “QC Rejected” in the database. Note in the batch comment field why all results were marked QC rejected.

## 12. Limitations of Method; Interfering Substances and Conditions

A. Interferences Addressed by This Method

1. Breakup of Argon Chloride ( $^{40}\text{Ar}^{35}\text{Cl}$ ) Interference on Arsenic ( $^{75}\text{As}$ ) Using DRC: The dynamic reaction cell of the ELAN ICP-DRC-MS is used in this method to break apart the argon chloride ( $^{40}\text{Ar}^{35}\text{Cl}$ ) interference on arsenic at m/z 75 [6] which is common to urine analysis by ICP-MS (see Section 1.b for an explanation of this process).
2. Correction & Elimination of Interferences ( $^{114}\text{Sn}$ ,  $^{98}\text{Mo}^{16}\text{O}$ ) on Cadmium ( $^{114}\text{Cd}$ ).

B. hematical Correction for Tin ( $^{114}\text{Sn}$ ) Interference:

The correction equation ( $-0.026826 \cdot \text{Sn}118$ ) is used in the “Equations” tab of the method to correct the counts observed as m/z 114 to exclude counts due to  $^{114}\text{Sn}$ .

C. Elimination of Molybdenum Oxide ( $^{98}\text{Mo}^{16}\text{O}$ ) Interference Using DRC:

The dynamic reaction cell of the ELAN ICP-DRC-MS is used in this method to eliminate interference from molybdenum oxide ( $^{98}\text{Mo}^{16}\text{O}$ ) onto cadmium at m/z 114 [7]. See Section 1.b for an explanation of this process.

D. Matrix Enhancement of Arsenic Signal:

Matrix induced signal enhancement in ICP-MS analysis from carbon on arsenic has been previously reported in the literature [17, 18]. In this method, ethanol (1.5% v/v) is added in the diluent and rinse solutions to “normalize” the arsenic signal enhancement in all blanks, calibrators, and samples.

E. Limitations of Method (Interferences Remaining in Method)

Calcium Chloride ( $^{40}\text{Ca}^{35}\text{Cl}$ ) Interference on Arsenic ( $^{75}\text{As}$ ):

It has been determined that a small interference remains at m/z 75 when the urine matrix contains **both** high chloride **and** high calcium levels [6]. Even at extreme calcium and chloride levels, this interference is has not been found to be significant (approximately 0.4  $\mu\text{g/L}$ ).

### 13. Reference Ranges (Normal Values)

In this method the 95% reference ranges (see Appendix, Table 10) for these elements in urine fall within the range of the calibrators.

### 14. Critical Call Results (“Panic Values”)

The collaborating agency with access to patient identifiers or the responsible medical officer is notified by FAX by the supervisor of any ferritin result that is  $<10$  ng/mL, which possibly represents a significant risk for iron deficiency. Copies of faxes sent concerning abnormal results are kept in a notebook by the supervisor for the duration of the study. For NHANES 1999+, Westat automatically notifies the NCHS survey physician because of several-times weekly electronic transmission of data.

### 15. Specimen Storage and Handling During Testing

Specimens may reach and maintain ambient temperature during analysis. Take stringent precautions to avoid external contamination by the metals to be determined.

### 16. Alternate Methods for Performing Test and Storing Specimens If Test System Fails:

If the analytical system fails, the analysis may be setup on other ELAN DRC instruments in the laboratory. If no other instrument is available, store the specimens at 4°C until the analytical system can be restored to functionality. If interruption longer than 4 weeks is anticipated, then store urine specimens at  $\leq -20^{\circ}\text{C}$ .

### **17. Test-Result Reporting System; Protocol for Reporting Critical Calls (If Applicable)**

Conduct test-result reporting as outlined in the DLS Policies and Procedures Manual. As stated in Section 14, the supervisor should notify the NCHS Medical Officer.

### **18. Transfer or Referral of Specimens; Procedures for Specimen Accountability and Tracking**

The designated person who receives specimen/samples delivered to Inorganic Toxicology and Nutrition Branch (ITN) sets up a "Specimen Folder." This person should fill out a tracking form and place it in the folder for the analyst performing the analysis. The form tracks location, status, and final disposition of the specimens. When the sample analysis is completed, place the tracking form in the Specimen Tracking Record Log Book located in the trace-metals library.


Use standard record-keeping means (e.g., electronic –Microsoft Access, optical disk or floppy disk) to track specimens. Maintain records for 3 years, including related Quality Assurance/QC data. Keep duplicate records in electronic or hardcopy format, and keep them off site if they are sensitive or critical. Use only numerical identifiers (e.g., case ID numbers). To safeguard confidentiality, provide only the medical supervisor (MS) or project coordinator (PC) with personal identifiers.

## 19. Summary Statistics and QC Graphs

### A. Urinary Barium

Summary Statistics for Urinary Barium by Lot


Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	162	1/25/2007	2/23/2009	0.7508	0.0505	6.7
HU-04311_UMP	162	1/25/2007	2/23/2009	4.9819	0.2168	4.4

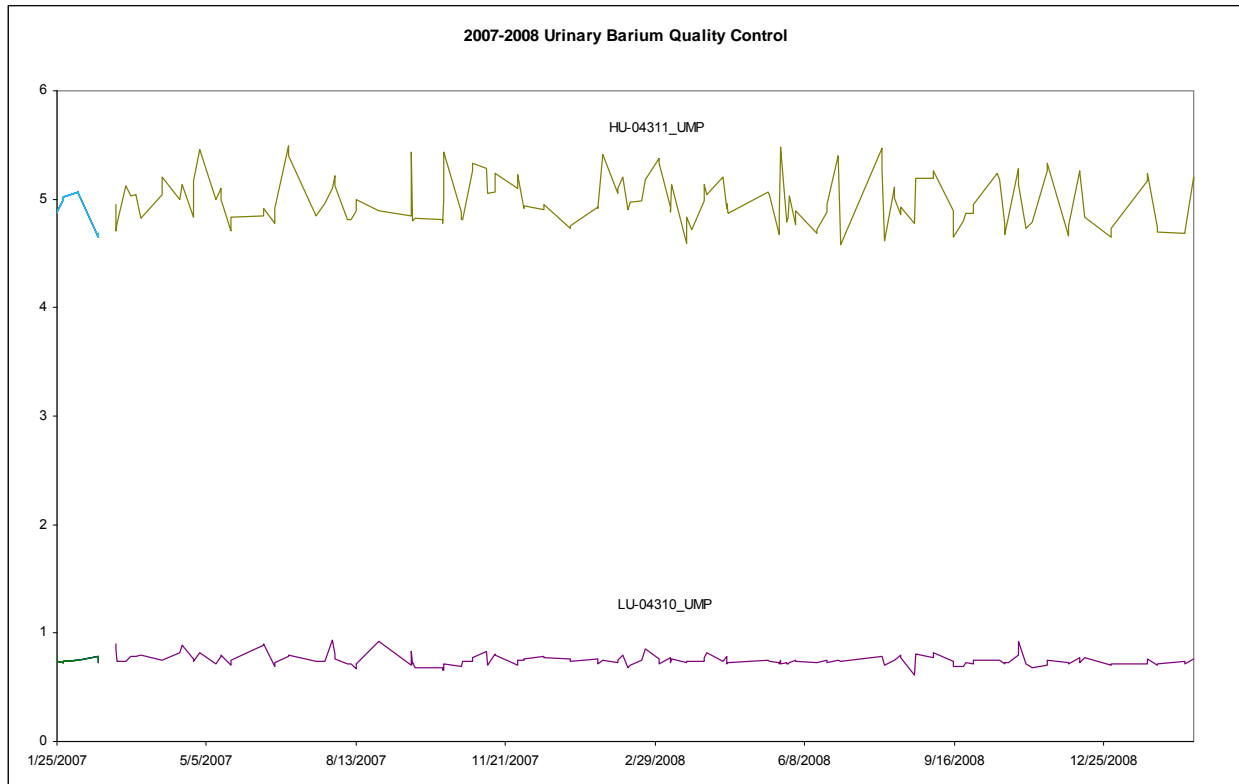

## B. Urinary Beryllium

### Summary Statistics for Urinary Beryllium by Lot


<b>Lot</b>	<b>N</b>	<b>Start Date</b>	<b>End Date</b>	<b>Mean</b>	<b>Standard Deviation</b>	<b>Coefficient of Variation</b>
LU-04310_UMP	153	1/25/2007	2/23/2009	0.6877	0.0419	6.1
HU-04311_UMP	152	1/25/2007	2/23/2009	5.3225	0.3816	7.2

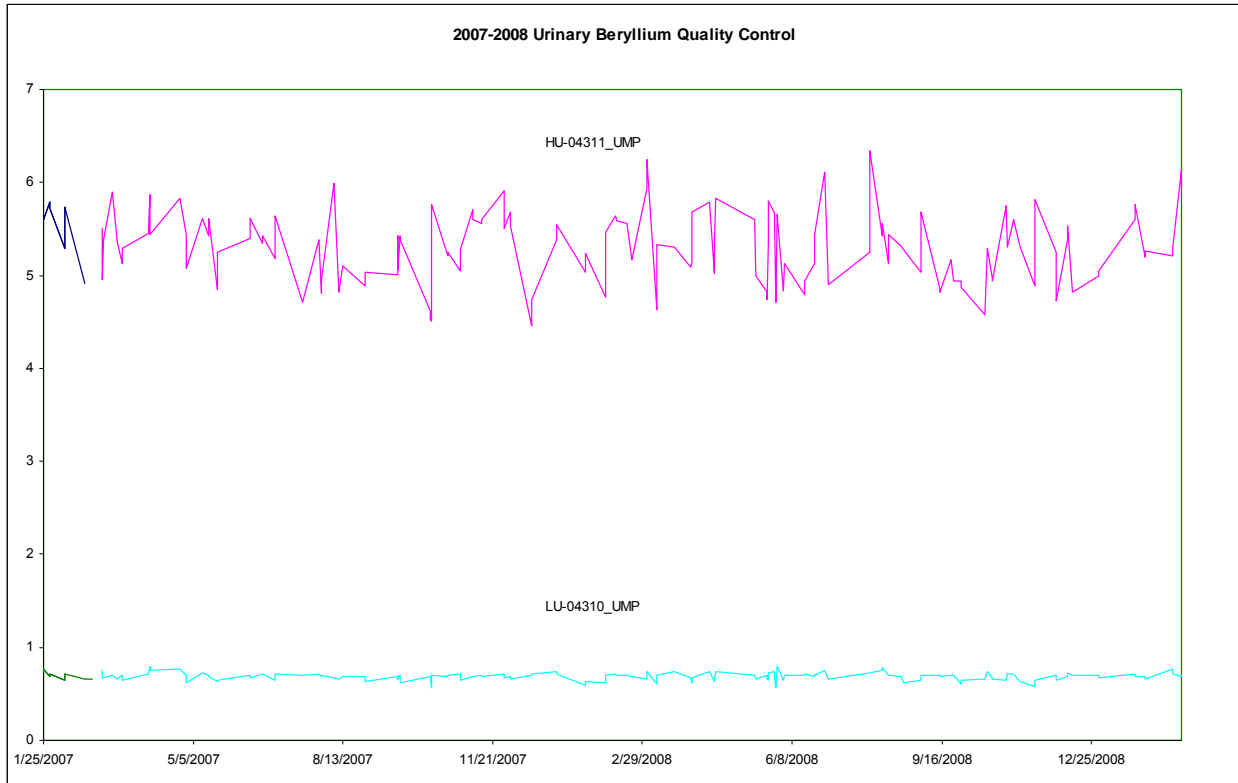

### C. Urinary Cadmium

#### Summary Statistics for Urinary Cadmium by Lot

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	160	1/25/2007	2/23/2009	0.3207	0.0152	4.8
HU-04311_UMP	160	1/25/2007	2/23/2009	1.6012	0.0666	4.2

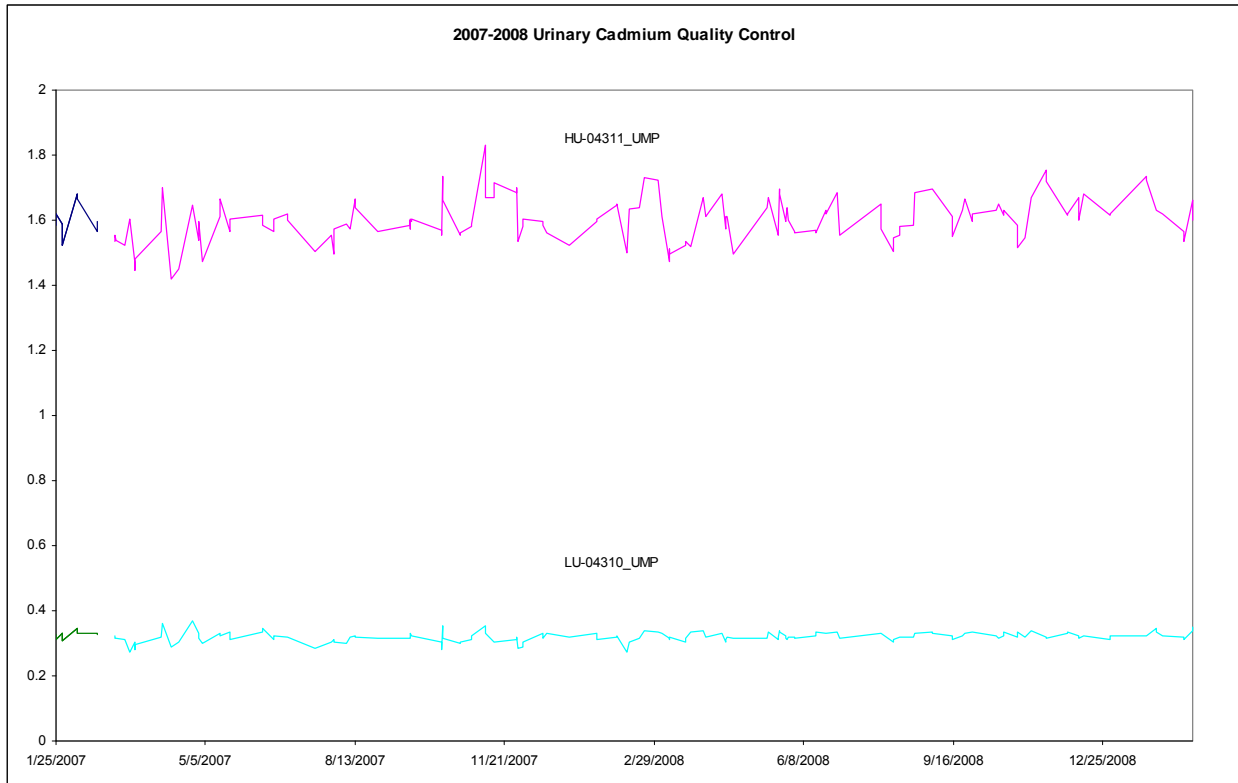




## D. Urinary Cobalt

### Summary Statistics for Urinary Cobalt by Lot


Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	162	1/25/2007	2/23/2009	0.4148	0.0249	6.0
HU-04311_UMP	161	1/25/2007	2/23/2009	1.8442	0.1067	5.8

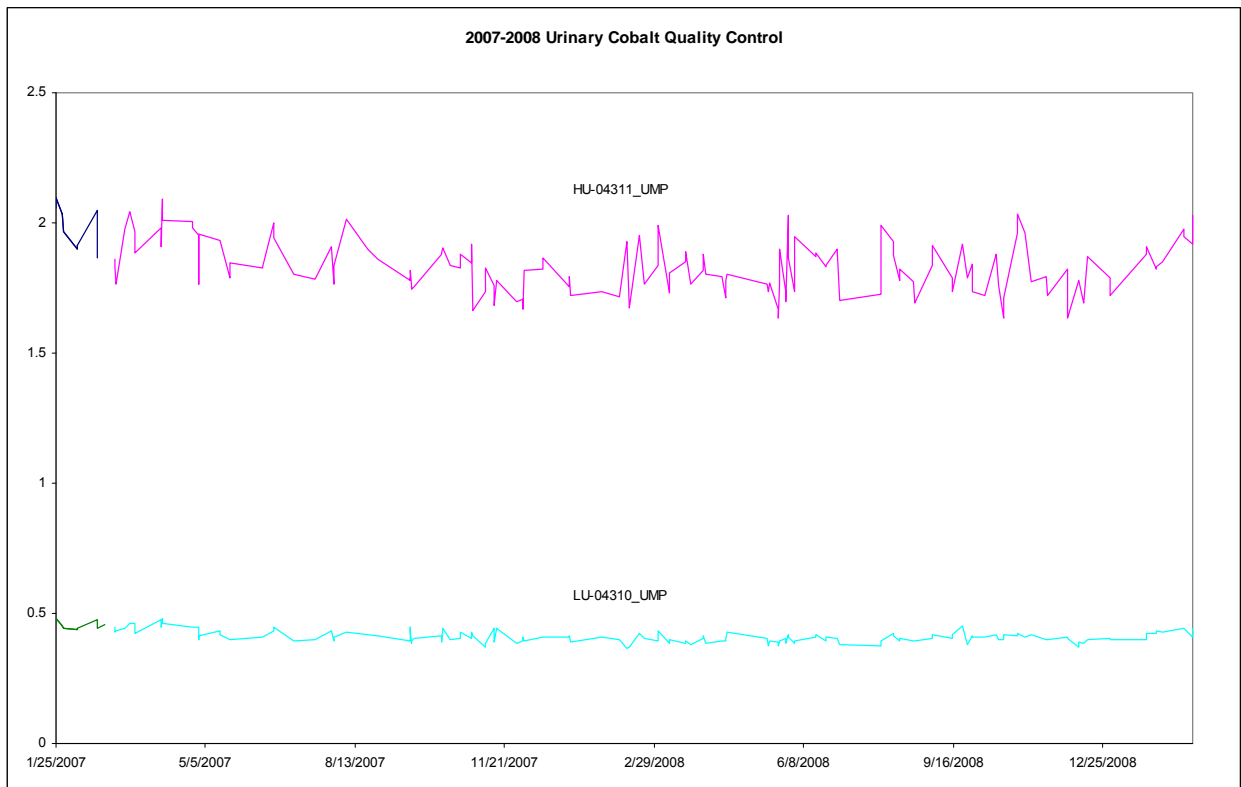

## E. Urinary Cesium

### Summary Statistics for Urinary Cesium by Lot


Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	164	1/25/2007	2/23/2009	2.3753	0.0569	2.4
HU-04311_UMP	164	1/25/2007	2/23/2009	9.7459	0.3870	4.0

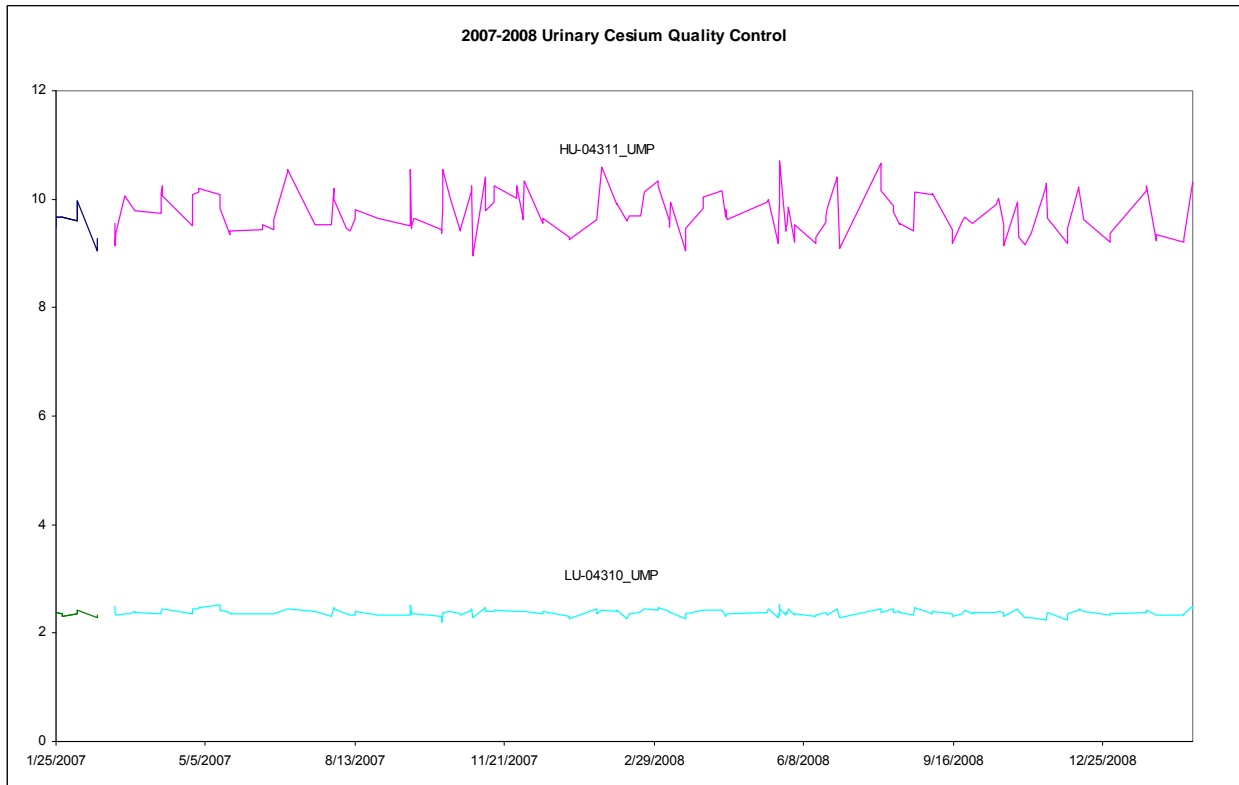

## F. Urinary Molybdenum

### Summary Statistics for Urinary Molybdenum by Lot

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	165	1/25/2007	2/23/2009	19.5361	0.3774	1.9
HU-04311_UMP	165	1/25/2007	2/23/2009	137.5939	2.5463	1.9

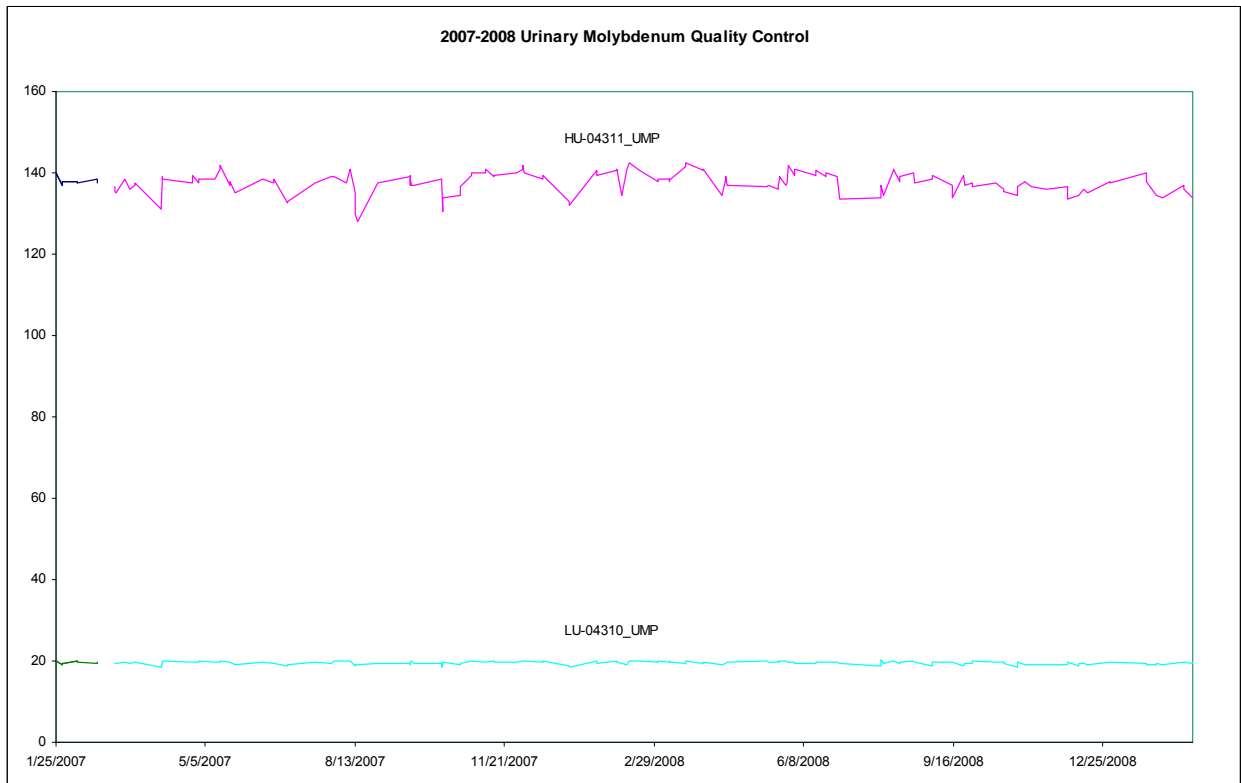

## G. Urinary Lead

### Summary Statistics for Urinary Lead by Lot


Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	151	1/25/2007	2/23/2009	0.4161	0.0269	6.5
HU-04311_UMP	151	1/25/2007	2/23/2009	2.9359	0.0670	2.3



## H. Urinary Platinum

### Summary Statistics for Urinary Platinum by Lot


Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	153	1/25/2007	2/23/2009	0.0983	0.0151	15.4
HU-04311_UMP	153	1/25/2007	2/23/2009	0.8537	0.1036	12.1

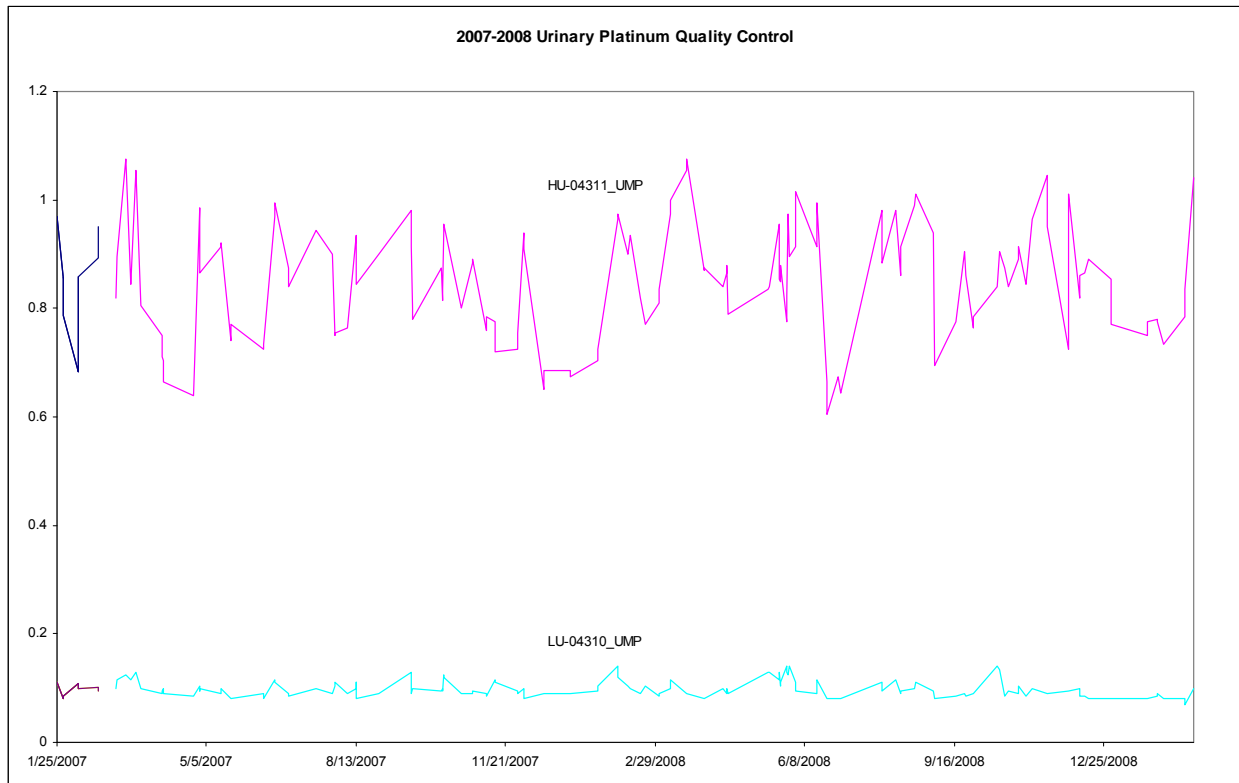

# I. Urinary Antimony

## Summary Statistics for Urinary Antimony by Lot


Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	159	1/25/2007	2/23/2009	0.1943	0.0104	5.4
HU-04311_UMP	159	1/25/2007	2/23/2009	0.6513	0.0266	4.1

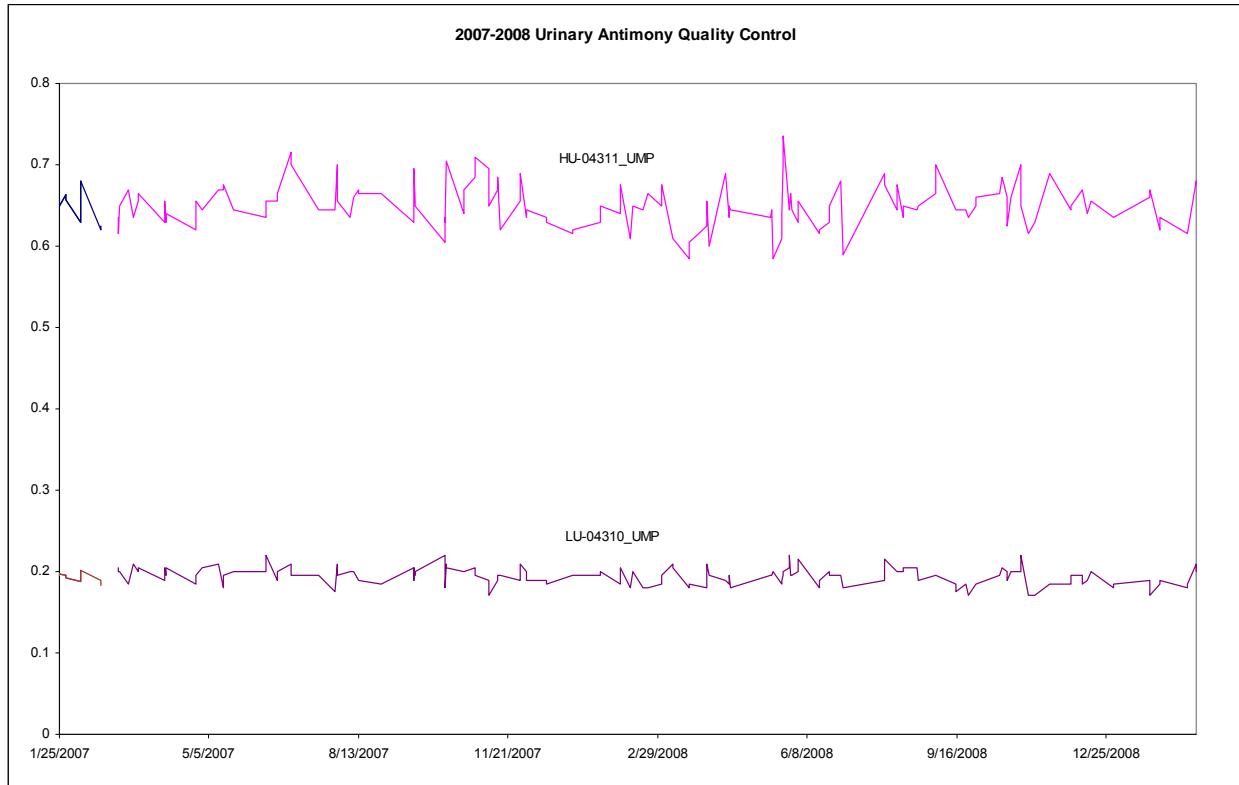

## J. Urinary Thallium

### Summary Statistics for Urinary Thallium by Lot


Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	167	1/25/2007	2/23/2009	0.1786	0.0061	3.4
HU-04311_UMP	167	1/25/2007	2/23/2009	0.5772	0.0134	2.3

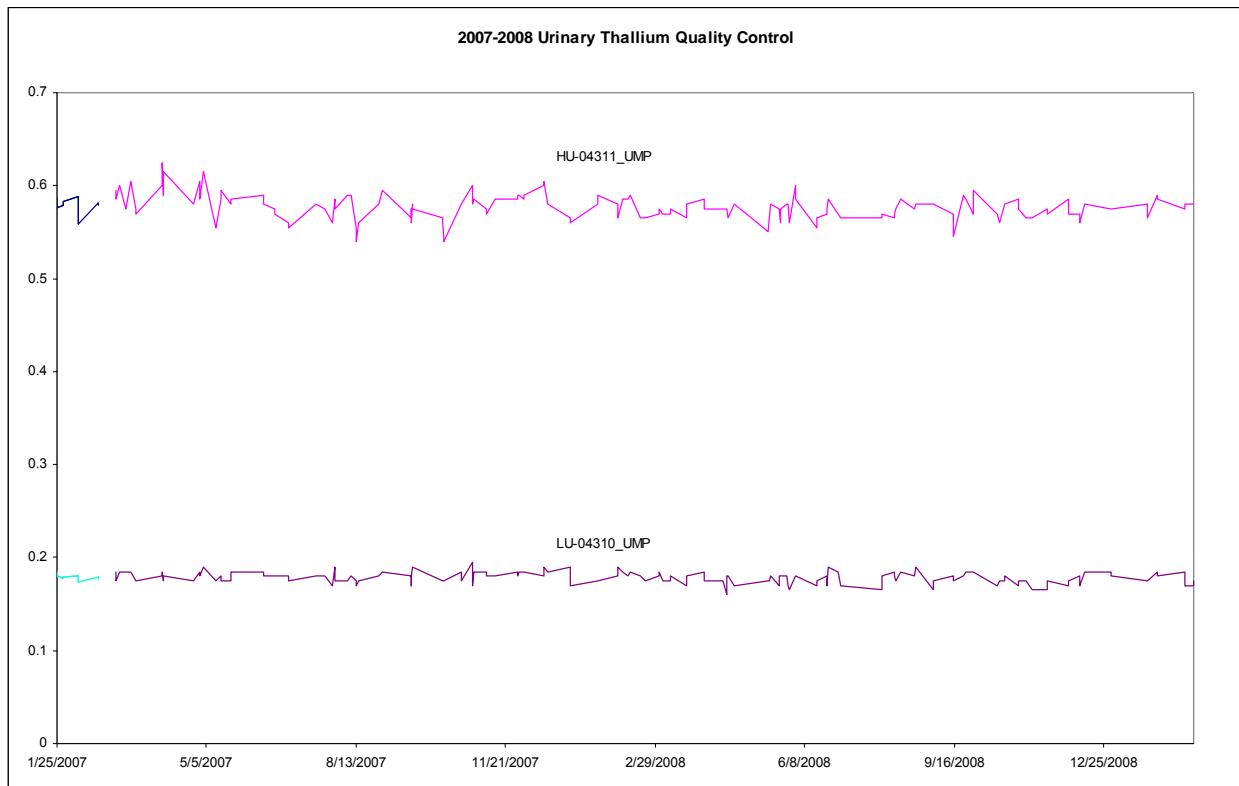

## K. Urinary Tungsten

### Summary Statistics for Urinary Tungsten by Lot

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	158	1/25/2007	2/23/2009	0.2181	0.0129	5.9
HU-04311_UMP	158	1/25/2007	2/23/2009	0.9448	0.0228	2.4

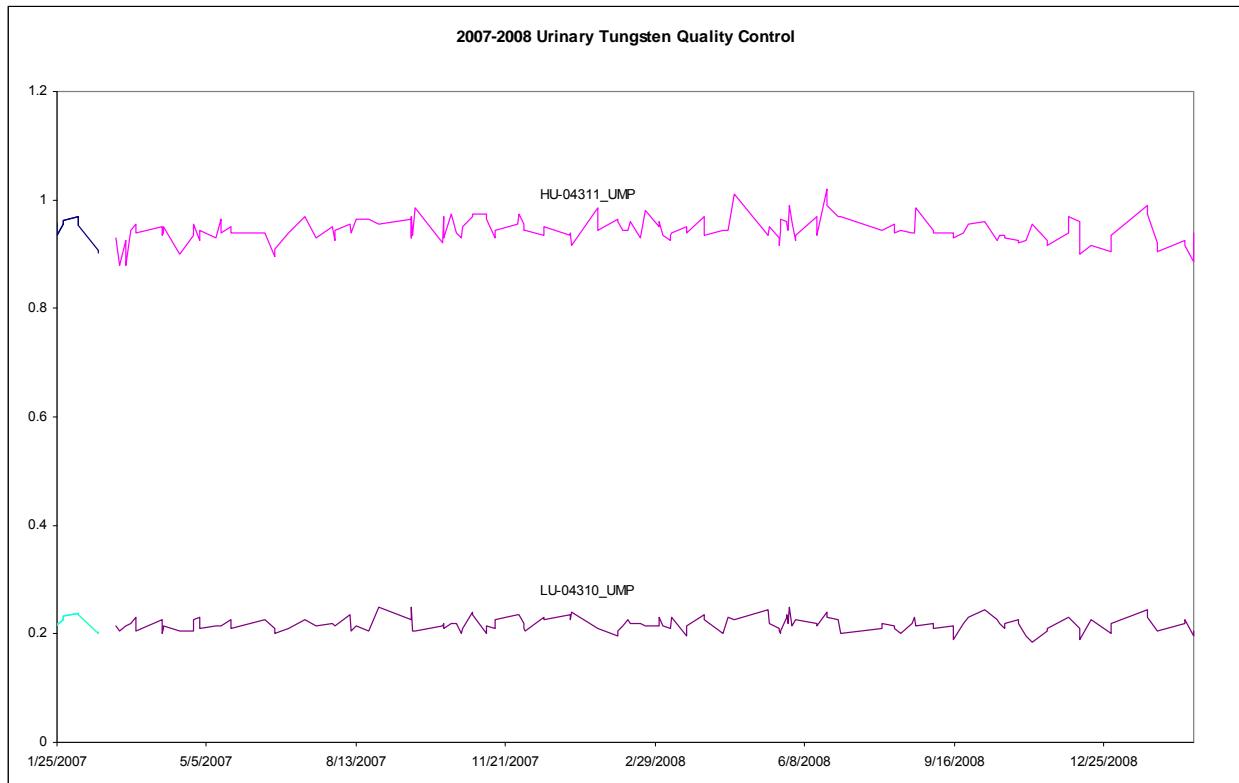




# L. Urinary Uranium

## Summary Statistics for Urinary Uranium by Lot

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
LU-04310_UMP	152	1/25/2007	2/23/2009	0.0133	0.0009	6.9
HU-04311_UMP	152	1/25/2007	2/23/2009	0.1275	0.0066	5.2

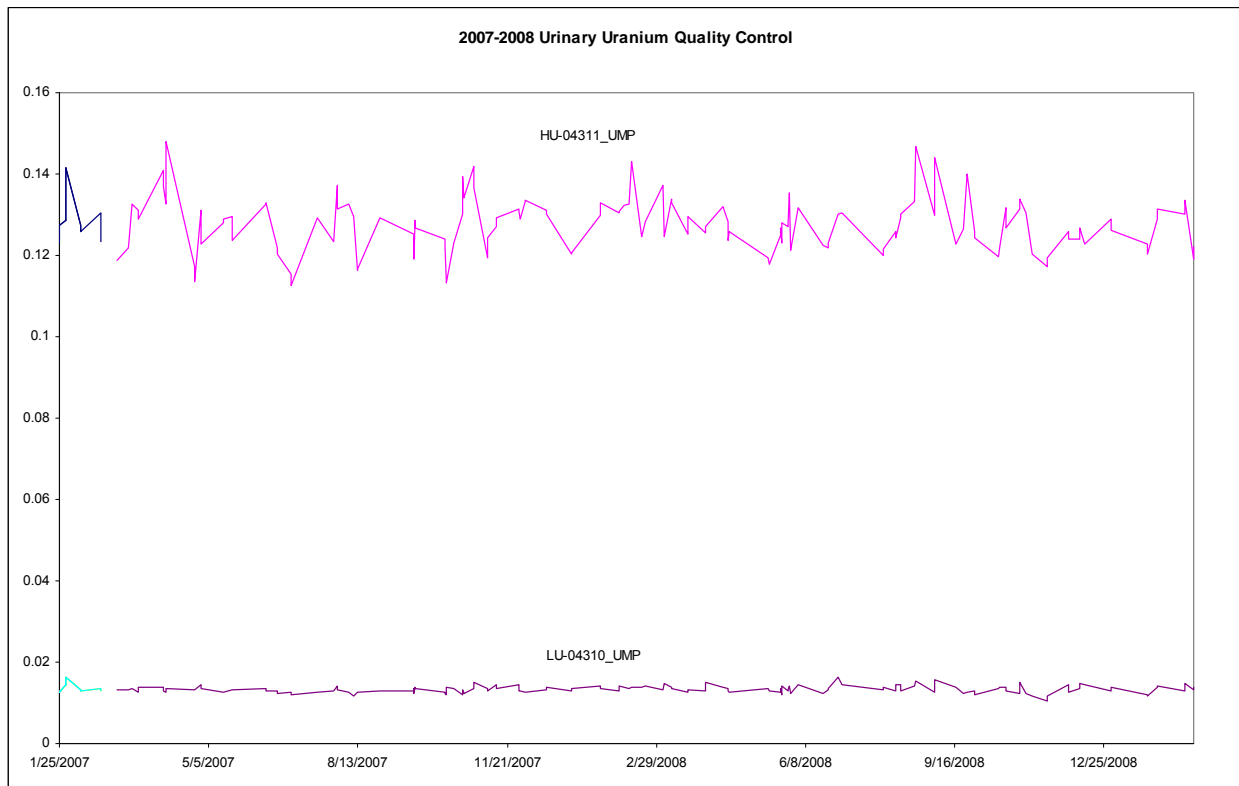

## Appendix

<b>Table 1. Instrument and Method Parameters</b>	
<b>Instrument:</b> PerkinElmer ELAN DRC <sup>Plus</sup> or DRC II ICP-MS CETAC ASX 500 series autosampler (tray B)	
<b>Optimization Window Parameters</b>	
RF power	1.45 KW
Plasma Gas Flow (Ar)	15 L/min
Auxiliary Gas Flow (Ar)	1.2 L/min
Nebulizer Gas Flow (Ar)	~0.90 – 1.0 L/min (optimized as needed for sensitivity)
Ion Lens Voltage(s)	AutoLens (optimized as needed for sensitivity)
QRO, CRO, CPV, Discriminator Threshold	Optimized per instrument by service engineer, or advanced user.
Parameters of x-y alignment, nebulizer gas flow, AutoLens voltages, mass calibration, and detector voltages are optimized regularly. Optimization file name = default.dac.	
<b>Configurations Window Parameters</b>	
Cell Gas Changes Pause Times	Pressurize Delay (From Standard to DRC mode) = 60 Exhaust Delay (From DRC to Standard mode) = 60 Flow Delay (Gas changes while in DRC mode) = 60 Channel Delay (Gas channel change in DRC mode) = 60
<b>File Names &amp; Directories</b>	
Method file names	<i>For calibration curve (programmed for urine blank)</i> 13 elem ur_methITU001B_HPS2107-003_urblk.mth  <i>For QC &amp; patient sample analysis (programmed for aqueous blank)</i> 13 elem ur_methITU001B_HPS2107-003_aqblk.mth
Dataset	Create a new dataset subfolder each day. Name as “2005-0718” for all work done on July 18, 2005
Sample File	Create for each day’s work
Report file name	<i>For sample results printouts</i> cdc_quant comprehensive_13 element.rop  <i>For calibration curve information</i> CDC_Quant Comprehensive (calib curve info).rop
Tuning	Default.tun
Optimization	Default.dac
Calibration	N/A
Polyatomic	elan.ply
Report Options Template (transferring results to the database)	CDC_Database Output.rop <i>Report Format Options: select only “Use Separator”</i> <i>File Write Option: Append</i> <i>Report File name: include date, instrument, and group being analyzed in file name (i.e. 2005-0311b_DRC2A_HM-0364.txt)</i>
<b>Method Parameters</b>	
<b>Method Parameters: Timing Page</b>	

Sweeps/reading	70
Readings/replicate	1
Replicates	3
Enable QC Checking	Off
Isotopes Monitored and Internal Standard Associations (Exact Mass)	<p><u>Group 1</u> (use <math>^{103}\text{Rh}</math> as an internal standard)  <math>^9\text{Be}</math> (9.0122), <math>^{59}\text{Co}</math> (58.9332), <math>^{98}\text{Mo}</math> (97.9055), <math>^{103}\text{Rh}</math> (102.905), <math>^{121}\text{Sb}</math> (120.904), <math>^{133}\text{Cs}</math> (132.905), <math>^{138}\text{Ba}</math> (137.905)</p> <p><u>Group 2</u> (use <math>^{193}\text{Ir}</math> as an internal standard)  <math>^{184}\text{W}</math> (183.951), <math>^{193}\text{Ir}</math> (192.963), <math>^{195}\text{Pt}</math> (194.965), <math>^{205}\text{Tl}</math> (204.975), <math>^{208}\text{Pb}</math> (207.977), <math>^{238}\text{U}</math> (238.05)</p> <p><u>Group 3</u> (use <math>^{103}\text{Rh}</math> as an internal standard)  <math>^{114}\text{Cd}</math> (113.904), <math>^{103}\text{Rh}</math> (102.905)</p> <p><u>Group 4</u> (use <math>^{71}\text{Ga}</math> as an internal standard)  <math>^{71}\text{Ga}</math> (70.9249), <math>^{75}\text{As}</math> (74.9216)</p>
Dwell Times	<p>30 ms for <math>^{59}\text{Co}</math>, <math>^{98}\text{Mo}</math>, <math>^{103}\text{Rh}</math> in Standard mode, <math>^{121}\text{Sb}</math>, <math>^{133}\text{Cs}</math>, <math>^{138}\text{Ba}</math>, <math>^{184}\text{W}</math>, <math>^{193}\text{Ir}</math>, <math>^{205}\text{Tl}</math>, and <math>^{208}\text{Pb}</math></p> <p>50 ms for <math>^{71}\text{Ga}</math> and <math>^{75}\text{As}</math></p> <p>100 ms for <math>^9\text{Be}</math>, <math>^{195}\text{Pt}</math>, <math>^{238}\text{U}</math>, <math>^{103}\text{Rh}</math> in DRC mode, and <math>^{114}\text{Cd}</math></p>
Scan Mode	Peak Hopping for all isotopes (1 MCA channel)
DRC channel A Gas Flow Rate	<p>10% hydrogen / 90% argon (5-7 psig delivery pressure)</p> <p>Typically 0.2 – 0.7 L/min *</p> <p>*(optimized per instrument, and periodically verified)</p>
DRC channel B Gas Flow Rate	<p>Oxygen (5-7 psig delivery pressure)</p> <p>Typically 1.1 – 1.8 L/min *</p> <p>* (optimized instrument, and periodically verified)</p>
RPa	0 for all isotopes
RPq	<p><i>Standard Mode</i>: 0.25 for all standard mode isotopes</p> <p><i>DRC Mode (Cd group)</i>: Typically* 0.65 - 0.75 for <math>^{114}\text{Cd}</math> (113.904), and <math>^{103}\text{Rh}</math> (102.905) in DRC mode. Use the same RPQ for each.</p> <p><i>DRC Mode (As group)</i>: Typically* 0.65 - 0.75 for <math>^{71}\text{Ga}</math> (70.9249), <math>^{75}\text{As}</math> (74.9216). Use the same RPQ for each.</p> <p>(* Optimize per instrument, and periodically verified)</p>
<b>Method Parameters: Processing Page</b>	
Detector mode	Pulse
Process Spectral Peak	N/A
AutoLens	On
Isotope Ratio Mode	Off
Enable Short Settling Time	Off
Blank subtraction	After internal standard
Measurement units	Cps
Process Signal Profile	N/A
<b>Method Parameters: Equations Page</b>	

Equations	On <sup>208</sup> Pb, use "+ Pb 206 + Pb 207" On <sup>238</sup> U, use "+ U 235" On 114Cd, use "- 0.027250 * Sn 118"
<b>Method Parameters: Calibration Page</b>	
Calibration Type	External Std.
Curve type	Simple Linear
Sample units	"µg/L" or "ppb"
Calibration Standard Concentrations (µg/L)	Be: 0.1, 0.3, 1, 3, 10 Co: 0.075, 0.225, 0.75, 2.25, 7.5 Mo: 3, 9, 30, 90, 300 Sb: 0.08, 0.24, 0.8, 2.4, 8 Cs: 0.2, 0.6, 2, 6, 20 Ba: 0.2, 0.6, 2, 6, 20 W: 0.06, 0.18, 0.6, 1.8, 6 Pt: 0.025, 0.075, 0.25, 0.75, 2.5 Tl: 0.04, 0.12, 0.4, 1.2, 4 Pb: 0.1, 0.3, 1, 3, 10 U: 0.005, 0.015, 0.05, 0.15, 0.5 Cd: 0.08, 0.24, 0.8, 2.4, 8 As: 2, 6, 20, 60, 200
<b>Method Parameters: Sampling Page</b>	
"Peristaltic Pump Under Computer Control"	On
Sample Flush	~90s at 10 rpm (optimize time so that solution reaches nebulizer before Read Delay begins)
Read Delay	20s at 10 rpm
Wash	120s at 24 rpm
Autosampler Locations of Blanks and Standards	<i>For calibration curve (points to urine blank)</i> cdc_13 elem ur_methITU001B_HPS2107-003_urblk.mth Urine Blank and Calibration Stds 1 – 5 in autosampler positions 11 – 16.  <i>For QC &amp; patient sample analysis (points to aqueous blank)</i> cdc_13 elem ur_methITU001B_HPS2107-003_aqblk.mth Aqueous Blank in autosampler position 19.

Appendix (continued)

<b>Table 2. Suggested maximum analyte concentrations for base urine.</b>	
<b>Analyte</b>	<b>Concentration (µg/L)</b>
Be	0.5
Co	0.25
Mo	30
Sb	0.2
Cs	3
Ba	2
W	0.2
Pt	0.25
Tl	0.2
Pb	0.75
U	0.03
Cd	0.25
As	5

<b>Table 3. Concentrations of Analytes in the Multi-Element Intermediate Stock Standard from High Purity Standards.</b>	
<b>Analyte</b>	<b><i>Intermediate Stock Standard Concentrations (mg/L)</i></b> High Purity Standards Item # SM-2107-003
Be	10
Co	7.5
Mo	300
Sb	8
Cs	20
Ba	20
W	6
Pt	2.5
Tl	4
Pb	10
U	0.5
As	200
Cd	8

Appendix (continued)

<b>Table 4. Preparation of Multi-element Intermediate Working Standards</b>					
<b>Standard #</b>	<b>1</b>	<b>2</b>	<b>3</b>	<b>4</b>	<b>5</b>
<b>Volume of Flask (mL)</b>	500	200	100	100	100
<b>Volume Spike of Int. Stock Std. (mL)</b>	0.050	0.060	0.100	0.300	1.00
	<b>Concentrations (ug/L) †</b>				
Be*	1	3	10	30	100
Co*	0.75	2.25	7.5	22.5	75
Mo*	30	90	300	900	3000
Sb*	0.8	2.4	8	24	80
Cs*	2	6	20	60	200
Ba*	2	6	20	60	200
W†	0.6	1.8	6	18	60
Pt†	0.25	0.75	2.5	7.5	25
Tl†	0.4	1.2	4	12	40
Pb†	1	3	10	30	100
U†	0.05	0.15	0.5	1.5	5
Cd*	0.8	2.4	8	24	80
As‡	20	60	200	600	2000
* Rh-103 used as internal standard † Ir-193 used as internal standard ‡ Ga-71 used as internal standard ‡ A further 1:10 dilution occurs when added to base urine. Enter the table 4 concentrations divided by 10 into the ELAN software (method window, calibration page).					

Appendix (continued)

<b>Table 5. Acceptable ways to perform two consecutive analytical runs, bracketing with bench quality control samples.</b>	
<b>Setup 1</b>	<b>Setup 2 (typical)</b>
<p><i>Run #1</i>            Calibration Standards            Low Bench QC            High Bench QC            patient samples            Low Bench QC            High Bench QC</p> <p><i>Run #2</i>            Low Bench QC            High Bench QC            patient samples            Low Bench QC            High Bench QC</p>	<p><i>Run #1</i>            Calibration Standards            Low Bench QC            High Bench QC            patient samples            Low Bench QC            High Bench QC</p> <p><i>Run #2</i>            Calibration Standards            Low Bench QC            High Bench QC            patient samples            Low Bench QC            High Bench QC</p>

Appendix (continued)

<b>Table 6. A typical SAMPLE/BATCH window.</b>			
<u>AS Location*</u>	<u>Sample ID</u>	<u>Measurements Action</u>	<u>Method</u>
219	DRCstability1	Run sample	.. .003_urblk.mth
219	DRCstability2	Run sample	.. .003_urblk.mth
219	DRCstability3	Run sample	.. .003_urblk.mth
219	DRCstability4	Run sample	.. .003_urblk.mth
Continue DRC stability samples . . .			
219	DRCstability9	Run sample	.. .003_urblk.mth
219	DRCstability10	Run sample	.. .003_urblk.mth
100	UrBlkChk1	Run blank, standards, and sample **	.. .003_urblk.mth
101	UrBlkChk2	Run sample	.. .003_urblk.mth
20	Aq Blk Check	Run blank and sample †	.. .003_aqblk.mth
102	L Bench QC	Run sample	.. .003_aqblk.mth
103	H Bench QC	Run sample	.. .003_aqblk.mth
21	Sample 1	Run sample	.. .003_aqblk.mth
22	Sample 2	Run sample	.. .003_aqblk.mth
23	Sample 3	Run sample	.. .003_aqblk.mth
104	L Bench QC	Run sample	.. .003_aqblk.mth
105	H Bench QC	Run sample	.. .003_aqblk.mth
<p>* The exact autosampler positions of QCs and patient samples do not have to be those shown above, but the order in which these are run should be as shown above.</p> <p>** When executing this row, the ELAN will first analyze the urine blank at AS position 11, then standards 1-5 at autosampler positions 12-16, <u>then</u> the “UrBlkChk1” sample at A/S position 100. The sampling information about AS positions 11-16 are stored in the “urblk” method file.</p> <p>† When executing this row, the ELAN will first analyze the aqueous blank at AS position 19, then the “Aq Blk Check” at AS position 20. The sampling information about AS positions 19 is stored in the “urblk” method file.</p>			



Appendix (continued)

<b>Table 7. Preparation of Multi-element Intermediate Working Standards</b>					
<b>Dilution ID</b>	<b>Water (μL)</b>	<b>Base Urine (μL)</b>	<b>AQ Intermediate Working Standard (μL)</b>	<b>Patient or QC urine sample (μL)</b>	<b>Diluent (μL)</b>
AQ Blank	1000	-	-	-	9000 *
Urine Blank and UrBlkChk	100	900	-	-	9000 *
Working Calibration Standards	-	900	100	-	9000 *
Patient Urine or Urine-Based QC	-	-	-	500	4500
Patient Urine <i>2x Dilution</i> <sup>†</sup>	500	-	-	500	9000 *
Patient Urine <i>10x Dilution</i> <sup>†</sup>	900	-	-	100	9000*
* 9000 μL diluent is best dispensed from the Digiflex™ as 2 4500-μL portions (i.e.- When preparing a Working Calibration Standard dilution, dispense 4500 μL diluent + 100 μL water in one cycle of Digiflex™, then 4500 μL diluent + 900 μL base urine in the next cycle of the Digiflex™ to prepare a 10 mL total volume dilution.)					
<sup>†</sup> Extra dilution is performed on urine samples whose concentration is greater than the concentrations listed in Table 8 in the Appendix (linearity of the method has been documented up to these concentrations). Any extra level of dilution can be prepared as long as the 9:10 ratio of diluent to total dilution volume is maintained. Use of the lowest possible dilution level is preferred because matrix differences may lead to different observed concentration results as the sample dilution becomes greater (i.e. 2x dilution) is preferred over 10x if 2x is sufficient to dilute analyte into the documented linearity range).					

Appendix (continued)

<b>Analyte</b>	<b>Highest Conc. (<math>\mu\text{g/L}</math>) Verified in Calibration Verification (“Range of Linearity Tested”, or “RLT”) *</b>
Be	100
Co	150
Mo	1800
Sb	100
Cs	1000
Ba	300
W	100
Pt	50
Tl	25
Pb	250
U	40
Cd	100
As	3000

\* If observed results are not within 10% of target, investigate the problem with the involvement of the lab supervisor.

Appendix (continued)

<b>Table 9. Boundary Concentrations for Urine Concentrations (µ/L).</b>			
<b>Analyte</b>	<b>1<sup>st</sup> Upper Boundary (“1UB”) *</b>	<b>2<sup>nd</sup> Upper Boundary (“2UB”) **</b>	<b>Range Maximum (“Lim Rep Delta”) †</b>
Be	0.2	0.4	0.3
Co	2.83	5.66	0.3
Mo	293.5	587	4.0
Sb	0.8	1.6	0.2
Cs	16.5	33	0.5
Ba	17.1	34.2	0.4
W	1.38	2.76	0.2
Pt	0.2	0.4	0.2
Tl	0.62	1.24	0.2
Pb	7.8	15.6	0.3
U	0.277	0.554	0.03
Cd	2.54	5.08	0.3
As	100	200	10

\* Typically, the 1<sup>st</sup> upper boundary (1UB) is the 99<sup>th</sup> percentile of non-weighted, non-creatinine corrected concentration results from the NHANES 1999-2000 subset groups. Concentrations observed greater than the “first upper boundary” (defined in the laboratory database as the “1UB”) should be confirmed by repeat analysis of a new sample preparation. The concentration assigned to the 1UB for an element is determined by study protocol but default concentrations are listed in this table. Report the original result, as long as the confirmation is within 10% of the original. Continue repeat analysis until a concentration can be confirmed.

\*\* Typically the 2<sup>nd</sup> upper boundary (2UB) is set to 2x the 1UB. At the discretion of the supervisor, the 1UB may vary per study according to the concerns of the study. Regardless of the study, report patient results confirmed to be greater than the 2UB to the QC reviewer as an “elevated result”.

† Range maximum is the range of the three replicate readings for a single sample analysis. This value is also called the “Lim RepDelta” in the database which handles data for the Inorganic Toxicology and Nutrition Branch. If the range of replicate readings is greater than the range maximum, and represents greater than a 10% relative standard deviation for the measurement, do not use the measurement for reporting.

Appendix (continued)

<b>Table 10. Reference Ranges for Urine Concentrations (from the Third National Report on Exposure to Environmental Chemicals [12] ). All results in <math>\mu\text{g/L}</math>.</b>							
<b>Analyte</b>	<b>Survey Years</b>	<b>Geometric Mean</b>	<b>50<sup>th</sup></b>	<b>75<sup>th</sup></b>	<b>90<sup>th</sup></b>	<b>95<sup>th</sup></b>	<b>N</b>
Be	99-00	$\leq 0.13$ *	$\leq 0.13$	$\leq 0.13$	$\leq 0.13$	$\leq 0.13$	2465
	01-02	$\leq 0.13$	$\leq 0.13$	$\leq 0.13$	$\leq 0.13$	$\leq 0.13$	2690
Co	99-00	0.375	0.400	0.630	0.940	1.32	2465
	01-02	0.379	0.410	0.610	0.930	1.27	2690
Mo	99-00	45.9	50.7	84.9	134	178	2257
	01-02	45.0	52.4	83.3	124	165	2690
Sb	99-00	0.132	0.130	0.210	0.330	0.420	2276
	01-02	0.134	0.130	0.180	0.260	0.340	2690
Cs	99-00	4.35	4.80	7.10	9.60	11.4	2464
	01-02	4.81	5.49	7.91	10.4	12.6	2690
Ba	99-00	1.50	1.50	3.00	5.40	6.80	2180
	01-02	1.52	1.63	3.12	5.22	7.48	2690
W	99-00	0.093	0.090	0.180	0.320	0.500	2338
	01-02	0.082	0.060	0.150	0.300	0.450	2652
Pt	99-00	$\leq 0.04$ **	$\leq 0.04$	$\leq 0.04$	$\leq 0.04$	$\leq 0.04$	2465
	01-02	$\leq 0.04$	$\leq 0.04$	$\leq 0.04$	$\leq 0.04$	$\leq 0.04$	2690
Tl	99-00	0.176	0.200	0.280	0.400	0.450	2413
	01-02	0.165	0.180	0.270	0.360	0.440	2653
Pb	99-00	0.766	0.800	1.30	2.10	2.90	2465
	01-02	0.677	0.600	1.20	2.00	2.60	2690
U	99-00	0.008	0.007	0.013	0.026	0.046	2464
	01-02	0.009	0.008	0.014	0.029	0.046	2690
Cd	99-00	0.193	0.232	0.475	0.858	1.20	2257
	01-02	0.210	0.229	0.458	0.839	1.20	2690
As †	See Table 11 p.60 for urine arsenic reference values.						
* Results were lower than the method detection limit of 0.13 $\mu\text{g/L}$ .							
** Results were lower than the method detection limit of 0.04 $\mu\text{g/L}$ .							
† Urine As was not included in the Third National Report on Exposure to Environmental Chemicals.							

Appendix (continued)

<b>Table 11. References to Total Urine Arsenic Concentrations</b>		
<u>Reference</u>	<u>Concentration (µg/L)</u>	<u>Group Type Sampled</u>
<b>Normal</b>		
Stokinger, 1981 [13]	<100	
Fowler, 1977 [14]	15	
Iffland, 1994 [15]	1 – 80 (Generally < 10)	
<b>Elevated</b>		
Iffland, 1994 [15]	300	After seafood consumption
	200	Copper smelter workers
	5 - 952 vs. 5 - 365	Wood treatment workers vs. comparison group (11)
	25.9 – 667	Seafood-preferring population
	74.1 378.1	Low-Inhalation exposure High-Inhalation exposure (12)
Gerhardsson et al., 1996 [16]	50 – 100	High intake of seafood or increased exposure of inorganic arsenic from food or air

## References

1. Thomas, R., *Practical Guide to ICP-MS (Practical Spectroscopy)*. 2003, New York, NY: Marcel Dekker. 336.
2. Tanner, S.D., Baranov, Vladimir I, *Theory, Design, and Operation of a Dynamic Reaction Cell for ICP-MS*. Atomic Spectroscopy, 1999. **20**(2): p. 45-52.
3. Tanner, S.D., V.I. Baranov, and D.R. Bandura, *Reaction cells and collision cells for ICP-MS: a tutorial review*. Spectrochimica Acta Part B-Atomic Spectroscopy, 2002. **57**(9): p. 1361-1452.
4. PerkinElmer SCIEX Instruments, *ELAN DRC II Hardware Guide*. 2001, Canada.
5. Mulligan, K.J., T.M. Davidson, and J.A. Caruso, *Feasibility Of The Direct Analysis Of Urine By Inductively Coupled Argon Plasma Mass-Spectrometry For Biological Monitoring Of Exposure To Metals*. Journal Of Analytical Atomic Spectrometry, 1990. **5**(4): p. 301-306.
6. Jarrett, J.M., *Total Urine Arsenic Biomonitoring Using Inductively Coupled Plasma Mass Spectrometry with a Dynamic Reaction Cell*. 2005, Centers for Disease Control and Prevention.
7. Jarrett, J.M., *Elimination of Molybdenum Oxide Interference In Urine Cadmium Analysis Using Inductively Coupled Plasma Reaction Cell Mass Spectrometry*. 2004, Centers for Disease Control and Prevention.
8. Office of Health and Safety in the Division of Laboratory Sciences, *Policies and Procedures Manual*. 2002, Division of Laboratory Sciences (DLS), National Center for Environmental Health, Centers for Disease Control and Prevention, Public Health Service, Department of Health and Human ServicesCenters for Disease Control and Prevention.
9. Centers for Disease Control and Prevention (CDC) Radiation Safety Committee, *CDC/ATSDR Occupational Health and Safety Manual (Radiation Safety chapter)*. Centers for Disease Control and Prevention, Public Health Service, Department of Health and Human ServicesCenters for Disease Control and Prevention.
10. ASTM G 94 – 92 (1998), "Standard Guide for Evaluating Metals for Oxygen Service," *ASTM International*, For referenced ASTM standards, visit the ASTM website, [www.astm.org](http://www.astm.org), or contact ASTM Customer Service at [service@astm.org](mailto:service@astm.org). For Annual Book of ASTM Standards volume information, refer to the standard's Document Summary page on the ASTM website.
11. U.S. Nuclear Regulatory Commission, *Regulatory guide 8.22 (revision 1). Bioassay at uranium mills*. 1988: Atlanta, GA.
12. Centers for Disease Control and Prevention, *Third National Report on Human Exposure to Environmental Chemicals*, <http://www.cdc.gov/exposurereport>. 2005.
13. Stokinger, H.E., *The metals*, in *Patty's industrial hygiene and toxicology*, G. Clayton and F. Clayton, Editors. 1981, John Wiley and Sons: New York. p. 1493-2060.
14. Fowler, B.A., in *Toxicology of trace elements*, R. Goyer and M. Mehlman, Editors. 1977, John Wiley and Sons: New York. p. p. 79.
15. Iffland, R., *Arsenic*, in *Handbook on metals in clinical and analytical chemistry*, H. Seiler, A. Sigel, and H. Sigel, Editors. 1994, Marcel Dekker, Inc.: New York. p. 238-250.
16. Gerhardsson, L. and S. Skerfving, *Concepts on biological markers and biomonitoring for metal toxicity*, in *Toxicology of metals*, L. Chang, Editor. 1996, CRC Press: Boca Raton, Florida. p. 98.

17. Amarasiriwardena, C.J., et al., *Determination of the total arsenic concentration in human urine by inductively coupled plasma mass spectrometry: a comparison of the accuracy of three analytical methods*. *Analyst*, 1998. **123**(3): p. 441-445.
18. Larsen, E.H. and S. Sturup, *Carbon-enhanced Inductively Coupled Plasma Mass Spectrometric Detection of Arsenic and Selenium and Its Application to Arsenic Speciation*. *Journal Of Analytical Atomic Spectrometry*, 1994. **9**: p. 1101-1105.