

# Laboratory Procedure Manual

*Analyte:* **Thyroglobulin Antibodies**

*Matrix:* **Serum**

*Method:* **Access 2 (Beckman Coulter)**

*Method No:*

*Revised:*

*as performed by:*

University of Washington Medical Center  
Department of Laboratory Medicine  
Immunology Division  
Director: Mark Wener M.D.  
Supervisor: Kathleen Hutchinson M.S., M.T. (ASCP)  
Authors: Michael Walsh, MT (ASCP), September 2006

*contact:*

Dr Mark Wener, M.D.

## *Important Information for Users*

University of Washington periodically refines these laboratory methods. It is the responsibility of the user to contact the person listed on the title page of each write-up before using the analytical method to find out whether any changes have been made and what revisions, if any, have been incorporated.

### Public Release Data Set Information

This document details the Lab Protocol for testing the items listed in the following table:

<b>File Name</b>	<b>Variable Name</b>	<b>SAS Label</b>
THYROD_F	LBXATG	Thyroglobulin Antibodies (IU/mL)

## 1. SUMMARY OF TEST PRINCIPLE AND CLINICAL RELEVANCE

The Access thyroglobulin antibody assay is a sequential two-step immunoenzymatic "sandwich" assay. A sample is added to a reaction vessel with paramagnetic particles coated with the thyroglobulin protein. After incubation, materials bound to the solid phase are held in a magnetic field while unbound materials are washed away. The thyroglobulin-alkaline phosphatase conjugate is added and binds to the TgAb. After a second incubation, the reaction vessel is washed to remove unbound materials. A chemiluminescent substrate, Lumi-Phos\*\* 530 is added to the reaction vessel and light generated by the reaction is measured with a luminometer. The light production is proportional to the concentration of Thyroglobulin antibody in the sample. The amount of Thyroglobulin antibody in the sample is determined by means of a stored, multi-point calibration curve.

Thyroglobulin is the precursor of the thyroid hormones tri-iodothyronine (T3) and thyroxine (T4). The presence of anti-Tg antibodies in patients has been reported in patients with Hashimoto's, Grave's disease (chronic primary hyperthyroidism) as well as other variants of chronic primary hypothyroidism such as myxedema and asymptomatic thyroiditis. The role of Tg antibodies is not clear. They do not appear to be pathogenic, but may be simply indicators of disease. Coexisting antibodies to TPO are common.

Autoantibodies to the two major thyroid antigens (Tg and TPO) are sometimes used as exclusionary tests in patients with thyroid disease because over 95% of patients with thyroiditis have autoantibodies to either or both of these antigens. However, a positive result cannot be used to make the diagnosis, since the antibodies can be found in a variety of other conditions. Although the presence of high titers of antibodies to thyroglobulin is often indicative of thyroid autoimmune disease, the data must be considered in light of other clinical and laboratory findings. Some individuals may have high levels of Tg antibodies with little or no evidence of clinical disease. By contrast a few clinically defined patients may have undetectable levels of this antibody. Low titers of Tg antibody may be found in apparently healthy individuals. The clinical significance of these results is undetermined.

This test is also needed to interpret the results of the tumor marker, thyroglobulin. Thyroglobulin antibodies often interfere with the measurement of the thyroglobulin. All samples ordered for thyroglobulin are screened for the presence of antibodies.

## 2. SAFETY PRECAUTIONS

Consider all samples received for analysis potentially positive for infectious agents including HIV and the hepatitis B virus. Observe universal precautions. Wear gloves, lab coat, and safety glasses when handling all human blood products and infectious viruses. Place disposable plastic, glass, paper, and gloves that contact blood in a biohazard bag or discard pan to be autoclaved. Disinfect all work surfaces with Vesphene solution. Dispose of all biological samples and diluted specimens in a biohazard bag at the end of the analysis.

Do not pipette by mouth. Do not eat, drink or smoke in designated work areas. Wash hands thoroughly after removal or personal protective devices used in handling specimens and kit reagents.

Material safety data sheets for all reagents used in the performance of this assay are kept

in the Immunology Division, University of Washington Medical Center (UWMC).

### 3. COMPUTERIZATION; DATA SYSTEM MANAGEMENT

- A. Each shipment of specimens received from the NHANES mobile unit arrives with a corresponding transmittal sheet and an electronic version of the shipping/resulting file. The file structure is determined by NHANES and is described in the National Health and Nutrition Examination Survey (NHANES) Contract Laboratory Manual.
- B. After the testing is completed results from the Access 2 are transferred to the laboratory server system, which is backed up daily. This instrument file contains the following information for each sample, control and calibrator tested.

- Patient ID
- Sample ID
- Rack
- Verify
- Test Name
- Interpretation
- Result
- Units
- Comp. Time
- Flags
- LIS
- Instrument
- RLU
- Pipettor
- Sample Type
- Sample Priority
- Test ID
- Reagent Pack Lot #
- Reagent Pack Serial #
- Dilution
- Calibrator level
- Comments
- Load Date/Time

- C. QC results are transferred to an Excel file using laboratory-developed software. This file calculates the QC statistics, plots Levey-Jennings charts, displays relevant instrument flags, tracks reagent lots and recent calibrations. QC results are reviewed prior to resulting samples.
- D. Sample results are transferred to an Excel file using laboratory-developed software that enters results after matching sample identifiers from the instrument file with those provided in the NHANES shipping/resulting file. This Excel file is formatted to match the NHANES shipping/resulting file and the program uses the conventions outlined in the NHANES Contract Laboratory Manual.
- E. Data entry is checked for errors.

- F. After the TGAB testing has also been completed, resulted, and checked, the result file is transmitted electronically to NHANES WESTAT. Electronic and hard copies of the files are kept in the laboratory.
- G. Technical support for this system is provided by Westat, Rockville, MD (1-301-294-2036)

**4. SPECIMEN COLLECTION, STORAGE, AND HANDLING PROCEDURES; CRITERIA FOR SPECIMEN REJECTION**

- A. No special instructions such as fasting or special diets are required.
- B. Serum is the preferred specimen type. Heparin plasma is acceptable. If testing is to be done within 48 hours, samples can be refrigerated at 2 to 8°C. Freeze at -20°C or colder for longer storage.
- C. Blood should be collected aseptically and the serum separated by standard laboratory techniques. Specimens may be collected by using regular or serum-separator Vacutainers. Serum should be separated from the cells within 2 hours of collection.
- D. The requested sample volume for the assay is 1.0 mL, and the minimum sample volume is 0.3 mL.
- E. Specimens may be stored in glass or plastic vials, as long as the vials are tightly sealed to prevent desiccation of the sample.
- F. Turbid samples or those with particulate matter should be centrifuged prior to assay.
- G. More than two freeze-thaw cycles is not recommended.

**5. PROCEDURES FOR MICROSCOPIC EXAMINATIONS; CRITERIA FOR REJECTION OF INADEQUATELY PREPARED SLIDES**

Not applicable for this procedure.

**6. PREPARATION OF REAGENTS, CALIBRATORS (STANDARDS), CONTROLS, AND ALL OTHER MATERIALS; EQUIPMENT AND INSTRUMENTATION**

A. Instrumentation

- 1. Beckman Access or Access II Immunoassay System (Beckman Coulter, Fullerton, CA.)

The Beckman Access is a fully automated, random access, instrument that features on-board storage of reagent packs in a refrigerated compartment; an ultrasonic probe tip for level sense detection, sample and reagent delivery, mixing, and probe cleaning to minimize carryover; barcode identification of specimens and reagent packs; temperature controlled reaction reactions; and

measurement and analysis of the light signal generated by the chemiluminescent reaction (RLU) using a smoothing spline math model.

The Thyroglobulin Antibody II assay parameter settings for the instrument are as follows:

Parameter	Setting
Sample Volume Requirements	
Minimum sample volume	175 ul
Sample volume used for testing	10 ul
No. of Standard Points	6
Calibration curve calculation	smoothing spline math model
Standard Curve Measuring Range (At initial dilution; approximate values, range is dependent upon standard value)	0 – 2500 IU/mL

2. Hewlett Packard DeskJet printer (Hewlett Packard, Boise, ID)
3. Computers (Dell Computer Corporation, Round Rock, Texas).
4. Centrifuge (Jouan Inc., Winchester, VA)

B. Equipment

1. Reaction Vessels (Beckman Coulter, Fullerton, CA)
2. Sample Cups (Fisher Scientific, Pittsburgh, PA)
3. Latex gloves, disposable (Any manufacturer).
4. Pipettes and tips (Rainin, Emeryville, CA)

C. Reagents

All reagents are purchased from Beckman Coulter, Fullerton, CA.

1. R1: Access Thyroglobulin Antibody reagent packs:  
Cat. No. A32898: 100 determinations, 50 tests/pack

Provided ready to use. Store upright and refrigerate packs at 2 to 10°C. Packs must be refrigerated at 2 to 10°C for a minimum of two hours before use on the instrument. Stable until the expiration date stated on the label when stored at 2 to 10°C. After initial use, the pack is stable at 2 to 10°C for 56 days. Signs of possible deterioration are a broken elastomeric layer on the pack or control values out of range.

R1a: Dynabeads Paramagnetic particles coated with streptavidin coupled with biotinylated human thyroglobulin, suspended in a Tris buffer with protein (bovine) <0.1% sodium azide, and 0.1% ProClin\*\* 300.

R1b: Human Thyroglobulin -alkaline phosphatase (bovine) conjugate in a Tris buffer, with protein (bovine), <0.1% sodium azide 0.1% ProClin\*\* 300.

R1c: Tris buffer with <0.1% sodium azide and 0.1% ProClin 300

2. Access Substrate      Cat. No. 81906: 4 x 130 ml

Lumi-Phos\*530 (buffered solution containing dioxetane), Lumigen\* PPD, fluorescer, and surfactant. Bring to room temperature (15 – 30 °C) at least 18 hours before use. Stable for 14 days at room temperature or after bottle has been opened.

3. Access Wash Buffer II: Cat # A16792

Provided ready to use. Store at room temperature (15 – 30 °C), stable until expiration date on label.

4. Access Thyroglobulin antibody Sample Diluent: Cat # 33866

HEPES buffer with bovine serum albumin (BSA), < 0.1% sodium azide, and 0.5% ProClin 300. Contains 0.0 IU/mL Thyroglobulin antibody. Provided ready to use. Allow the contents to stand for 10 minutes at room temperature and mix by gently inverting prior to use. Avoid bubble formation. Stable until the expiration date when stored at 2-8 °C.

D. Standards/Calibration Preparation

Access Thyroglobulin antibody Calibrators Cat. No. A36920: 2.5 ml/vial

Provided ready to use. Store at 2 to 10°C. Mix contents by gently inverting before use. Avoid bubble formation. Stable until the expiration date stated on the vial labels when stored at 2 to 10°C. Control values out of range are a sign of possible deterioration.

Human thyroglobulin antibody in human serum at levels of approx. 50, 250, 500, 1000 and 2500 IU/mL. Zero calibrator is human serum with <0.1% sodium azide, 0.5% ProClin\*\* 300. S0 Std contains 0.0 IU/mL Thyroglobulin antibody. Refer to calibration card for exact concentrations. Calibration Card (with actual calibrator concentration information)

E. Preparation of Quality Control Materials

Three different levels of serum controls are run with each run. The controls are purchased from BioRad Laboratories (Hercules, CA) or prepared in-house. Commercial controls are stored and used according to the manufacturer's recommendations. In house controls are stored frozen (-20°C or colder). Once thawed, the controls are stored at 2-8 °C. All controls are used within their stated expiration dates.

## 7. CALIBRATION AND CALIBRATION VERIFICATION PROCEDURES

TgAB concentrations are calculated by using a calibration curve. This method utilizes a smoothing spline standard curve with a direct relationship of measured light produced (RLU) to concentration of TgAB in the serum sample. Serum results are expressed as IU/ mL.

Calibrators are standardized to the WHO 65/93 international standard.

An active calibration curve is required for all tests. For the Access Thyroglobulin Antibody II assay, calibration is required every 56 days or whenever new lot numbers of reagents are placed into use. Refer to the Operator's Guide and Reference Manual for complete instructions on calibration procedures.

### B. Verification

1. Three levels of control are run for each test series. If, within a testing series, these controls do not conform to specifications as defined in the quality control manual, the entire series is invalidated.
2. New lot numbers of calibrator are verified by running 100 or more samples tested on the previous lot number. The correlation is analyzed using one or more linear regression formulas.

## 8. PROCEDURE OPERATING INSTRUCTIONS; CALCULATIONS; INTERPRETATION OF RESULTS

### A. Preliminaries

1. Bring all controls and patient specimens to room temperature before use. Mix any specimens or controls that have been frozen. Centrifuge samples with particulate matter prior to testing.
2. Prime system: pipettor - 1 time, and substrate - 4 times
3. Check reagent, substrate, wash buffer, and reaction vessel status. Load any needed supplies onto the instrument. Mix reagent pack contents by gently inverting pack several times before loading on the instrument. Do not invert open (punctured) packs — mix reagents by swirling gently.

### B. Instrument Operation (see operator's manual for details).

1. Check sample volume to make sure that there is sufficient volume to perform testing. Gently mix, uncap and load specimens into specimen racks, with the barcode in the open slot. Make sure there are no bubbles. Alternately, use the barcode wand to identify the specimens. If the barcode is not reading properly, sample IDs can be entered manually. Load the racks onto the instrument.
2. Select the ThgAB test. Note: if other thyroid testing is also ordered, the entire 8 test panel can be ordered as a group. Testing is done in singlicate. Select the sample(s) to be used for the random repeat testing.



3. The instrument automatically calculates all results. After testing is completed, results are printed and review by the technologist.
4. Remove specimens and controls soon after the instrument finishes pipetting from the sample. Return controls to the refrigerator and refreeze specimens.
5. Perform scheduled instrument maintenance (daily, weekly, and monthly) as outlined on the maintenance log. See the operator's manual for specific instructions.

C. Recording of Data

1. Using a lab developed program, specimen results are transferred from the instrument data file into the assay specific results table created from the send file corresponding to the specific sample box. The file format is Excel (Microsoft Corporation, Redmond WA). A copy of this file is printed out and checked for accuracy of data entry.
2. Control results are entered to the Assay Specific QC/Levy-Jennings Table using the Excel program. Compliance with the Westgard rules is evaluated. A copy of this table is printed out and checked for accuracy of data entry.

D. Replacement and Periodic Maintenance of Key Components

1. Daily Maintenance:

**Start-up:**

Inspect fluidics module.  
Check system supplies and replace as needed.  
Clean exterior of substrate, dispense, and aspirate probes.  
Prime pipettor – 1X and substrate - 4X.  
Verify temperature.

**Shut-down:**

Check waste containers, empty if needed  
Perform clean

2. Weekly Maintenance:

Change probes and clean them  
Clean exterior of the analyzer  
Clean upper portion of the main pipettor with alcohol wipe  
Inspect waste filter bottle for fluid  
Run system check

3. Periodic Maintenance to be performed by the manufacturer's service engineer.

E. Calculations

Patient test results are determined automatically by the system software. The amount of analyte in a sample is determined from the measured light production by means of a stored nonlinear calibration curve. Patient test results can be reviewed using the Sample Results screen. Refer to the Operator's Guide for complete instructions on reviewing results.

## 9. REPORTABLE RANGE OF TEST RESULTS

Results are reported to the nearest tenth (0.1). The lowest reportable TgAB result is 0.9 IU/mL. Results above the top standard (generally near 2500 IU/mL) are reported as "greater than" (e.g. >2500 IU/mL). Estimates of imprecision can be generated from long-term quality control pool results.

## 10. QUALITY CONTROL (QC) PROCEDURES

- A. Bench quality controls are used in this analytical method. Bench quality control specimens are tested with each analytical run (a set of consecutive assays performed without interruption) so that judgments may be made on the day of analysis. The data from these materials are then used in estimating methodological imprecision and in assessing the magnitude of any time-associated trends.
- B. The bench controls are purchased in sufficient quantity to provide serum samples for all the assays for approximately 1 year. Ranges are established after 20 parallel runs with previously established controls. The quality control pools comprise three levels of concentration spanning the low, borderline and high ranges for TgAB.
- C. Bench quality controls are placed at the beginning of each analytical run. After analysis, the long-term quality control charts (Levey-Jennings) for each control material are consulted to determine if the system is "in control." The Levey Jennings chart plots the quality control material observations on the y-axis and the date of the observation on the x-axis. Quality control material observations are compared with the 95% and 99% confidence limits as well as with the center line (the overall mean of the characterization runs) prior to reporting any results. The system is out of control if any of the following events occur for any one of the quality control materials:
  - The observation from a single pool falls outside the 99% confidence limits.
  - The observations from two pools fall either both above or both below the 95% confidence limits.
  - The observations from eight successive runs for one pool fall either all above or all below the center-line and the current result is above or below the 95% confidence limits.

## 11. REMEDIAL ACTION IF CALIBRATION OR QC SYSTEMS FAIL TO MEET ACCEPTABLE CRITERIA

If the run is declared "out of control", the system (instrument, calibration standards, etc.) is investigated to determine the root of the problem before any results are released. Consult with the supervisor for appropriate actions.

## 12. LIMITATIONS OF METHOD; INTERFERING SUBSTANCES AND CONDITIONS

- A. The upper reportable value is approximately 2500 IU/mL, and is dependent upon the assigned value of the current lot of calibrator.
- B. The lowest reportable value is approximately 0.9 IU/mL. According to the manufacturer, this is the lowest detectable level of TgAB distinguishable from zero with 95% confidence.

- C. The TgAB results should be interpreted in light of the total clinical presentation of the patient, including: symptoms, clinical history, data from additional tests, and other appropriate information.
- D. This assay does not demonstrate "high dose hook" below 350,000 IU/mL.
- E. According to the manufacturer the following substances do not interfere with the assay:
  - Hemoglobin up to 500 mg/dL
  - Bilirubin up to 40 mg/dL
  - Triglycerides up to 3000 mg/dl
  - Human albumin levels of 6 g/dL
- F. For assay employing antibodies, the possibility exists for interference by heterophile antibodies in the patient sample. Patients that have been regularly exposed to animals or have received immunotherapy or diagnostics procedures utilizing immunoglobulins or immunoglobulin fragments may produce antibodies, e.g. HAMA, that interferes with immunoassays. Additionally, other heterophile antibodies such as human anti-goat antibodies may also be present in patient samples. Such interfering antibodies may cause erroneous results. Carefully evaluate the results of patients suspected of having these antibodies.
- G. **WARNING:** The concentrations of TgAB in a given specimen determined with assays from different manufacturers can vary due to differences in assay methods and reagent specificity. Values obtained with different assay methods cannot be used interchangeably.

### **13. REFERENCE RANGES (NORMAL VALUES)**

0-4.0 IU/mL

Based on manufacturer's studies following the criteria outlined by the National Academy of Clinical Biochemists (NACB) for establishing a normal reference range for thyroid antibody testing. Testing was done on males under the age of 30, with TSH levels between 0.5-2.0 mIU/L, no goiter, no personal or family history of thyroid disease, and absence of non-thyroid autoimmune disease. After completing the screening, 137 samples were tested generating a 95% non-parametric upper reference limit below 9 IU/mL. Additionally, 519 normal samples were collected from males and females in the United States ranging in age from 18-74 years old. The screening criteria included serum TSH levels between 0.5-2.0 mIU/L, no goiter, no personal or family history of thyroid disease, and absence of non-thyroid autoimmune disease. After completing the screening 519 samples were testing. 96% of these samples fell below 4 IU/mL.

### **14. CRITICAL CALL RESULTS ("PANIC VALUES")**

Not applicable to this procedure.

### **15. SPECIMEN STORAGE AND HANDLING DURING TESTING**

Specimens should be maintained at 20-25 °C during testing. After testing, the samples are stored at -70 °C or colder.

**16. ALTERNATIVE METHODS FOR PERFORMING TEST OR STORING SPECIMENS IF TEST SYSTEM FAILS**

There are no acceptable alternative methods of analysis. Specimens may be stored at 4-8 °C for no longer than 2 days. Otherwise, specimens should be stored -70 °C or colder until the system is returned to functionality.

**17. TEST RESULT REPORTING SYSTEM; PROTOCOL FOR REPORTING CRITICAL CALLS (IF APPLICABLE)**

Not applicable to this procedure.


**18. TRANSFER OR REFERRAL OF SPECIMENS; PROCEDURES FOR SPECIMEN ACCOUNTABILITY AND TRACKING**

Standard record keeping should be used for tracking specimens. Samples are inspected upon arrival and new boxes are added to an Excel worksheet (sample log) used to track boxes. This sample log is used to track the status of testing and resulting.

The residual serum is stored at  $\leq -70$  °C for 6 months after analysis, then it is returned to the NHANES Repository for long-term storage.

19. Summary Statistics and QC Graphs

		Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
41542	5	23JAN09	29JAN09	173.640	12.039	6.9
40701	47	23JAN09	29JUL09	15.260	1.049	6.9
40703	47	23JAN09	29JUL09	35.749	2.357	6.6
41572	78	03FEB09	27JAN10	172.081	12.327	7.2
40731	36	05AUG09	27JAN10	42.053	3.356	8.0

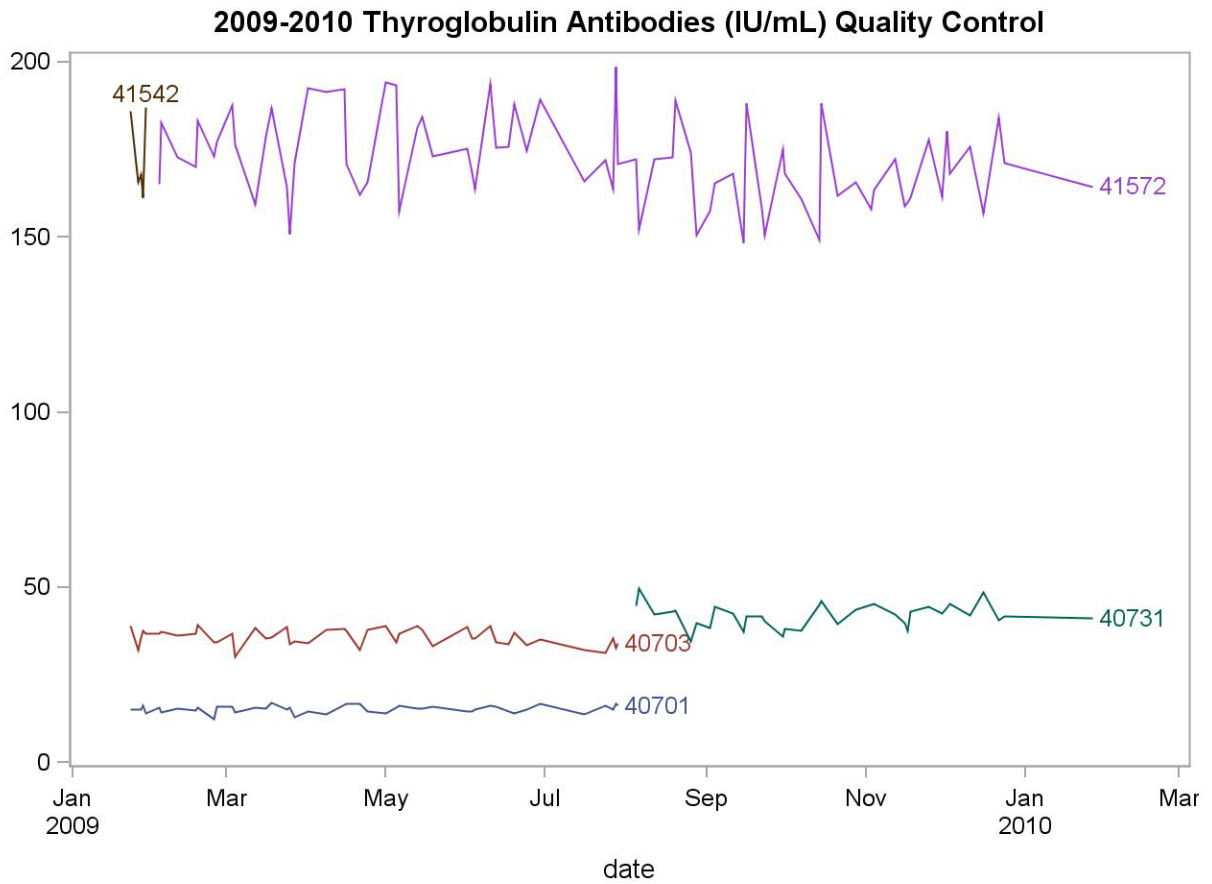

## REFERENCES

### Manufacturer Information:

Beckman Access Immunoassay System Operator's Guide and Reference Manual  
Thyroglobulin antibody II kit inserts, Beckman Coulter, Inc. 2006 Beckman Coulter

National Health and Nutrition Examination Survey (NHANES) Contract Laboratory Manual.  
September 2006

Burtis CA, Ashwood ER, editors. Tietz textbook of clinical chemistry. 2<sup>nd</sup> ed. Philadelphia: WB Sanders, 1994.

Larsen, PR, Ingbar, SH. The thyroid gland. In: Wilson JD, Foster DW, editors. Williams textbook of endocrinology. 8<sup>th</sup> ed. Philadelphia: W.B. Sanders, 1992: 357-487

Demers LM, Specer CA. Laboratory medicine practice guidelines – Laboratory support for the diagnosis and monitoring of thyroid disease. The National Academy of Clinical Biochemistry.

Feldt-Rasmussen, U. Analytical and clinical performance goals for testing autoantibodies to thyroperoxidase, thyroglobulin and thyrotropin receptor. Clinical Chemistry; 42:1 pp. 160-3 (1996)