

Laboratory Procedure Manual

Analyte: **Insulin**

Matrix: **Serum**

Method: **Tosoh Bioscience AIA-900**

As performed by: Diabetes Diagnostic Laboratory
University of Missouri
1 Hospital Dr. Columbia, MO 65212

Contact: Dr. Randie Little/Curt Rohlfing/Rhonda Howard

Important Information for Users

The University of Missouri, Diabetes Diagnostic Laboratory periodically refines these laboratory methods. It is the responsibility of the user to contact the person listed on the title page of each write-up before using the analytical method to find out whether any changes have been made and what revisions, if any, have been incorporated.

Public Release Data Set Information

This document details the Lab Protocol for testing the items listed in the following table:

File Name	Variable Name	SAS Label
INS_J	LBXIN	Insulin (uU/mL)
	LBDINSI	Insulin (pmol/L)

1. SUMMARY OF TEST PRINCIPLE AND CLINICAL RELEVANCE

The AIA-PACK IRI is a two-site immunoenzymometric assay which is performed entirely in the AIA-PACK. This assay is designed for IN VITRO DIAGNOSTIC USE ONLY for the quantitative measurement of Insulin (IRI) in human serum on Tosoh AIA System analyzers. Insulin present in the test sample is bound with monoclonal antibody immobilized on a magnetic solid phase and enzyme-labeled monoclonal antibody in the AIA-PACK. The magnetic beads are washed to remove unbound enzyme-labeled monoclonal antibody and are then incubated with a fluorogenic substrate, 4-methylumbelliferyl phosphate (4MUP). The amount of enzyme-labeled monoclonal antibody that binds to the beads is directly proportional to the insulin concentration in the test sample. A standard curve is constructed, and unknown sample concentrations are calculated using this curve.

Insulin, the antidiabetic hormone, is produced in the pancreatic β cell as^{1,2} a large preproinsulin containing 109 amino acid residues with a molecular weight of approximately 11,500. This peptide is rapidly converted by cleavage to proinsulin consisting of 86 amino acid residues with a molecular weight of approximately 9,000 and is stored within the B cell secretory granules.³ Equimolar quantities of insulin (51 amino acids with a molecular weight of approximately 6,000) and C-peptide (31 amino acids, molecular weight of approximately 3,000) are produced through proteolytic cleavage and then secreted along with a small amount of proinsulin.² Insulin exists in polymeric forms depending on the pH and zinc content. The monomeric insulin molecule is composed of two polypeptide chains, "alpha" and "beta" which are connected by two interchain disulfide bridges of cystine.³ Insulin is degraded in most tissues and has a plasma half-life of 7 - 15 minutes in man. It is also rapidly and completely inactivated in the gastrointestinal tract. Daily production of insulin in a healthy adult is 40 - 50 units. After binding to its specific receptors on target cell membranes, insulin acts as an anabolic and anticatabolic hormone, influencing the rates of carbohydrate, lipid, and protein and electrolyte metabolism.² Insulin release is stimulated by glucose. A failure to respond to this glucose stimulus may be one of the fundamental defects in human diabetes.^b The factors stimulating or inhibiting insulin release or the factors decreasing the tissue response to insulin have been well documented.

2. SAFETY PRECAUTIONS

All specimens, control materials, and calibrator materials are considered potentially infectious. All staff are required to use personal protective equipment including gloves, lab coats, and eye protection when handling bodily fluids (substances). Staff not following proper PPE measures will receive a verbal warning. If a second offence occurs, this staff member will be written up. If this occurrence continues, the lab director and manager will schedule a meeting with the staff member to discuss and prevent additional occurrences. All technical staff are required to read and follow assigned Lab Safety Procedures. Safety Data Sheets (SDSs) and the Chemical Inventory Lists are located at the bench for each reagents used in the lab. All lab staff is required to complete mandatory on-line safety trainings provided by MU Environmental Health & Safety (EHS) department.

Disposal of all waste material should be in accordance with local guidelines.

All Sharps must be placed in rigid puncture proof containers. Each container should be labeled with the biological hazard system. These containers should be kept in an upright position. Training provided by EHS on Biohazard Containers should be performed by each testing personnel at least annually.

3. COMPUTERIZATION; DATA SYSTEM MANAGEMENT

Data are maintained on a secured Microsoft Access / Microsoft SQL server client-server system in a 128-bit authenticated Windows domain environment.

- A. Laboratory services are requested through the Westat system operations via an email notification containing a unique manifest list of the samples and sample analysis type (e.g. FGL), which confirms that specimens have been shipped to DDL.
- B. Each Manifest Form should include and be verified against each sample received:
 - 1) Patient Sample ID #
 - 2) Test Name
 - 3) Date Collected
 - 4) Shipment ID #
 - 5) Shipment Date

- 6) Lab Name
 - 7) Lab ID
 - 8) Survey Year
- C. Once specimens are received and verified the corresponding file is imported electronically into the SQL server database via secure transfer.
 - D. After analysis the results, date analyzed and tech initials are imported from the instrument into the SQL server database via secure transfer.
 - E. Data check sheets are printed out and checked against the instrument printouts by the supervisor.
 - F. After results are cleared by the supervisor a results file in the specified format is exported and uploaded to Westat via secure transfer.

**4. SPECIMEN COLLECTION, STORAGE, AND HANDLING PROCEDURES;
CRITERIA FOR SPECIMEN REJECTION**

- A. Specimen Collection – Serum collected in Red Top Tubes is required for the assay. EDTA, heparinized and citrated plasmas SHOULD NOT BE USED. SST or gel tubes have not been validated. No special patient preparation is necessary. A venous blood sample is collected aseptically without additives (Red top tube).
- B. Sample and Reagent Volumes (N/A) – 50 µL required for analysis
- C. Sample Minimums – min volume required for analysis directly from collection tube is 300 µL (ideally 1 mL is preferred).
- D. Storage
Serum insulin in Red Top Tube: 6 hr at room temperature
Serum transferred : 7 days at 2-8 °C, ≤ 60 days at - 20 °C, 10 years at – 70 °C
- E. Specimen Handling – Specimen are received under frozen conditions on dry ice. Each specimen must arrived in the laboratory labeled with two unique accession identifiers generated by NHANES. Each Manifest Form should include be verified against each samples received. Once checked all samples should be storage at -70°C until analysis. Refer to the NHANES study handling and reporting procedure for additional details. Specimens should be thawed at room temperature prior to analysis and after analysis samples should be returned to -70 °C for long-term storage.

- F. Specimen Rejection;
 - a. Unlabeled specimens.
 - b. Quantity not sufficient samples.
 - c. Incorrect lab name.
 - d. Incorrect test name.
 - e. Wrong lab ID.
 - f. Other specimen types such as urine, whole blood. Only accept serum samples
 - g. The below must be followed for unacceptable specimens;
 - i. Document all details on the sample manifest form and the Pre-analytical Problem Resolution Log.

Immediately contact NHANES for a corrective action on handling the issue.

5. PROCEDURES FOR MICROSCOPIC EXAMINATIONS; CRITERIA FOR REJECTION OF INADEQUATELY PREPARED SLIDES

Not applicable for this procedure

6. EQUIPMENT AND INSTRUMENTATION, MATERIALS, REAGENT PREPARATION, CALIBRATORS (STANDARDS), AND CONTROLS

A. Instrumentation

Automated Enzyme Immunoassay Analyzer (AIA-900); ST IRI Test Code 058.

Refer to Operator's Manual for Startup, Operation, and Shut down procedures.

B. Materials

REF	CONTENT
022930	AIA-900 System
020968	AIA-PACK Substrate Set II
020360	AIA-PACK IRI Calibrator Set
020560	AIA-PACK IRI Sample Diluting Solution
020955	AIA-PACK Wash Concentrate Set
020956	AIA-PACK Diluents Concentrate
020970	AIA-PACK Detector Standardization Test Cups
020971	AIA-PACK Sample Treatment Cup
018581	AIA-PACK Sample Treatment Cup

C. Reagent Preparation

Refer to Package Insert for instruction for use details.

D. Calibrators (Standards) – Calibrators are handled and prepared following manufacturer’s instructions.

E. Controls

- 1) At least two levels of controls are purchased Bio-rad Laboratories; Lyphochek® Immunoassay Plus Controls. Refer to manufacturer’s package insert for processing and handling details.
- 2) In-house controls once thawed are ready for use – Controls were prepared by collecting one unit each of whole blood from three or more non-diabetic volunteers. Transfer 500 µL aliquots to polypropylene storage tubes. Cap tightly and freeze at -70°C or colder.

7. CALIBRATION AND CALIBRATION VERIFICATION PROCEDURES

A. When using new calibrator lots, enter the calibrator concentration values and lot number into the software test file (Refer to the Tosoh AIA System Operator’s Manual for details).

- 1) Calibrators for ST AIA-PACK IRI are lyophilized and require reconstitution prior to use.
- 2) Load the appropriate amount of ST AIA-PACK IRI test cups on the instrument.
- 3) Add the appropriate amount of each calibrator to sample cups.
- 4) The ST AIA-PACK IRI Calibrator Set contains assigned concentrations of insulin. The assigned value is determined on a lot-by-lot basis and is designed to provide an assay calibration range of 1.0 to 320 µU/mL. of insulin.
- 5) Verify that both the calibrator lot and concentration numbers have been correctly entered into the software.
- 6) Calibrators should be run in triplicate measurements.
- 7) Press the CALIB button to order calibration.
- 8) Press the UP or DOWN button to select #IRI on the screen.
- 9) Display the Numeric keypad screen in order to enter the lot No. or use the barcode scanner reading the lot No. from the barcode label.
- 10) Move the cursor to CAL.

- 11) If the correct lot No. and concentration values are displayed, press OK button. The calibration program for the analyte will appear on the ORDER (NON-BARCODE) screen. Each calibrator is run in triplicate.
- 12) Press EXIT button back to HOME screen.
- 13) After placing the sample racks loaded with calibrators and test cups on the sample loader, press the (ASSAY START-NON) button on the HOME screen. The END OF REQUEST will display the number of calibrator.
- 14) Press START button to start the assay. The display will return to the HOME screen.

Calibration verification is performed at a minimum when a system is first placed in service and every 90 days. Recalibration is required (regardless of the length of time since last performed) immediately if any of the following occurs:

- a. A change of reagent lots for chemically or physically active or critical components, unless the laboratory can demonstrate that the use of different lots does not affect the accuracy of patient/client test results or the range used to report patient/client test data.
 - b. If QC materials reflect an unusual trend or shift or are outside of the laboratory's acceptable limits, and other means of assessing and correcting unacceptable control values fail to identify and correct the problem.
 - c. After major maintenance or service. The Laboratory Director must determine what constitutes major maintenance or service.
 - d. When recommended by the manufacturer.
- A. Calibration Acceptability Criteria
- a. The mean rate for the zero calibrator should be <math><3.0\text{ nM/sec}</math>.
 - b. Since there is a direct relationship between concentration and rate, the rates should increase as the concentration increases.
 - c. The replicate values should be within a 10% range.
- B. Calibration Review and Acceptance
- a. Upon completion of the assay, review the calibration curve carefully using the criteria listed above.

- b. Edit the calibration if necessary, then accept the calibration curve by following the below steps:
 - Press the (SUB MENU) button on the HOME screen to display the SUB MENU screen.
 - Press the (CALIBRATION) button to display the CALIBRATION REVIEW (PENDING) SCREEN.
 - Select the calibration to be reviewed and press the REVIEW button.
 - Once the assay results have been confirmed, press the (CALCULATE) button, the Calibration Curve Graph will appear on the screen.
 - Press the (ACCEPT) button to accept the calibration. If the calibration curve of the same lot already exists, the previous curve will be overwritten.
 - Press the (PRINT) button to print the calibration curve.

8. **PROCEDURE OPERATING INSTRUCTIONS; CALCULATIONS; INTERPRETATION OF RESULTS**

A. Preliminaries

- 1) Allow frozen samples and controls to reach room temperature. Invert gently to mix.
- 2) Bring the substrate reagent, calibrators if necessary and test cups to room temperature.
- 3) Label sample cups (provided by Tosoh) using a fine-point permanent marker with sample ID or control ID.
- 4) Check the Diluent solution and wash solution bottle level. Fill the bottles as needed.
- 5) Empty the waste bottle and solid waste container.

B. Instrument set-up – (Refer to Chapter 5 of the Operator’s Manual for **start-up** procedure, Chapter 8 for the **Operations** procedure, and Chapter 9 for **SHUT Down** procedure).

- 1) Place Substrate in substrate compartment. Power on. Wait for Log On screen.
- 2) Log on by pressing OK to retain current operator ID. For new operator, press (OPERATER) button then (MODIFY) button to enter your ID. After enter your

ID, press OK. Your ID will be saved on user list. Press OK again and the following screen will appear with your user ID.

- 3) Press OK again and daily check Screen will appear. Following the screen to do the daily, weekly and monthly maintenance.
- 4) Place a cup adapter rack with a Standardization Cup (STD) in position 2 to the instrument. Press OK. Automated maintenance will begin.
- 5) Record results of the Substrate Background Measurement on the Substrate Background form and keep the printed result for troubleshooting.
- 6) Press (RT.OPEN) button (orange color) to refill tip rack. Then press OK.

C. Processing Samples:

- 1) Using an Eppendorf pipette, transfer the required volume of controls and samples into the corresponding sample cups. Remove any bubbles. For each run one set of controls is placed in front and 2nd set is placed at the end.
- 2) Put sample cup into the sample rack with test cup. Place the rack on the instrument using a maximum of eight racks at a time. Put end marker tube in the last hole of the last rack so that the instrument will stop.
- 3) Press ASSAY START (NON) button on HOME screen. The number 0 will be displayed in the box of the END OF REQUEST.
- 4) Press START to start the assay.
- 5) To run more samples repeat step 2 and 4.

D. Procedural Notes

- 1) If other number other than 0 displayed in box of the END OF REQUEST, press the cell to display the NUMERIC KEYPAD screen. Enter 0 than press OK.
- 2) The AIA-900 is provided with a camera which reads the analytic name and lot number printed on the test cup. It enables the Instrument to automatically distinguish one specimen from another and also to recognize the assay to be done for each specimen.
- 3) Any changes to procedure must be documented. Major changes to the SOP may include the way a procedure is performed or calculations and requires the approval from the Medical Director. Minor changes include typo graphical errors or other minor corrections that do not change the

way the procedure or calculation is performed and do not require approval of the laboratory director. Major SOP changes must be reviewed by the Medical Director prior to approval.

- 4) Any changes to the SOP will be communicated to technical staff via verbal communication, email notification and changes to procedures are tracked in Navex. Technical Staff after reading the changes made to the SOP will sign and date their SOP review.

E. Printing after Calibration

The AIA System performs all sample and reagent handling operations automatically. The AIA Systems read the rate of fluorescence produced by the reaction and automatically converts the rate to insulin concentration in $\mu\text{U/mL}$ then both the rate and concentration are printed out.

- 1) If the calibration curve is undetermined before the analysis, use the steps below to recalculate results after a calibration curve is determined.
 - a. Press RESULTS button in HOME screen.
 - b. Place the cursor line under the first sequence number desired.
 - c. Press SELECT button.
 - d. Use the down arrow key to move the cursor to the last sequence number desired.
 - e. Press SELECT again. A green > sign will be on left of the result selected.
 - f. Press FUNCTION button, the screen will appear with 4 options. Select option RECALC then press OK.
 - g. Repeat step 3 to 6 but select option PRINT instead of RECALC then press OK. The result will be printed on the printer tape.

F. Replacement and periodic maintenance of key components

- 1) Perform the Daily maintenance on the day of assay by following the DAILY MAINTENANCE schedule form. Discard the used sample cups when the analysis is complete. Turn off the instrument once every 24 hours.
- 2) Clean the B/F Wash Probe tip weekly with a cotton ball moistened with 70% ethanol.
- 3) Wash the Substrate line weekly.

- 4) Replace the wash probe tip on the end of the B/F wash probe monthly.
- 5) Clean diluents and wash tanks tri-monthly with 1:100 dilution of hypochlorous acid solution then rinse reservoirs with DI water
- 6) AIA-900 performs a substrate background measurement each time daily maintenance is run and the results are automatically printed out. If the substrate background measurement is within specifications, an OK will be displayed next to 4MU Background. If the substrate background is too high a "BH" (blank high) error flag will be printed and Substrate Replacement will be incomplete. Prime or replace the substrate and repeat daily maintenance. If the lamp intensity level is within specifications an OK will be displayed next to Lamp Intensity Level. If the lamp intensity level is too low an "LL" (lamp low) error flag will be printed. The "LL" is warning that the lamp will need to be replaced soon.

G. Calculations

- 1) The AIA Systems perform all sample reagent handling operations automatically. The AIA Systems read the rate of fluorescence produced by the reaction and automatically converts the rate of insulin concentration to $\mu\text{U/mL}$.

H. Recording of Data

- 1) Quality Control Data—All replicate values of quality control data plus all pertinent assay information are recorded in the Tosoh Insulin Assay Log Database located on the network drive. Print out the Tosoh Insulin Diary Sheet.
- 2) Analytical Results—When the assay is accepted, record the results.
- 3) Refer to the NHANES study sample handling and reporting procedure in Navex for additional details.
- 4) During the data entry process, check the lab accession number. Enter any comments associated with the specimen in the comment field.
- 5) Print a data check sheet with test result. Test result is checked against the instrument print out by the supervisor. A copy of the data check sheet is kept in the appropriate book at the Diabetes Diagnostic Laboratory.

9. REPORTABLE RANGE OF RESULTS

A. Analytic Range – 1 – 320 $\mu\text{U}/\text{mL}$

- 1) Limit of Detection (LOD) – 1 $\mu\text{U}/\text{mL}$
- 2) Accuracy & Precision – Refer to the validation studies
- 3) Analytical Specificity – As per Package Insert details
- 4) Linearity Limits – 1 – 320 $\mu\text{U}/\text{mL}$
- 5) Reportable/Non-reportable values:

Results less than 1.0 $\mu\text{U}/\text{mL}$ are reanalyzed for confirmation and then reported as "<1.0" (less than 1.0). Insulin results that exceed 320 $\mu\text{U}/\text{mL}$ should be 1:10 (or 1:100 depending on the re-assayed result) diluted with the ST AIA-PACK IRI Sample Diluting Solution, and reassayed so the diluted specimen reads between 1.0 to 320 $\mu\text{U}/\text{mL}$. Insulin results that exceed 320 $\mu\text{U}/\text{mL}$ are reported as ">320.0" (greater than 320.0) with the added below comment, "Insulin concentration above the linear range for this assay, sample was re-analyzed as a 1:10 (or 1:100 [add the correct dilution factor] that was used) dilution and result was multiplied by 10 (or 100) to yield a final result of (add final result) $\mu\text{U}/\text{mL}$."

10. QUALITY CONTROL (QC) PROCEDURE

- A. QC is used to ensure that a test system is performing accurately. QC aliquots are tested in the same manner as patient specimens and by the same personnel performing patient testing. QC is used to monitor the analytic performance of the testing being performed at DDL and to certify that all results reported by the lab are accurate.
- B. Batch QC: The bench quality control consists of at least three levels of controls, which cover the spectrum of insulin ranges for both normal and diabetic populations. At least two levels of commercial controls are purchased from Bio-Rad Laboratories. Other matrix appropriate controls are prepared in-house and stored in $-70\text{ }^{\circ}\text{C}$ or colder. One vial of each level of control is thawed and used in each assay.
- C. At least twenty runs (in duplicate measurements beginning and at the end of each analytical run) are used to establish target control values and tolerance limits (ranges). For manufacturer controls, the established acceptable limits should fall within manufacturer's defined limits for each level of control. Defined

- acceptability limits specific for each assay are approved by the Lab Co-Director prior to implementation.
- D. Controls are tested with every run. Control results must fall within the specified ranges prior to the release of results. Control results are reviewed daily by the bench tech and lab supervisor/delegee (prior to the release of results) and examined weekly by the Lab Supervisor or delegee (for trends).
 - E. QC chart plots are plotted for the daily means and ranges of the replicates (the difference between the highest and the lowest value of a single control within a run) for each control and compares them to the established target ranges. These are approved (monthly) by the Lab Director or delegee.
 - F. If control results are outside of the acceptable limits, patient results should not be reported. Patient samples should be return to optimal storage conditions until all QC issues are resolved. QC Corrective actions could include repeating of QC (of the same lot), performing a recalibration procedure or changing out the appropriate reagents.
 - G. QC results must be acceptable prior to releasing patient results. If QC results are still unacceptable after performing corrective actions immediately, call the appropriate technical hotline for additional troubleshooting advice and alert the supervisor or delegate immediately.
 - H. All Corrective Action/Preventive Action (CAPA) events concerning out of range controls, calibrators, or rejected runs that occur at the time of testing, must be documented (in the daily diary sheet for that particular assay) and immediately in the communicated (verbally and by email) to the Lab Director or Supervisor. The QC policy found in the Quality Management Program should be used as a guide for additional details.
 - I. For the insulin assay two types of quality control (QC) systems are used in this analytical method: 1) "sample QC" and 2) "batch QC." For sample QC, 2% of specimens are randomly selected and analyzed either within-assay or between-assay for quality assurance purposes. If the difference between duplicates is greater than 20%, the specimen is reanalyzed. Batch QC specimens are assayed at the beginning and end of each analytical run. If the stock of these controls becomes low, another batch is ordered or prepared in time to analyze it concurrently with the current QC materials. The new controls are used only after

- their means and the ranges have been established by performing 20 characterization runs. Daily means and ranges of the controls are calculated from 20 interassay determinations. The bias ranges of the daily means are set at ± 1 SD or the 67% confidence interval (CI); the warning limits (WL) are the ± 2 SD or the 95% CI and the control limits (CL) are the ± 3 SD or the 99% CI. For the daily ranges, the bias limit is the mean + 1 SD with warning and control limits set at the mean +2 SD and the mean + 3 SD, respectively.
- J. After each assay run, all controls are recorded on the Daily Diary Log Sheet. The QC guideline for this insulin assay declares a system as "out-of-control" if any of the following events occur:
- 1) The mean from a single run for a single control falls outside 99% confidence limits (3sd).
 - 2) The means from a single run for both controls fall outside 95% confidence limits (2sd).
 - 3) The means from eight successive runs for a single control fall either all above or all below the mean line. Runs for which the mean falls within 1sd of the established mean are not counted in this trend.
 - 4) The range from a single run for a single control falls above 99% confidence limits.
 - 5) The ranges from a single run for both controls fall above 95% confidence limits.
 - 6) The ranges from eight successive runs for a single control fall above the mean line.
- K. Refer to the Quality Controls' Manual for the current QC materials' used with acceptable ranges. QC acceptable ranges are also located at the bench where the assay is being performed. Please refer to the current QC acceptable ranges if there are concerns prior to releasing patient results and the lab supervisor (or delegate) if you have questions.
- L. Quality Control—Frequency
- 1) In order to monitor and evaluate the precision of the analytical performance, QC is tested in the same manner as patient samples and by the same personnel performing patient testing.

- 2) QC corrective action could include retrieving new controls of the same lot, reconstituting of new QC, performing a re-calibration, or performing analytical precision studies to assess the integrity of the QC in used against the acceptable limits.
- 3) If it is determined that there is a systemic issue with the AIA-system it is the responsibility of the tech to immediately call the technical service hot line at 1-800-248-6764 for additional troubleshooting advise.
- 4) Once all analyzer issues are resolved and QC is acceptable all patient samples from the rejected run should be reanalyzed.
- 5) QC should also be performed after calibration; at least two levels of controls are run in order to accept the calibration curve.
- 6) Controls are also performed after a PM (sampling mechanism changes, maintenance of the wash probe or detector lamp adjustment or change).

11. REMEDIAL ACTION IF CALIBRATION OR QC SYSTEMS FAIL TO MEET ACCEPTABLE CRITERIA

- A. When the QC results fail to meet the acceptable criteria, check the sample cup containing the QC specimen for bubbles and reanalyze the QC specimen.
- B. If QC is still unacceptable perform a recalibration of the assay and reanalyze the QC. QC must be considered acceptable prior to releasing patient results.
- C. If steps above do not result in correction of the "out-of-control" values for QC materials, troubleshoot the instruments and reagents until the system is back "in control". Refer to the above QC procedure section for additional details.
- D. Samples should be placed under optimal conditions during instrument and reagent troubleshooting.
- E. All corrective actions should be documented.
- F. QC must be considered acceptable prior to releasing patient results.

12. LIMITATIONS OF METHOD; INTERFERING SUBSTANCES AND CONDITIONS

- A. Free Bilirubin (up to 17 mg/dL) and conjugated Bilirubin (up to 19 mg/dL) do not interfere with the assay
- B. Lipemia, as indicated by triglyceride concentrations (up to 1,660 mg/dL), does not interfere with the assay. ³

- C. Ascorbic acid (up to 20 mg/dL) does not interfere with the assay.
- D. Protein, as indicated by human albumin concentrations (up to 5g/dL), does not interfere with the assay.
- E. Hemoglobin may interfere with the assay, particular at the low range. Hemolyzed sample should not be used for this assay because hemolysis may falsely lower the value due to insulin degrading enzyme in red blood cells. ^{1,3}
- F. For patients that are being treated with Asfotase Alfa (Strensiq®), elevated test results are likely to occur using this method. Results should be not be used for patients receiving Asfotase Alfa treatment.

Asfotase Alfa is a prescription medicine used for treatment of patients with hypophosphatasia (HPP), a rare genetic metabolic disorder characterized by the abnormal development of bones and teeth because of defective mineralization.

About HPP, from NORD (National Organization for Rare Disorders);

ref: <http://raredisease.org/rare>

An alternative test method should be used that does not utilize alkaline phosphatase technology. All results should be interpreted with respect to the clinical picture of the patient.

Refer to Product Insert for additional details.

13. REFERENCE RANGES (NORMAL VALUES)

- A. Tosoh suggested mean fasting levels for healthy individuals is below 17 uU/mL.
- B. Reference ranges for insulin were updated at the Diabetes Diagnostic Laboratory in January 2013 by combining results from two volunteer groups. A 360 caloric standard meal (Boost™) challenge was performed in February 2009 for fasting and 120 minutes on non-overweight, non-diabetic subjects (n=15, mean age=34, M:F=8:7). A 360 caloric standard meal (Boost™) challenge was performed in October 2012 for fasting, &120 minutes on non-obese, non-diabetic subjects (n=29, mean age =39, M:F=20:9). All participants fasted overnight for at least 10 but no more than 15 hours. Any subject with a BMI greater than 25 kg/m² and fasting glucose greater than 100mg/dl were excluded from the calculation. Refer to validation study binder.³ The means and observed ranges are:

Insulin Reference Range $\mu\text{U/mL}$

	Fasting	120 min
N	44	44
Mean	6	14
Range	1-14	1-35

Proper interpretation of the Insulin results can be difficult: values are affected by many factors, such as body mass index, age and state of nutrition.

Results that are outside of these reference ranges do not necessarily mean the abnormal test result is of clinical significance. This should only be determined by a physician after careful evaluation of the individual person's health record.

14. CRITICAL CALL RESULTS ("PANIC VALUES")

Not Applicable for this Procedure

15. SPECIMEN STORAGE AND HANDLING DURING TESTING

Specimens are stored at -70°C until analyzed. On the day of analysis, thaw the specimens. Mix thoroughly. Upon completion of analysis, refreeze at -70°C .

16. ALTERNATE METHODS FOR PERFORMING TEST OR STORING SPECIMENS IF TEST SYSTEM FAILS

If the instrument is unable to perform the test, the specimens are stored at -70°C until testing is available. NHANES should be notified if the platform will be down for an extended period of time.

17. TEST RESULT REPORTING SYSTEM; PROTOCOL FOR REPORTING CRITICAL CALLS (IF APPLICABLE)

NHANES data files with results are exported from the NHANES database in the specified format and uploaded to Westat via secure transfer weekly.

18. TRANSFER OR REFERRAL OF SPECIMENS; PROCEDURES FOR SPECIMEN ACCOUNTABILITY AND TRACKING

All shipments are recorded on the NHANES Shipping Log upon receipt. Actions taken during the course of analysis, result reporting, and specimen retention are also recorded on the log. Specimens are stored frozen at -70°C or colder after analysis;


specimen locations are recorded according to sequential DDL accession number and box number. After one year specimens may be discarded.

19. SUMMARY STATISTICS AND QC GRAPHS

See following pages

2017-2018 Summary Statistics and QC Chart for Insulin ($\mu\text{U/mL}$)

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
40301	7	19JAN17	28FEB17	18.320	0.642	3.5
40302	7	19JAN17	28FEB17	56.196	1.111	2.0
40303	7	19JAN17	28FEB17	182.167	2.397	1.3
IH13	100	19JAN17	23JAN19	78.098	2.382	3.0
40331	93	07MAR17	23JAN19	17.187	0.692	4.0
40332	93	07MAR17	23JAN19	58.183	2.036	3.5
40333	93	07MAR17	23JAN19	177.749	4.111	2.3

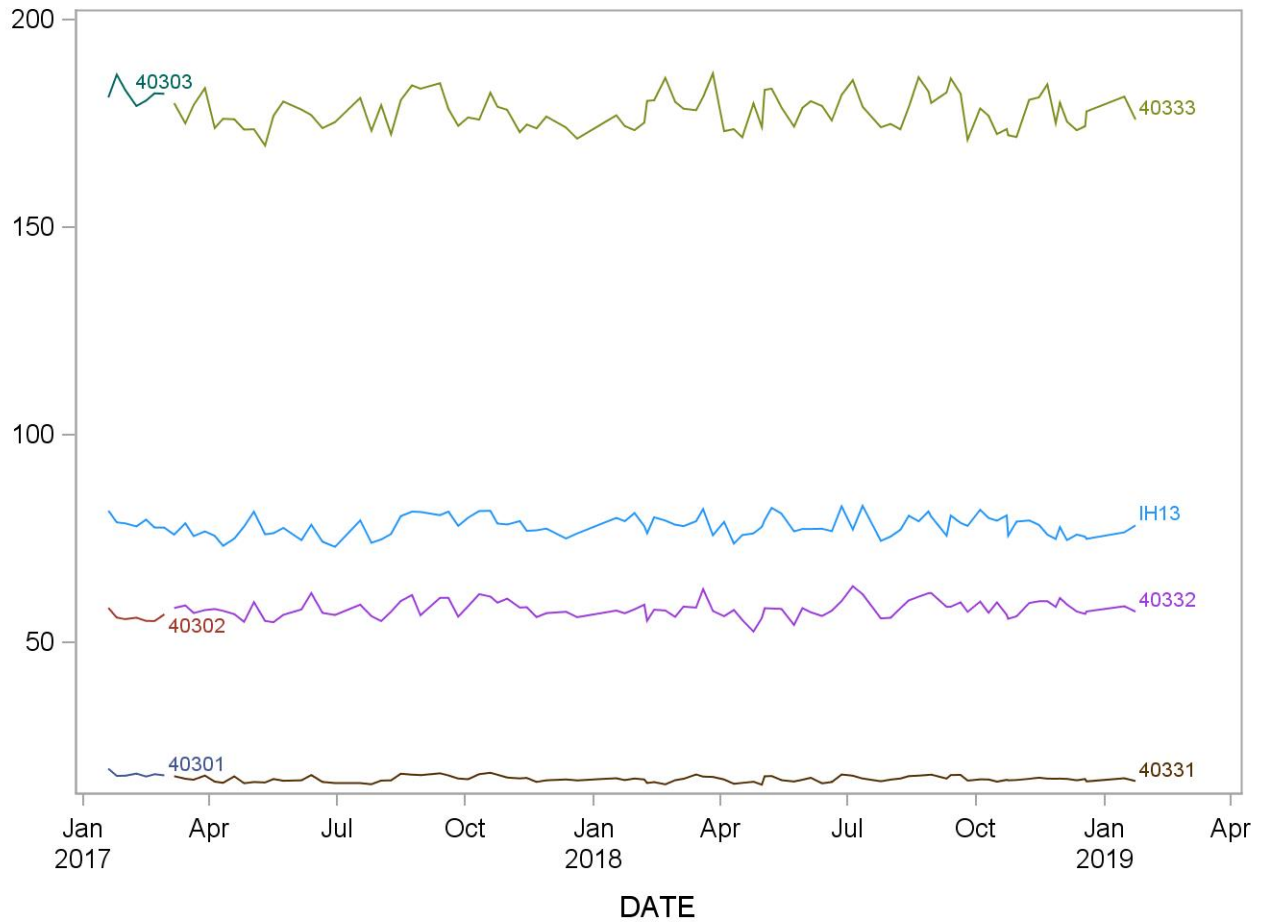

REFERENCES

1. Wrenshall, G.A., et al, 1962, The Story of Insulin. Indiana University Press. Bloomington.
2. Lavy, P.L., 1967, The Pancreatic Beta Cell Structure and Function. *New Eng. J. Med.* 276, 187.
3. Steiner, D.F., 1970, Proinsulin and the Biosynthesis of Insulin. *New Eng. J. Med.* 283, 522.
4. Levine, R., 1970, Mechanisms of Insulin Secretion. *New Eng. J. Med.*, 283, 522.
5. Rieser, P., 1967, Insulin, Membranes and Metabolism. Williams and Wilkins, Baltimore.
6. Ashmore, J., Weber, G. 1968, Hormonal Control of Carbohydrate Metabolism in Liver in "Carbohydrate Metabolism and its Disorders" (ed. Dickens, Randle and Wholan) Academic Press, London, New York.
7. Samols, E., 1965, Immunological Aspect of Insulin on the Nature and Treatment of Diabetes Mellitus. (ed. Leibel, Wrenshall) Excerpta Medica Foundation, Amsterdam.
8. Fajans, SS., J.W. Conn., 1959, The Early Recognition of Diabetes Mellitus. *Ann. N. Y. Acad. Sci.* 82:208.
9. Cerasi, E., R. Luft, 1967, Plasma Insulin Response to Glucose Infusion in Healthy Subjects and in Diabetes Mellitus. *Acta Endocr. (kbh)* 55:278-304.
10. Cerasi, E., et al, 1972, Decreased Sensitivity of the Pancreatic beta-Cells to Glucose in Prediabetic and Diabetic Subjects. *Diabetes* 21: 224-234.
11. Young, D., 1990, Effects of Drugs on Clinical Laboratory Tests. 3rd Edition, Washington, DC, American Association for Clinical Chemistry Press.